


**CADHAM PROVINCIAL LABORATORY**

# **GUIDE to SERVICES**

**2010 Edition**

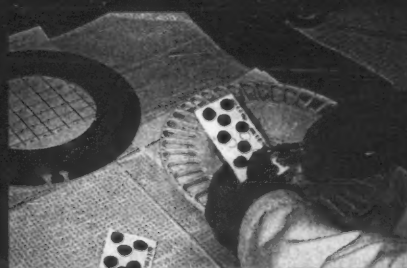

***Serving Manitoba since 1897***

**Manitoba** 

**CADHAM PROVINCIAL LABORATORY**

**MANITOBA HEALTH**

**Location:** Cadham Provincial Laboratory  
750 William Avenue  
Winnipeg, Manitoba

Telephone: (204) 945-6123  
Fax: (204) 786-4770  
Email: [cadham@gov.mb.ca](mailto:cadham@gov.mb.ca)  
Website: [www.gov.mb.ca/health/publichealth/cpl](http://www.gov.mb.ca/health/publichealth/cpl)

**Mailing Address:** Cadham Provincial Laboratory  
P.O. Box 8450  
Winnipeg, Manitoba R3C 3Y1

Revised: January 2008  
February 2009  
February 2010  
(Electronic version only)

## TABLE OF CONTENTS

<b>Responsibilities .....</b>	<b>4</b>
<b>Senior Staff .....</b>	<b>5</b>
<b>Abbreviations Used .....</b>	<b>6</b>
<b>General Guide to Laboratory Use .....</b>	<b>7</b>
<b>Biohazard Response Team .....</b>	<b>11</b>
<b>Infection Control Program .....</b>	<b>11</b>
<b>Outbreak Response Support .....</b>	<b>12</b>
<b>Sexual Assault Protocol .....</b>	<b>15</b>
<b>Needlestick Injury Protocol .....</b>	<b>16</b>
<b>1.0 Technical Support Services .....</b>	<b>17</b>
1.1 Specimen Submission Requirements .....	17
1.1.1 Requisition Requirements .....	17
1.1.2 Specimen Labeling Requirements .....	18
1.1.3 Newborn Screening and Maternal Serum Screening .....	18
1.1.4 Transport of Specimens .....	19
1.2 Specimen Rejection Policy .....	19
1.3 Packaging and Transport of Specimens .....	20
1.3.1 Surface Transport .....	20
1.3.2 Air Transport .....	21
1.4 Transporting Specimens to CPL .....	21
1.5 Transport Supplies .....	22
1.6 Transport Media (TM) .....	23
1.7 Transshipping of Specimens by CPL .....	23
<b>2.0 Clinical Microbiology .....</b>	<b>24</b>
2.1 Specimen Collection .....	25
2.1.1 Abscesses .....	26
2.1.2 Blood for Culture .....	26
2.1.3 Body Fluids .....	26
2.1.4 Bullae, Cellulitis, Petechiae, Vesicles .....	27
2.1.5 Cerebrospinal Fluid .....	27
2.1.6 Cervix and Endometrium .....	28
2.1.7 Conjunctiva .....	28
2.1.8 Nasopharynx .....	29
2.1.9 Nose .....	29
2.1.10 Pus .....	29
2.1.11 Skin .....	30
2.1.12 Sputum .....	30

2.1.13	Stool, Feces, Rectal .....	30
2.1.14	Throat .....	31
2.1.15	Transtracheal Aspirate .....	31
2.1.16	Vagina .....	31
2.1.17	Urethra .....	31
2.1.18	Mid-Stream Urine .....	32
2.1.19	Wounds.....	32
2.2	STI Bacteriology .....	33
2.2.1	Gonorrhea .....	33
2.2.2	Chlamydia.....	34
2.3	Reference Microbiology .....	35
2.4	Antibiotic Susceptibility Testing .....	36
2.5	Reporting .....	37
2.6	Clinical Microbiology Turnaround Times.....	37
<b>3.0</b>	<b>Serology - Parasitology .....</b>	<b>39</b>
3.1	Serology Test List .....	40
3.2	Sample Requirements .....	43
3.3	Requisitions .....	45
3.4	Transport .....	45
3.5	Referred Out Serology Tests .....	46
3.6	Parasitology Test List.....	46
<b>4.0</b>	<b>Virus Detection .....</b>	<b>52</b>
4.1	Specimen Requirements .....	53
4.2	Specimen Collection .....	53
4.3	Virology Test List.....	56
<b>5.0</b>	<b>Newborn Screening and Public Health Chemistry .....</b>	<b>58</b>
5.1	Newborn Screening Program .....	59
5.2	Maternal Serum Screening (Quad Testing) Program .....	66
<b>6.0</b>	<b>Information Management .....</b>	<b>67</b>
<b>7.0</b>	<b>Alphabetical Index of Testing .....</b>	<b>71</b>
<b>8.0</b>	<b>Forms and Requisitions .....</b>	<b>122</b>



## **RESPONSIBILITIES**

Cadham Provincial Laboratory (CPL) is responsible for several province-wide public health, reference and diagnostic services.

It is the central public health microbiology reference laboratory for Manitoba and operates an Infection Control Program in support of rural hospitals, long-term care facilities and Manitoba Health disease control programs.

CPL is directly linked to the Communicable Disease Control (CDC) Branch of the Public Health Division in surveillance of communicable diseases and is the principal laboratory participant in outbreak investigations.

CPL is also the sole centre for laboratory services in virology, chlamydiology and infectious diseases serology and serves patients, practitioners and public health units in Manitoba, and some of Northwest Ontario and Saskatchewan.

CPL is the sole provider of neonatal metabolic disease screening and maternal serum screening for Manitoba and Northwest Ontario.

CPL co-ordinates and conducts the evaluation of suspicious packages and substances for biohazardous materials for Manitoba.

CPL also participates in the training of physicians, nurses and graduate students.

## SENIOR STAFF

### Administration

Laboratory Director: P. Van Caeseele, MD FRCPC	945-6456
A/Administrative Director	945-6302
Privacy Officer	945-6456
Provincial Epidemiologist Nurse Co-ordinator	945-6685
Administrative Officer	945-6337
Education Co-ordinator	945-6230
Outbreak Co-ordinator	945-7473
Safety Officer	945-8021

### Clinical Microbiology

Chief Technologist	945-7184
Scientist; Classical	945-7278
Scientist; Molecular	945-7473

### Information Management

Data Entry Supervisor	945-8001
Information Technologist	945-6865

### Newborn Screening and Public Health Chemistry

Chief Technologist	945-7980
Scientist	945-8021

### Serology and Parasitology

Chief Technologist	945-7582
Scientist	945-7545

### Technical Support Services

Chief Technologist	945-6230
--------------------	----------

### Virology

Chief Technologist	945-6858
Scientist; Virus Detection	945-6878
Scientist; Virus Research	945-7136

## ABBREVIATIONS USED

2 SP	= Chlamydia transport medium	NAAT	= nucleic acid amplification test
CTM		NAD	= nucleic acid detection
Ab	= antibody	NAT	= nucleic acid testing
AD	= antigen detection	NML	= National Microbiology Laboratory
ADB	= anti-DNAase B	NPA	= nasopharyngeal aspirate
Agg	= agglutination	NPS	= nasopharyngeal swab
AIDS	= acquired immunodeficiency syndrome	NT	= neutralization
ALC	= 70 percent alcohol	PCR	= polymerase chain reaction
ASOT	= antistreptolysin O titre	PFGE	= pulsed field gel electrophoresis
BAD	= bacterial antigen detection	PHA	= passive hemagglutination
C & S	= culture and sensitivity	PHI	= Public Health Inspector
CMV	= cytomegalovirus	PHIN	= personal health information number
CONV	= convalescent	PHN	= Public Health Nurse
CPL	= Cadham Provincial Laboratory	QA	= quality assurance
CQI	= continuous quality improvement	QC	= quality control
CSF	= cerebrospinal fluid	RHA	= Regional Health Authority
CT	= cytotoxicity or chlamydia trachomatis	RDA	= RNA/DNA amplification
DFA	= direct fluorescent antibody	RFLP	= restriction fragment length polymorphism
EDC	= expected date of confinement	RPHA	= reverse passive hemagglutination
EIA	= enzyme immunoassay	RPR	= rapid plasma reagin assay
EM	= electron microscopy	SA	= sexual assault
ESBL	= extended spectrum beta lactamase producing	SARS	= severe acute respiratory syndrome
FBI	= food-borne illness	STAT	= immediately
FVT	= fecal verotoxin	STI	= sexually transmitted infection
GC	= gonorrhea	SAF	= sodium acetate acetic acid formalin
HA	= hemagglutination	TB	= tuberculosis
HAV	= hepatitis A virus	TDG	= transportation of dangerous goods
HBV	= hepatitis B virus	TDGR	= transportation of dangerous goods regulations
HCV	= hepatitis C virus	TI	= (1-2%) tincture of iodine
HI	= hemagglutination-inhibition	TM	= transport medium
HIV	= human immunodeficiency virus	VHF	= viral hemorrhagic fever
HSV	= herpes simplex virus	VISA	= vancomycin-intermediate <i>Staphylococcus aureus</i>
HTLV	= human T-lymphotrophic virus	VRE	= vancomycin resistant enterococcus
ID	= immunodiffusion	VRSA	= vancomycin resistant <i>Staphylococcus aureus</i>
IFA	= indirect fluorescent antibody	VT	= verotoxin
IHA	= indirect hemagglutination	VTM	= viral transport medium
KOH	= potassium hydroxide	WB	= Western blot
LA	= latex agglutination		
LCM	= lymphocytic choriomeningitis		
LGV	= lymphogranuloma venereum		
MoAb	= monoclonal antibody		
MOH	= Medical Officer of Health		
MRSA	= methicillin resistant <i>Staphylococcus aureus</i>		

## **GENERAL GUIDE TO LABORATORY USE**

### **Services**

Cadham Provincial Laboratory provides public health laboratory services that include Microbiology, Virology, Parasitology, Serology, Newborn Screening, Public Health Chemistry and Quality Assurance. Reference services for identification and typing of microorganisms are available to all medical and veterinary laboratories in the Province.

Services are available to all registered medical practitioners and midwives, hospitals, health units, medical officers of health, public health inspectors and other recognized health practitioners. There is a charge for laboratory services to patients not insured by Manitoba Health.

Advice is provided by the senior staff on laboratory issues relating to communicable disease. Staff members may visit hospitals or places where outbreaks are occurring at the request of appropriate authorities.

### **Hours of Operation**

The regular hours of service are 0800 to 1630 hrs., Monday to Friday. The Laboratory is partially staffed on Saturdays, Sundays and statutory holidays.

### **Patient Inquiry Services**

Results are provided to authorized personnel for all telephone inquiries, Monday to Friday, 0800-1630; Saturday, limited service 0900-1400 (urgent requests outside of the regular operating hours will be responded to by the medical staff on-call).

### **STAT Testing:**

**Please ask yourself these questions before making a STAT request:**

- 1. Why is test required "STAT"? Could this wait till the next regular shift in the Laboratory?*
- 2. Will doing this test after hours "STAT" alter the management of the patient?*

STAT testing is available 24 hours a day as follows:

#### **Monday through Friday:**

- STAT testing must be arranged through the appropriate Section of CPL prior to shipment.
- A requisition with the appropriate information and clearly marked STAT (a colored sticker is optimum) must accompany the specimen.
- Prior approval from CPL's medical staff or an Infectious Disease Specialist must be obtained for STAT viral testing.

- Prior approval must be obtained from CPL's medical staff for all remaining STAT testing, except for organ donor emergencies.

**After 4:30 p.m., and on Weekends and Holidays (call back):**

- Call (204) 945-6655 and the on-site Security Guard will refer the call to the medical staff on call.

**Specimen Delivery**

Specimens may be delivered at any time, but may not be processed until the next business day if delivered after 1600 hrs.

**Specimen Hazards**

Specimens that break or leak during transport pose a serious physical and infection risk to staff that transport, receive or process them. All specimens sent to the laboratory must be properly packaged and transported. Refer to Packaging and Transport of Specimens (see section 1.3). If the shipper's location and/or the patient specimen can be identified without peril to staff, CPL will notify the sender.

**Reporting Procedure**

Positive results of epidemiological importance or which are likely to be required with urgency by the physician are telephoned with hard copy reports to follow.

Reports issued are for the information of medical or public health staff only. **Please do not have patients contact CPL for interpretation of results.**

**Alert/Critical Results Call Practice**

Preamble: Listed results will be telephoned and/or reported via fax or the CPL PDN system to the physician or other clinical personnel responsible for the patient's care.

**1.0 Virology:**

- All STAT Results
- CMV PCR positives
- Positive preliminary enterovirus results on patients with myocarditis, pericarditis, CNS symptoms, newborns or pregnant women.
- Outbreak results
- Positive results from immunocompromised patients

- Positive rapid RSV results
- Requests for telephone results

## 2.0 Clinical Microbiology:

- Methicillin resistant *Staphylococcus aureus*
- Vancomycin resistant Enterococcus
- *Legionella*
- *Bordetella pertussis* and *parapertussis*
- *Corynebacterium diphtheriae*
- From all sterile fluids and sterile sites (e.g. CSF, pleural fluid, etc.):  
all positive direct smears  
all isolates - preliminary and final.
- All enteric pathogens - preliminary and final.  
**EXCEPTION** - multiple simultaneous samples on the same patient.
- *Blastomyces dermatitidis* - presumptive from KOH and culture, and final confirmed results.
- *Clostridium difficile* toxin positive results
- Verotoxin positive results
- Any unusual or high profile isolates, e.g. a suspected risk level 3 organism.
- Any results specifically requested to be phoned or 'STAT' results.



### 3.0 Serology

- Needle stick on request
- Organ donor results
- Positive results for:
  - Measles or Rubella IgM
  - Hepatitis A - IgM
  - New syphilis (only during outbreaks)
  - Hanta IgM
  - SARS antibody

### 4.0 Newborn Screening and Public Health Chemistry

- Newborn Screening  
All samples exceeding critical limits are referred to a pediatric geneticist or endocrinologist.
- Maternal Serum Screening  
Amnios > 2.0 MoM

### 5.0 Reportable Diseases and Conditions After Hours to MOH On Call

- |   |                               |
|---|-------------------------------|
| • Anthrax   | • Plague                      |
| • Botulism  | • Polio                       |
| • Cholera   | • Rabies                      |
| • Diphtheria  | • <i>Salmonella typhi</i>     |
| • Food poisoning (caused by <i>Bacillus cereus</i> and other unspecified) | • SARS or SARS Ab             |
| • <i>Haemophilus influenza</i> , invasive                                 | • Smallpox                    |
| • Measles and Measles IgM   | • Tetanus                     |
| • Meningitis (other bacterial)  | • Typhoid Fever               |
| • <i>Neisseria meningitidis</i> , invasive                                | • Viral Hemorrhagic Fever     |
|   | • Viral Meningitis outbreaks  |
|   | • Western Equine Encephalitis |
|   | • Yellow Fever                |

## **BIOHAZARD RESPONSE TEAM**

Cadham Provincial Laboratory and the Office of the Chief Medical Officer of Health co-ordinate and conduct the evaluation of suspicious packages and substances for biohazardous materials for Manitoba.

Any spill or suspicious package/substance response first requires triage of the event through the regional or on-call Medical Officer of Health (MOH). **CPL only responds to requests from MOHs in this regard.** The MOH may be reached at the regional public health office or after hours at (204) 788-8666.

The Manitoba team members for Health Canada's Emergency Response Assistance Plan for biosafety level 4 material are located at CPL and may also be reached at (204) 788-8666.

## **INFECTION CONTROL PROGRAM**

Cadham Provincial Laboratory provides an Epidemiologist Nurse Co-ordinator to assist the rural hospitals and all long-term care facilities in Manitoba to establish and maintain infection control programs. Consultation on infection control issues, training for facility infection control practitioners and continuing education for facilities on infection control and infectious diseases are available. Liaison with public health, urban hospitals and community agencies and professional organizations at regional, provincial and national levels complements the program and provides a valuable interface.

## OUTBREAK RESPONSE SUPPORT

Cadham Provincial Laboratory provides laboratory support to Public Health and health care facilities in the investigation of outbreaks.

Note: When submitting specimens for outbreaks, include any history available and transport immediately and directly to CPL. Utilize CPL expertise to ensure appropriate specimens are collected.

An outbreak is the occurrence in a defined area of cases of an illness with a frequency clearly in excess of normal expectancy. The number of cases indicating presence of an outbreak will vary according to the infectious agent, size or type of population exposed to the disease, previous experience or lack of exposure to the disease, and time and place of occurrence. Therefore, the status of an outbreak is relative to the usual frequency of the disease in the same area among the same population, at the same season of the year.<sup>1</sup> Commonly, outbreaks involve gastrointestinal illness, with or without a food-borne component, respiratory illness or parasitic infestation. Antimicrobial-resistant organisms, primarily as colonization, may also be investigated as an outbreak.

Knowledge of circulating pathogens and early detection of new or re-emerging organisms is paramount to disease prevention and control. CPL contributes to this general surveillance by electronic reporting of reportable diseases to Public Health, by informal reports at Infectious Disease reviews, by surveillance programs and by consultations with section staff, the Epidemiologist Nurse Co-ordinator and/or the CPL Director. This listing of reportable diseases can be found on-line at: [www.gov.mb.ca/health/publichealth/cdc/index.html](http://www.gov.mb.ca/health/publichealth/cdc/index.html).

Reports of suspected or actual outbreaks come to CPL from a variety of individuals who may include infection control practitioners, public health nurses, medical officers of health (MOH), public health inspectors (in the case of food-borne illness) and occasionally concerned citizens. Response to this information will vary depending on the nature of the outbreak:

- 1) Additional or rapid testing may be done on specimens after discussion between laboratory and public health care facility staff.
- 2) When an outbreak involves or is anticipated to involve numerous individuals, the MOH or designate may request an outbreak code to be applied to all outbreak samples. This code enables samples to be traced more easily and provides phone reports of positives and written reports

of negatives and positives on a daily basis to the MOH. To obtain a code, the MOH or designate contacts the CPL Outbreak Co-ordinator with a summary of the outbreak and the type of testing desired i.e., virology, bacteriology, toxin detection, serology, or parasitology. Advice regarding appropriate investigations and specimens is available from CPL medical staff or section chief technologists. The MOH or designate then requests the facility, public health nurse or inspector to add the code to all requisitions related to the outbreak.

- 3) During an outbreak, clear communication amongst disciplines is essential. CPL may refer the health care facility/individual to Public Health, may call Public Health directly or may do both. Effort is made to expedite delivery of specimens and to attempt to have laboratory results available to the Public Health Outbreak Co-ordinator at the earliest opportunity. Diseases that are reportable by provincial regulation are reported electronically to the CDC Branch of Manitoba Health on a daily basis. In food-borne illness outbreaks, liaison between CPL and an environmental testing laboratory may occur when both environmental and human specimens are being tested.

Outbreak investigation usually involves testing approximately six affected individuals and not more than ten. Specimens should be taken from appropriate sites and placed in appropriate containers and transport medium when applicable (see sections 2, 3, 4 or 7). They should be sent as soon as possible to CPL in appropriate transport containers and under appropriate conditions (see Section 1). Specimens which leak, are damaged or lack identification (name and requisition number) cannot be processed. Requisitions must be filled out **completely** and where applicable, clearly indicate the outbreak code. Adding "outbreak" or "food-borne illness" if applicable is also helpful.

Supplies for collection and transportation of specimens may be obtained from CPL (see Section 1.5).

Outbreaks involving antimicrobial-resistant organisms e.g., MRSA, VRE, are usually managed by the infection control practitioner(s) of the health care facility in collaboration with the CPL Epidemiologist Nurse Co-ordinator. Arrangements may be made for tracking of specimens and reporting of results as well as molecular epidemiologic investigation. Protocols for the control of MRSA, VRE and ESBLs can be found online at <http://www.gov.mb.ca/health/publichealth/cdc/fs/reportable.pdf>.

Consultation regarding infection control issues during outbreaks is available from CPL by calling the Epidemiologist Nurse Co-ordinator.

In outbreaks where the Manitoba Provincial Outbreak Response Plan (ORP) is activated, CPL will usually be a participant on the team.

Details of pathogen specific outbreak response protocols may be found in the Manitoba Health Communicable Disease Management Protocol Manual, available online at <http://www.gov.mb.ca/health/publichealth/cdc/protocol/>.

1. Heymann, David L. (Editor). *Control of Communicable Diseases Manual, 18th Edition*. American Public Health Association, Washington DC, 2004.

## **SEXUAL ASSAULT PROTOCOL**

CPL is the laboratory-co-ordinating site for the investigation of infectious diseases transmitted during sexual assault.

Consult with local/regional protocols for detailed procedures.

The following infectious agents may be considered in the investigation of sexual assault:

- HIV
- HBV
- HCV
- Chlamydia
- Gonorrhea
- HSV
- HAV
- (others, depending on circumstances)

In all cases, investigating practitioners must take precautions to definitively label each requisition and patient specimen container with the patient name, date and site of collection. Requisitions must also be labeled:

### **SEXUAL ASSAULT**

Specimens/requisitions not labeled as above may be rejected or discarded within two months of testing, limiting evidence available to establish transmission of infection.

Requisitions may be folded to protect the patient's identity and privacy.

Delivery should follow a chain-of-custody protocol.

**NOTE: The abbreviations 'S.A.' or 'SAP' are not acceptable.**



## NEEDLESTICK INJURY PROTOCOL

Please see the Manitoba Health's Integrated Post-Exposure Protocol located at [www.gov.mb.ca/health/publichealth/cdc/index.html](http://www.gov.mb.ca/health/publichealth/cdc/index.html).

There is a small window of time available to exposed individuals to optimize their chances of preventing HIV or HBV transmission after needlestick injury. HIV, HBV and HCV related testing can be conducted on a STAT basis if required. This requires co-ordination with CPL.

Weekdays - tests performed on specimens received before 1430 hours.

Same-day testing - advise source laboratory to transport specimen to CPL STAT.

Weekend or after hours testing - if required, call CPL Security at 945-6655 to page physician on-call.

Results - source and exposed patient name, and CPL requisition number are required.

## 1.0 TECHNICAL SUPPORT SERVICES

The Technical Support Services Section is responsible for:

- Continuous Quality Improvement (CQI), Quality Assurance (QA) and Quality Control (QC) for CPL.
- verification of all diagnostic kits and microbial identification systems for CPL and the rural diagnostic laboratories.
- media preparation and QC.
- co-ordination of CPL's accreditation program
- specimen receiving and processing, mail room services, and wash-up / sterilization.
- acting as a resource for CQI, QA and QC activities for the medical community in Manitoba.

### 1.1 SPECIMEN SUBMISSION REQUIREMENTS

Specimens collected and transported to CPL require the following information:

#### 1.1.1 REQUISITION REQUIREMENTS

Each specimen must have its own requisition.

##### Patient Information (Mandatory)

1. Patient surname, given name.
2. Patient address (street and number), city/town, postal code.
3. Personal Health Identification Number (PHIN); if no PHIN use MB Health number or unique organizational identifier (for out-of-province), etc.
4. Date of birth.
5. Gender.

##### Ordering Practitioner's Information (Mandatory)

6. Ordering practitioner's last name, first name or first initial (full name preferred).
7. Ordering practitioner's reporting address (street and number), city/town, postal code.
8. Ordering practitioner's phone number.
9. Secure fax number.

##### Specimen Information

10. The source and type of specimen.
11. Specimen collection date and time.
12. Collected by and phone number

13. Requested test(s) or procedure(s).
14. Medical/Clinical information (symptoms, history, diagnosis, risk factors, etc.).
15. For referral isolates, suspected organism and previous identification test results are required.

#### **1.1.2 SPECIMEN LABELING REQUIREMENTS (SPECIFIED ON THE CONTAINER) (MANDATORY)**

- Patient's surname and given name.
- As well as one other unique identifier, i.e., PHIN, DOB, MB Health number, requisition number (if there is one). PHIN is preferred.
- Collection date and time.
  - **Do not use addressograph label on specimen**
- For non-nominal HIV and Retrovirus specimens, use the patient code and requisition number as the two unique identifiers.

**If specimens are received at CPL that are inappropriately labeled or packaged, the sender will be notified by report. Repeat errors will result in non-processing of specimens.**

#### **1.1.3 NEWBORN SCREENING AND MATERNAL SERUM SCREENING REQUIREMENTS**

Newborn Screening Specimen requirements: Specimen collection instructions are given on the back of the collection card.

Newborn Screening requires a blood spot card which **must** contain the following information:

- Mother's name
  - Infant's gender
  - Birth weight and
  - Collection date and time, (use 24-hour clock i.e. 4:00 p.m. = 16:00 hr.).
- For proper collection see Section 5.1.

Maternal Serum Screening Specimen Requirements: A minimum of 0.5 ml of serum is required at 16-18 weeks of gestation for testing. **See section 5.2.** A 10 ml serum separator tube is required.

Maternal Serum Screening additional information required for accurate interpretation of results:

- Gestational age (ultrasound information most accurate)
- Ethnicity (race)

- Patient weight at time of phlebotomy
- Insulin dependant diabetes mellitus, (IDDM) status
- Multiple gestation of current pregnancy.

#### **1.1.4 TRANSPORT OF SPECIMENS**

- Must comply with Transport of Dangerous Goods (TDG) Regulations for ground transport and the International Air Transport Association (IATA) Regulations for air transport. (see Section 1.3)
- When required, serum must be separated from the clot in properly processed serum separator tubes or by aliquoting serum into a separate tube.
- Check blue box transport instructions.
- For priority processing: clearly mark "STAT" on the outside of the package and on the requisition.

The laboratory physician on call must be consulted regarding after hour or weekend requests for STAT testing.

#### **1.2 SPECIMEN REJECTION POLICY**

Specimens recieved at CPL will be rejected for analysis for the following reasons:

##### **1. Specimens that cannot be safely processed:**

- i.e. Specimens with needle attached

##### **2. Improperly transported specimens:**

- Improperly packaged
- Specimens transported at incorrect temperatures
- When clinically relevant, collection time is missing
- Time-sensitive specimens received too late for processing

##### **3. Improperly labeled specimens:**

- Two unique identifiers are not available on sample, e.g., name and PHIN
- Patient's name on the specimen container does not match the requisition
- Insufficient information on requisition to send the report
- Specimens from more than one collection time are submitted without collection date and time on all specimen containers and requisitions.

##### **4. Improper specimen collection:**

- Specimen type or source is inappropriate for analysis
- Insufficient specimen for analysis
- Specimen collected in incorrect container(s) or preservative(s)

Where possible, the ordering practitioner will be notified of the reason(s) for rejection of a specimen.

### **1.3 PACKAGING AND TRANSPORT OF SPECIMENS**

All specimens submitted to the laboratory for testing must be packaged in such a manner as to prevent the spillage, breakage, or damage to the specimen itself, and/or to accompanying specimens. The safety of the environment, and the safety of all persons involved in the shipping, handling and receiving of these specimens must be ensured by preventing exposure to the contents of the shipment at any time.

In Canada, the Transport of Dangerous Goods Act regulates how dangerous goods, (including Class 6.2 Infectious Material) may be transported.

The packaging of diagnostic and infectious specimens should generally follow the triple packaging system consisting of a primary receptacle, a secondary packaging, rigid outer packaging and for liquid specimens, an absorbent material placed between the primary receptacle and the secondary packaging capable of absorbing the entire contents of the primary receptacle(s). The appropriate IATA packaging instruction should be consulted for the type of specimen(s) to be shipped and the mode of transportation.

#### **1.3.1 SURFACE TRANSPORT**

##### **Diagnostic Specimens**

A diagnostic specimen is any human or animal material, including but not limited to excreta, secreta, blood and its components, tissue and tissue fluids, that is handled, offered for transport or transported for the purpose of diagnosis. Specimens suspected of containing or known to contain an infectious substance can not be shipped as a diagnostic specimen.

##### **Infectious Substances**

Infectious substances are substances containing viable micro-organisms, including but not limited to, a bacterium, virus, rickettsia, parasite or fungus, or a recombinant, hybrid, or mutant thereof, that are known or reasonably believed to cause disease in humans or animals, and that are included in risk group 2, 3 or 4 of Division 2 of Class 6 dangerous goods.



### 1.3.2 AIR TRANSPORT

#### Category A

An infectious substance which is transported in a form that, when exposure to it occurs, is capable of causing permanent disability, life-threatening or fatal disease to humans or animals. Substances meeting these criteria must be assigned to UN 2814 or UN 2900. Assignment must be based on known medical history and symptoms, endemic local conditions or professional judgment.

#### Category B

An infectious substance that does not meet the criteria for inclusion in Category A. Infectious substances in Category B must be assigned to UN 3373 except for cultures that must be assigned to UN 2814 or UN 2900, as appropriate. The proper shipping name of UN 3373 is Diagnostic or Clinical specimens.

Refer to Blue Box Packaging Directions and CPL Diagnostic and Infectious Specimen Transport Guidelines (see section 8.0).

### 1.4 TRANSPORTING SPECIMENS TO CPL

- All specimens should be shipped by bus, courier or air, whichever is the fastest.
- Specimens sent by bus are picked up at the Winnipeg Bus Depot at 7:45 a.m. and 11:45 a.m. Monday through Friday, and 10:00 a.m., Saturday and Sunday. Urgent specimens may be picked up at other times by prior arrangements. Please call (204) 945-6805 prior to sending urgent specimens.
- Only where bus or other transport is not available should diagnostic specimens be sent via Priority Post. **Specimens must never be sent through the regular mail.**
- All specimens must be clearly labeled and the requisition completely and appropriately filled out.
- Frozen specimen packing should contain sufficient dry ice or frozen gel packs for package contents and distance traveled to maintain specimen integrity.
- Cold, not frozen, gel packs should be used on top or between specimen bags for refrigerated specimens.
- Because of possible contamination, gel packs will be discarded if not received in press and seal bag.
- Ship specimens with similar temperature requirements in the same shipping container.



## 1.5 TRANSPORT SUPPLIES

CPL will provide the following kits/supplies for the collection and transportation of specimens. Orders for supplies are to be sent to CPL via fax at 786-4770 (see form in Section 8.0). CPL will supply forms for ordering of supplies (call 945-6806). Only in emergent situations should this request be telephoned. Maternal Serum Screening (MSS) requisitions are obtained from Genetics at 787-2098.

### DESCRIPTION

### UNIT OF ISSUE

#### Forms

Address Labels - White	Role of 250
General Requisition (Micro, Serology, Virology)	Pkg. of 250
HIV (AIDS) Information Pkg.	Each
HIV Requisitions - General	Pkg. of 25
HIV Requisitions - Viral Load, for specialists or HIV caregivers only	Pkg. of 25
Newborn Screening Specimen Collection Card	Bundle of 50
Newborn Information Pamphlet	Each

#### Containers

Stool Specimen Container - 30 ml	Each
Urine - 60 ml	Each
Glass Slide	Each
Sputum Specimen Container - 70 ml	Each
Ziplock Bag	Each

#### Reagents

VTM	2 ml vial
SAF Stool Preservative	1 litre bottle – MSDS will be supplied
Gram's Crystal Violet	1 litre bottle
Gram's Safranin	1 litre bottle
Gram's Iodine	1 litre bottle
Decolorizer 50/50	1 litre bottle
10% KOH 10 ml	25 ml bottle
40% KOH 40 ml	50 ml bottle
Kovac's	50 ml bottle
Ehrlich's	50 ml bottle
10% Ferric Chloride	50 ml bottle
2SP Chlamydia Transport Media	2 ml vial (special request)

#### Swabs/Kits

Amies Charcoal Transport	Box of 50 or each
Dacron Swabs	Box of 100 or each
Chlamydia/GC GenProbe Aptima:	
Unisex Swab Collection Kit	Box of 50 or each
Urine Collection Kit	Box of 50 or each
Flocked Swabs	Each
Micro Trak	Box of 20 or each

Cardboard mailing tubes, Strep ID kits, styrofoam containers, blood culture bottles and urine dip slide kits are obtainable from Materials Distribution Agency - 945-6040.

## 1.6 TRANSPORT MEDIA (TM)

MEDIA	APPEARANCE	USE	STORAGE
Amies Charcoal TM (with swab)	Black	General T.M. suitable for routine <u>bacterial</u> cultures and sensitivities, especially good for sensitive pathogens ( <i>B. pertussis</i> , <i>N. gonorrhoeae</i> ). Substitute nasopharyngeal swab for included swab where necessary.	Store at room temperature. Do not freeze. Observe expiry date on package
STI Collection Kits for Chlamydia and Gonorrhea	Clear Diluent (1ml)	A rapid NAAT for detection of <i>N. gonorrhoeae</i> and <i>Chlamydia trachomatis</i> from endocervical, urethral and urine specimens.	Store at 2°C-30°C until expiration date on the kit. After specimen collection store at 2°C-30°C. Do not discard swab or buffer.
Sodium Acetate, Acetic Acid, Formalin (SAF) (Parasitology)	Clear no precipitate	2 parts SAF and 1 part stool - thoroughly emulsified at time of collection	Store at room temperature.
Viral Transport Media (VTM)	Clear, straw color	T.M. for swabs, tissues and aspirates requiring viral culture. Not suitable for blood, CSF, urine or stool.	Store frozen at -20°C for up to 12 weeks; 4°C for 1 week. If it is thawed during shipping, it should be refrozen at the earliest opportunity, and may be used for 6 weeks. After collection of specimen, transport to CPL ASAP at 4°C with a cold pack.
2SP Chlamydia Transport Media (Special Request)	Clear	TM for swab specimens requiring PCR for LGV, haemophilis ducreyi, Chlamydia pneumoniae or Chlamydia psittaci	Store frozen at -20°C until used – observe expiry date on vial. After specimen collection, transport to CPL at 4°C.

## 1.7 TRANS-SHIPING OF SPECIMENS BY CPL

In order for CPL to ensure specimens are trans-shipped and received at the appropriate testing site, the appropriate requisition must be completed accurately with the result reporting information (name and location where report is to be sent) clearly visible. **Referrals to the National Microbiology Laboratory (NML) must be processed through CPL.**

## 2.0 CLINICAL MICROBIOLOGY

Clinical Microbiology services involve the detection, isolation and epidemiological characterization of bacterial or fungal pathogens or toxins from clinical specimens. Procedures include, but are not limited to:

- isolation using enriched and selective culture media
- pathogen identification and characterization by biochemical, serological, microscopic and molecular techniques
- antimicrobial susceptibility testing
- toxin testing

The sub-sections include:

- Enteric bacteriology - enhanced foodborne-illness investigation and detection of enteric pathogens.
- Sexually transmitted infections - detection of *Chlamydia trachomatis*, *Neisseria gonorrhoeae*, and other sexually transmitted microorganisms.
- Respiratory - detection of *Bordetella pertussis*, *Legionella* spp., *Corynebacterium diphtheriae*, and other microorganisms.
- Mycology - isolation and identification of yeasts, dimorphic fungi, dermatophytes and some medically significant molds.
- Toxins - detection of verotoxins and *Clostridium difficile* toxins.
- Emerging antimicrobial resistance - screening, detection and characterization of current and emerging antibiotic resistance in organisms such as *Staphylococcus aureus* and *Enterococcus* spp.
- Miscellaneous bacteriology - isolation and identification of pathogens from a variety of clinical specimens and referred-in isolates.
- Molecular typing - molecular typing of bacterial pathogens.

### Service Hours:

Clinical Microbiology provides 7-day service, with some limitations on weekends and statutory holidays. Emergency on-call service is by special arrangement only and must be authorized by the Director or designate.

### Clinical Microbiology programs:

These include, but are not limited to:

- Enteric program including enteropathogenic *Escherichia coli* and provincial screening for verotoxins
- Outbreak and food-borne illness investigation
- Provincial STI screening for chlamydia and gonorrhea
- Prenatal screening for Group B *Streptococcus*
- Respiratory program including Diphtheria toxigenicity, Pertussis and *Legionella*

- Emerging resistance screening - e.g. MRSA, VRE, Carbapenemases
- Molecular diagnostics and epidemiology
- Support of the Provincial Mycobacteriology Program
- Provincial public health reference service
- Laboratory-based surveillance

#### **Referred out services**

- Acts as an intermediary for referring of isolates or specimens to the National Microbiology Laboratory, National Reference Centres, or other appropriate referral and reference laboratories.

#### **Other services**

- Development of and participation in externally and internally funded microbiology research projects.
- Participation in University of Manitoba under-graduate and graduate medical education programs.
- Inservices for Red River College students, Cadham Provincial Laboratory staff, and public health nurses.
- Participation provincially and nationally on issues of public health laboratory and program importance.
- Education and research.

**Note:** This is a general description of services and not meant to be exclusive.

## **2.1 SPECIMEN COLLECTION**

The following is a review of the steps necessary to secure the optimal sample for culture.

**Note:** Sterile gloves should be worn whenever specimens are collected.

**Standard Skin Antisepsis** (best method) for obtaining blood and body fluid specimens:

- After palpation, scrub the site with 70% ALC for a minimum of 30 seconds.
- Apply TI to the area, allowing contact for at least 2 minutes. Let air dry (do not blow).
- Remove the TI with ALC using increasing outward circular movement (2 minutes).
- For superficial lesions such as abscesses and bullae, a gentle disinfection with ALC, allowed to dry, is sufficient.

### **2.1.1 ABSCESSSES**

1. Prepare the surface as per Standard Skin Antisepsis.
2. Aspirate at least 0.5 ml and preferably 1.0 ml of purulent material.
3. Send the specimen immediately to the laboratory. If delay in transportation is anticipated, inject 1.0 ml or more of pus into blood culture medium.
4. Swabs of pus must be placed into transport medium, because of their tendency to dry.
5. Do not freeze. Keep swabs at room temperature (25°C) or 4°C, and keep fluids at room temperature.
6. Expel contents of syringe into sterile screw-capped container.
7. Submitting the contents in a syringe is not recommended, but if absolutely unavoidable, remove the needle and cap and replace with the sterile cap provided. Tape the plunger to avoid spillage.

### **2.1.2 BLOOD FOR CULTURE**

Blood culturing is no longer available.

All blood cultures should be directed to a regional or private laboratory and collected according to that laboratory's procedure.

Special blood culture medium is required for the isolation of mycobacterium. Contact the TB Laboratory at Health Sciences Centre. Telephone: (204) 787-7652.

### **2.1.3 BODY FLUIDS (EXCEPT URINE AND CEREBROSPINAL FLUID)**

1. Prepare the surface as per Standard Skin Antisepsis (See Section 2.1).
2. Use sterile needle and syringe.
3. Handle all specimens so as to ensure viability of potential anaerobic pathogens, i.e., collection into blood culture bottle or an anaerobic transport vial.
4. Send 5 to 10 ml to the laboratory. It is best to transport such specimens immediately to the laboratory, not only to maximize appropriate processing, but to ensure prompt results from immediately available laboratory procedures.



5. Where a delay in processing is anticipated, inject 1 ml into a blood culture bottle and make a smear by spreading one drop of fluid in the center of a clean microscope slide. The smear should be allowed to dry in air and then be fixed over heat. Any remaining fluid should be submitted along with the smear and the inoculated blood culture bottle (25°C).
6. When clotting is anticipated, dilute the sample with sterile saline.
7. Do NOT refrigerate or freeze samples, keep them at room temperature (25°C).

#### **2.1.4 BULLAE, CELLULITIS, PETECHIAE, VESICLES**

1. Prepare the surface as per Standard Skin Antisepsis (See Section 2.1). In the case of bullae and vesicles, care is exercised to avoid lesion rupture.
2. A sterile needle and syringe are used.
3. As much material as feasible is aspirated, and placed in appropriate transport media (bacterial or viral).
4. If no aspirate is available, non-bacteriostatic sterile saline may be injected and aspirated. It is best to attempt this at lesion edges.
5. Petechiae pose special problems, and some prefer excoriation of the skin with a needle tip after vigorous cleansing. In this event, TI should be removed with ALC prior to this exercise. A swab is then used to immediately inoculate chocolate agar plates at the bedside or to put material on a slide for Gram stain.
6. One half of one ml of sterile saline may be injected into the advancing edge of a cellulitis and subsequently aspirated. The aspirate may then be injected into a blood culture medium.
7. Collect at least 0.5 ml, preferably 1.0 ml or more.
8. Do NOT freeze, keep at 4°C.

#### **2.1.5 CEREBROSPINAL FLUID**

1. The physician wears sterile gloves, a gown, and a mask.
2. Prepare the surface as per Standard Skin Antisepsis (See Section 2.1).
3. Drape the surrounding skin with sterile linen.
4. In adults, a needle insertion is ideally followed by collection of more than 2 ml of cerebrospinal fluid into a sterile container for which a leakproof cap is available.
5. Ideally separate tubes are used to collect specimens for a cell count and biochemical analysis.
6. Transport the specimen immediately to the laboratory, as the organisms likely to be isolated are fastidious.



7. Where a delay in processing is anticipated, inject 1 ml into a blood culture bottle and make a smear by spreading one drop of fluid in the center of a clean microscope slide. The smear should be allowed to dry in air. Any remaining fluid should be submitted along with the smear and the blood culture bottle.
8. Store at 35°C or at room temperature (25°C).

### **2.1.6 CERVIX AND ENDOMETRIUM**

1. The patient is placed in the lithotomy position.
2. A speculum is inserted and the cervix is visualized. Excess mucus is removed with a cotton ball or a swab, before the specimen is collected.
3. For cervical cultures, the swab is inserted in the distal portion of the cervical os, and allowed to remain for 10 to 30 seconds. Specimens for chlamydial studies are taken as above, but require rotation of the swab to obtain the superficial layer of cells required for testing.
4. Endometrial cultures should be approached either by needle aspiration or by a double lumen catheter through the cervical os. A slight cut has been previously made in the advancing end of the catheter to allow egress of a number 8 infant feeding tube, which is fed through the lumen of the Foley catheter to obviate normal flora contamination.
5. Place swab in charcoal transport medium for C&S; place chlamydial/GC swab in GenProbe container (see section 2.2).
6. Do NOT freeze or refrigerate.
7. Store at room temperature (25°C).
8. Send as quickly as possible (within 48 hours for C&S).

Note: Only aspirated material is considered to be useful for anaerobic culture. Additional information on other collection kits is provided under the STI Bacteriology and Virology entries, see index.

### **2.1.7 CONJUNCTIVA**

1. Premoisten sterile swab with sterile saline and obtain secretions from inner aspect of eyelid.
2. Transport in charcoal transport medium.
3. For Chlamydia testing, see section 2.2.

### **2.1.8 NASOPHARYNX**

1. The patient is comfortably seated, preferably with the head tilted back.
2. A nasal speculum is gently inserted.
3. A nasopharyngeal swab, on a malleable wire with a Teflon coated non-toxic tip, is inserted parallel to the palate through the speculum into the nasopharyngeal area.
4. The swab is rotated gently and allowed to remain for 20 to 30 seconds. The swab is then removed and placed in transport medium.
5. It is important to stress the use of charcoal transport medium with these specimens because the swab tip is small and vulnerable to drying and the organisms likely to be present are rather fastidious.
6. The specimens should be transported promptly to the laboratory.
7. Store at room temperature (25°C).
8. Do NOT freeze or refrigerate.

### **2.1.9 NOSE**

1. Anterior nares cultures are easily taken with a regular Dacron swab. In small children this is best done with a swab such as described in section 2.1.8.
2. Place swab in Ames charcoal transport medium and send to CPL immediately.
3. Store at 25°C (room temperature).

### **2.1.10 PUS**

1. Aspirate a minimum of 0.5 ml by sterile syringe, if possible, and submit in sterile tube and/or on a swab well-soaked in pus. Send swab in transport medium.
2. To make thin smears, use the swab or by pressing a small spot of pus between two slides and then sliding them apart. Dry in air. Place slides in cardboard slide-mailer and secure with an elastic band.
3. For anaerobic culture, inject pus or other material into a blood culture tube, or into an anaerobic transport system.
4. Store at 25°C (room temperature).

### **2.1.11 SKIN**

1. See also sections: abscesses; bullae, cellulitis, petechiae, vesicles and wounds.
2. For fungi: Cleanse with ALC, scrape the advancing edge of the infected tissue and collect scrapings in dark paper, folded and properly packaged. Keep specimen dry. Do NOT place in transport medium.

### **2.1.12 SPUTUM**

1. Sputum is a very poor specimen unless patient co-operation is assured and unless special laboratory assessment is performed to determine adequacy of the specimen based on the numbers of squamous epithelial cells and leukocytes. Even optimal specimens fail to indicate the causal agent of pneumonia in up to 90% of cases.
2. The patient should be properly instructed as to what is desired. Early morning sputa from the lungs after rinsing the mouth out with water and gargling; removal of dentures and plates is desired. Keep the amount of saliva in the specimen to a minimum. A sterile wide-mouthed, screw capped, leak proof container is provided for the expectorated material.
3. If the patient is unable to produce sputum, induction may be effected by postural drainage, saline nebulization, or chest percussion. Please inform the laboratory by notation on the requisition when this type of a specimen is obtained. Otherwise, it may be mistaken for saliva, and be rejected.
4. Since some of the organisms are fastidious, the specimen should be transported promptly to the laboratory.
5. In infants, tracheal secretions should be submitted.
6. Collect at least 1.0 ml of sputum, but no more than 30 ml.
7. Store at 4°C. Do not freeze or use preservatives.

### **2.1.13 STOOL, FECES, RECTAL**

1. From a clean, urine-free receptacle, transfer stool into a sterile screw-capped container, until at least one-third full and no more than one-half full. Transport promptly.
2. Rectal swab is the specimen of choice for VRE and ESBL, and a suitable alternative for Shigella. Insert sterile swab 1 inch into the anal canal so that feces is evident on the swab. Transport in charcoal transport medium.
3. For C&S, toxin testing or FBI investigation, store at 4°C.
4. Keep from freezing or leaking.

5. If unusual pathogens are suspected (i.e. *Vibrio cholerae*, *Yersinia*, or *Plesiomonas shigelloides*), please indicate this on the requisition. Also indicate whether C&S, C. difficile toxin, or verotoxin testing is desired.

#### **2.1.14 THROAT**

1. In adults, a rayon swab is used with good visualization (use a tongue blade and a good light source). Vigorously swab both tonsillar fauces and the posterior pharynx, reaching up behind the uvula and culturing an ulceration, exudate, lesion or area of inflammation. Place in charcoal transport medium.
2. Store at room temperature (25°C) or 4°C.
3. Submit within 24-48 hours.
4. Indicate on requisition if diphtheria or gonorrhea is suspected.

#### **2.1.15 TRANSTRACHEAL ASPIRATE**

1. This technique is not a routine culture technique and is best done by an experienced individual.
2. Submit specimens in a blood culture media tube or anaerobic transport system, in order to preserve possible anaerobes.
3. Collect at least 0.5 ml, preferably more, in a sterile container.
4. Store at 4°C.
5. Do NOT freeze.

#### **2.1.16 VAGINA**

1. Wipe away excessive secretions; obtain secretions from mucosal membrane of the vaginal vault with sterile transport swab.
2. Intrauterine device may be sent in a sterile container; transport within 24 hours at room temperature (25°C).
3. Vaginal or vaginal-rectal swabs for Group B Streptococcus should be collected at 35-37 weeks gestation.

#### **2.1.17 URETHRA**

1. Use a sterile bacteriologic wire or disposable plastic loop to obtain the specimen from the anterior urethra by gently scraping the mucosa. An alternative to the loop is a sterile rayon urethral swab that is easily inserted into the urethra.
2. Transport of these cultures are as described for cervix and endometrium (see section 2.1.6).
3. For Chlamydia / GC testing, see section 2.2.

### 2.1.18 MID-STREAM URINE

1. It is best to obtain early morning specimens whenever possible. The urine of patients who are receiving forced fluids may be sufficiently diluted to reduce the colony count below  $10^8$  per L.
2. A properly collected mid-stream urine is ideal. An intermittently inserted catheter or suprapubic aspirate is also appropriate. Always use sterile containers.
3. Collect a mid-stream urine as follows: Clean area thoroughly with soap and water and rinse with wet gauze pads; holding the labia apart or the foreskin retracted, begin voiding; after several ml have passed, collect sample without stopping flow.
4. Urine culture - inoculate a urine dip slide with the collected urine as follows: Dip agar slide into the freshly collected urine, allow to drain, replace in container, screw lid tightly. NOTE: Actual urine for culture will not be processed unless the urine arrives at CPL before 4:30 p.m. on the day of the collection.
5. Urine for chlamydia (see Chlamydia in section 2.2.2).
6. Urine for *Legionella* Antigen (see Legionnaires disease in section 7.0).

### 2.1.19 WOUNDS

1. For the closed wound technique, see sections abscesses, bullae, cellulitis, petechiae and vesicles (2.1.1 and 2.1.4).
2. For open wounds:
  - a) Clean the sinus tract opening or the wound surface with normal saline.
  - b) These areas frequently yield "normal flora" organisms. Therefore, it is important to attempt to culture the base or edges of the wound.
  - c) Swab specimens of sinus tracts may be acceptable, but aspiration material obtained by needle or catheterization is preferable. Curettings obtained from the lining of the sinus tract also provide excellent culture material. For ulcerations or open wounds, curettings or biopsy specimens are best. These are placed into a sterile transport container. Tissue or aspirated material provides the greatest yield. Do not freeze. Keep at room temperature (25°C). Transport ASAP.



## **2.2 STI BACTERIOLOGY**

### **2.2.1 GONORRHEA**

#### **Collection**

1. See section on cervical specimens (see section 2.1.6).
2. A rectal swab can be obtained without an anoscope, by inserting a dacron swab approximately 1 inch into the anal canal. The swab is then moved from side to side to sample the crypts and left for 10 to 30 seconds to allow absorption of organisms onto the swab.
3. In the male, urethral culture is usually obtained either with a rayon swab or a plastic loop (see section 2.1.17). Both smear for gram stain and culture are indicated.
4. Throat swabs are occasionally helpful in the diagnosis of gonorrhea, but vaginal swabs generally have a poor yield of positive results.
5. Urine NAAT is available routinely for males, females without a cervix (due to hysterectomy), or those refusing a complete examination. See section 2.2.2 Chlamydia for information on collection and transportation of specimens.

**Note:** Culture for this organism is useful when the specimen can be processed in the laboratory within 24 hours of procurement. It is one way of detecting in the cervix, urethra or eye and the only satisfactory way at present for detection from other sites, i.e. throat, vagina or rectum.

**Note:** If purulent material is present in the urethra or if other lesions are present these should be cultured, in addition to the cervix.



### **Adults**

A swab in Amies charcoal transport medium from anus, throat, eye, vagina, cervix or urethra for culture, where culture can be started same day. Dry swabs are unsuitable.

**or**

NAAT kits (Aptima unisex swab collection kit for endocervical and urethral specimens and Aptima urine specimen collection kit); follow manufacturer's instructions. **Use NAAT where any delay in testing is anticipated.** Use only the swab provided in the NAAT kit, any other type of swab will not be processed. The same kit can also be used for Chlamydia detection. Be sure to leave the extraction fluid in the tube supplied.

### **Prepubertal Children**

Culture is as above, irrespective of length of delay of processing. Unfortunately many *N. gonorrhoeae* infections will be missed if delay is greater than 24 hours. A smear should accompany these specimens. (One streak about 0.5-1.0 inch or 1.5-2.0 cm long is made on a clear glass slide. Air dry). Vaginal swabs in Amies charcoal TM are suitable samples from prepubertal females, but must be cultured.

### **Conjunctivitis in Newborn**

A swab for culture in Amies charcoal TM should be taken. If submitted for NAAT, results will be reported as "for investigational purposes only."

A smear should be submitted.

## **2.2.2 CHLAMYDIA**

### **Collection**

See section on cervical swabs (section 2.1.6) and urethral swabs (section 2.1.17). For urine specimens patient should not have urinated one hour prior to collection. The first 20-30 ml (NOT midstream) should be collected into a sterile plastic preservative-free container, then transferred to the Aptima urine specimen transport tube ASAP. Store at 2°C-30°C until transportation to CPL is available.

### **Tests Available**

The two detection methods readily available at CPL are:

1. NAAT used for urines from males, cervical and urethral specimens.
2. Direct fluorescent antibody (Microtrak) used for all other specimens (e.g. throat, rectal, nasopharyngeal). Culture of *C. trachomatis* requires prior arrangement with the laboratory. This should be used in the instance of treatment failure.

### **Adults**

Collect a cervical swab using the GenProbe Aptima unisex swab collection kit from adult and legally consenting adolescent females.

Collect first void urine using Aptima urine collection kit only on females without a cervix (due to hysterectomy) or those refusing a complete examination.

Collect first void urine (recommended) or a urethral swab using the GenProbe Aptima unisex swab collection kit or Aptima urine collection kit from males.

All other sites (throat, anus) require fluorescent antibody testing (Microtrak). If Microtrak kits are not available, make a smear approximately 1/4 inch (10 mm) in diameter on a clean glass slide and air dry.

See manufacturer's instructions (GenProbe Aptima or Microtrak) for handling swabs, urines and slides.

### **Prepubertal Children**

For vaginal and urethral sites, NAAT testing is preferred; Microtrak testing is acceptable but not preferred. For all other sites, e.g. throat, rectal, etc., Microtrak testing is appropriate. First void (not midstream) urine may also be collected. **Please indicate boldly on requisition that these specimens are from young children.**

### **Eyes**

Eyes can be tested using Microtrak testing or NAAT. Results for eyes tested by NAAT will be reported as "for investigational purposes only."

### **Newborn**

Pulmonary, tracheal secretions and nasopharyngeal aspirates should be submitted in sterile containers. When in doubt as to procedure please phone for advice (945-7204).

## **2.3 REFERENCE MICROBIOLOGY**

Cadham Provincial Laboratory provides a range of specialized reference activities which may assist in the identification or typing of microbes isolated in other clinical laboratories. Material to be examined must be submitted in pure culture in an acceptable manner, bacterial or fungal cultures being on slants of appropriate medium in a tightly capped or stoppered small bottle or stoutly constructed tube. Suspect level 3 organisms should be appropriately packaged and labelled (see section 1.3).

It is essential that the following details accompany the material to be examined: the patient's name, home address, physician, age, sex, details of the original specimen and the illness being investigated, and results of tests already performed in the submitting laboratory.

**The services offered include:**

1. Toxin testing
2. Serotyping
3. Biotyping
4. Molecular typing
5. Nucleic acid amplification techniques

Specimens submitted for tests that are not performed at CPL will be forwarded to an appropriate reference laboratory. Some delay may be expected when this occurs. Acknowledgement is always given when the Laboratory reports results of tests performed elsewhere. E.g., Streptococcal Serotyping, *Clostridium botulinum* I.D., *Chlamydia pneumoniae* I.D., Mycoplasma and Ureaplasma Culture or PCR.

A variety of molecular diagnostic and molecular typing techniques are utilized in the Molecular division of the Clinical Microbiology section of CPL. These tests are typically based on methods developed by other laboratories. The protocols for these methods have either been provided directly by these laboratories or have been published in the peer-reviewed literature and validated in-house prior to being put in service. Any clients of CPL who require more detailed information on these tests (e.g., primer sequences used, PCR reaction conditions) should phone the Scientist in the Clinical Microbiology section for further information.

## **2.4 ANTIBIOTIC SUSCEPTIBILITY TESTING**

1. Antibiotic susceptibility tests are done routinely only on organisms considered to be significant. If testing with particular antibacterial agents is desired, the request should be clearly noted on the requisition. It is not practical to test every available antibiotic.

Routine susceptibility tests are not needed when resistance has not been described to the antibiotic of choice (e.g. *Streptococcus pyogenes* to penicillin).

2. Determination of the concentration of antibiotics in blood or other specimens may be arranged in consultation with the Laboratory.

## **2.5 REPORTING**

Results are sent out by fax or mail. Positive test results which are likely to be required with urgency by the physician are telephoned with hard copy results to follow. Refer to the "Alert/Critical Results Call Practice" in the General Guide to Laboratory Use section.

## **2.6 CLINICAL MICROBIOLOGY TURNAROUND TIMES**

Turnaround time is dependent on a variety of factors including:

- purity of submitted isolate
- fastidiousness and growth requirements of organism
- unusual phenotypic traits of the isolate
- complexity of testing methods required for workup
- amount of isolate information provided by the submitting laboratory

### **General Turnaround Time for a Positive Culture**

- **Bacterial cultures and yeasts** - rapidly growing cultures can often be identified after overnight incubation, and preliminary results may be available within 24 hours of receipt of specimen.
- **Antibiotic susceptibility tests** - 48 hours (because the organism must be isolated in culture before being tested).
- **Fungi, culture and identification** - 2 to 3 weeks for hair, skin and infections; 2 to 6 weeks for others.

**Note:** Referral to reference centres will increase turnaround time.

**Laboratory Tests****Minimum Turnaround Time  
For A Negative Test (working days)**

Actinomyces request	8 days
Anaerobic Culture – miscellaneous specimen	3 days
<i>C. difficile</i> culture & toxin (colon tissue)	5 days
<i>C. difficile</i> toxin (stool)	5 days for tissue culture/same day for rapid test
Chlamydia (by NAAT)	2 days
Chlamydia DFA	2 days
Diphtheria Culture	3 days
Diphtheria Toxigenicity	4 days
Direct Microscopy	24 hrs. or same day
Ear C&S	2 days
ESBL Screen	2 days
GC (NAAT)	2 days
GC Culture	3 days
Genital C&S	2–3 days
<i>Helicobacter pylori</i> Culture	8 days
Legionella Antigen Detection	24 hrs. or same day
Legionella Culture	15 days
Miscellaneous Specimen (wound, ulcer, skin)	3 days
MRSA Screen	3 days
Mycology Culture – CSF	42 days
Mycology Culture – routine specimen, not CSF	21 days
Peritoneal dialysate fluid C&S	6 days
Pertussis Culture	8 days
Sputum C&S	2 days
Stool C&S	3 days
Stool FBI	3 days
Throat, Nasal, Mouth C&S	2 days
Verotoxin	3 days
VRE Screen	3 days



### 3.0 SEROLOGY-PARASITOLOGY

Serological tests involve the detection and determination of antigen or antibody using a variety of laboratory procedures.

#### The Serology Section performs procedures for:

- Screening and diagnostic purposes - To detect acute or chronic infections due to viral, bacterial, fungal or parasitic agents.
- Immune Status Assessment - To detect past exposure or to evaluate response to immunization.
- Quantitative and qualitative testing of viral agents.
- Comprehensive Parasitology testing.
- Molecular-based typing and detection.
- Viral load and genotyping for patient management and surveillance purposes.
- Outbreak support.

#### Serology Programs

- Prenatal screening program includes Rubella, HBsAg, Syphilis testing. It is recommended that all expectant women be tested for HIV.
- Manitoba Bone Marrow Transplant Program for donors and recipients.
- Organ donor and eye bank donor screening.
- Dialysis Program consists of regular screening of all dialysis patients for hepatitis and HIV if requested.
- Sentinel surveillance for vector-borne diseases such as Western Equine Encephalitis.
- Outbreak support services.

#### Referred out Services

- The Serology Section acts as an intermediary for the referral of specimens to appropriate serology reference centres.

#### Emergency On-Call Service (Paging system)

Serology provides 24-hour call back services. Call 945-6655 and the on-site Security Guard will refer the call to the medical staff on call. Emergency Needlestick exposure testing at the approval of CPL Physician. Most after-hour STAT requests must be authorized through the 'on-call' physician. The only exception is organ donor transplant screening serology.

**Note:** This is a general description of services and not meant to be exclusive.

### 3.1 SEROLOGY TESTS

Serological tests involve the detection and quantification of specific antibody titers, using a variety of laboratory procedures. These procedures may be used for:

- **Screening or diagnostic** - To detect acute or chronic infections due to viral, bacterial, fungal or parasitic agents.
- **Immune Status Assessment** - To detect past exposure or to evaluate response to immunization.
- **Prenatal testing** includes HBsAg, Rubella, Syphilis and HIV (HIV opt-out available) tests. Additional tests may be requested (e.g., toxoplasmosis, Varicella Zoster, Parvovirus and HCV). Toxoplasmosis testing will only be performed if specifically requested and history of contact with cats is stated on the requisition.
- **Quantitative Testing of Viral Agents** - to monitor response to/or direct therapy.
- **Viral Genotyping** - as part of an epidemiologic investigation or patient management.
- **Hepatitis** - Indicate clinical condition or reason for requesting the test. HBsAg and HCV Ab will be tested on adults and HAV IgM on children and seniors.
- **Post-exposure protocol** - Follow instructions developed by Manitoba Public Health at: [www.gov.mb.ca/health/publichealth/cdc](http://www.gov.mb.ca/health/publichealth/cdc) or telephone 788-6722 during work hours, or 945-0183 after hours and on holidays/weekends.
- **Molecular** - Molecular biology techniques are used to meet the growing demands of the medical community to develop quantitative and qualitative assays to assist in patient management. NOTE: For tests not performed locally, samples will be referred to a reference laboratory. Referral institution is always indicated on the report.

#### Requisition requirements for Hepatitis B testing

There are several different protocols for appropriately testing for HBV, each depending on the clinical scenario. It is therefore, extremely important to include patient clinical information on the requisition when requesting HBV tests.

The following are some important facts to consider when ordering:

#### HBs Ag (HBV surface antigen)

This is a marker of current active HBV infection. It cannot differentiate between acute and chronic infection. Screen test of choice.

**HBe Ag (HBVe antigen)**

A marker of highly infectious active HBV infection. May be absent sometimes in certain mutated forms (core mutant)

**HBs Ab (HBV surface antibody)**

A marker of immunity to HBV. Cannot differentiate between immunity acquired from vaccine and natural infection.

**HBe Ab (HBVe antibody)**

Another marker of immunity, but not protective and not indicative of resolved infection.

**HBc IgM (HBV core IgM antibody)**

A marker of acute HBV infection. The first reliable marker to appear in HBV infection.

**HBc Ab (Total HBV core antibody)**

Another marker used to differentiate naturally acquired immunity from immunization acquired immunity.

**HBV DNA (quantitative or qualitative)**

Only to be used under special circumstances, a marker of active HBV replication.

LABORATORY TESTS	SPECIMEN REQUIRED	MIN. VOL. REQ. FOR TESTS (ml)	FREQUENCY OF TESTING	TURN-AROUND TIME IN WORKING DAYS
ADB	Serum	.1-.5	Weekly	6
ASOT	Serum	.1-.5	Weekly	6
CHLAMYDIA IgM Ab	Serum	.5-1	Weekly	6
CMV IgG	Serum	.1-.5	Weekly	6
CMV IgM	Serum	.1-.5	Weekly	6
EBV IgM & IgG	Serum	.1-.5	Weekly	6
HB Core IgM & IgG	Serum	.1-5	Daily	2
HBsAG	Serum	1-5	Daily	2
HCV GENOTYPING	Plasma	5	Batched	21
HCV-RNA	Plasma	1-2	Batched	7
Hbs Ab	Serum	1-5	Daily	2
HBV Viral Load	Plasma	5	Batched	7
HCV Ab	Serum	1-5	Daily	2
HCV Viral Load	Plasma	5	Batched	21
HEP A IgG	Serum	1-5	Daily	2
HEP A IGM	Serum	1-5	Daily	2
HIV 1/2 Ab	Serum	1-5	Daily	2
HIV Viral Load	Plasma (EDTA)	5	Batched	5
H. PYLORI Ab	Serum	1-5	Bi-weekly	5
HSV IgM & IgG	Serum	.1-.5	Weekly	6
HTLV 1/2	Serum	1-5	Weekly	6
LEGIONELLA Ab	Serum	.1-.5	Monthly	25
LYME Ab	Serum	.1-.5	Weekly	6
MEASLES IgM & IgG	Serum	.1-.5	Weekly	6
MUMPS IgM & IgG	Serum	.1-.5	Weekly	6
MYCOPLASMA IgM	Serum	.5-1	Weekly	6
PARVO IgM & IgG	Serum	.1-.5	Weekly	6
PNEUMOCOCCAL Ab (pre and post-vaccine sample required)	Serum	.5-1	Monthly	25
RUBELLA IgG	Serum	.5-1	Daily	2
RUBELLA IgM	Serum	.1-.5	Weekly	6
SYPHILIS (DFA)	Slide		As required	2
SYPHILIS (RPR)	Serum	.5-1	Daily	2
SYPHILIS (VDRL)	Serum or CSF	.5-1	Daily	2
TOXOPLASMA IgM & IgG	Serum	.5-1	Weekly	6
VZV IgM	Serum	.1-.5	Weekly	6
VZV IgG	Serum	.1-.5	Daily	2
WEST NILE Ab	Serum	1-2	Weekly (min.)	6
WNNAT	Plasma (EDTA)	5	Batched	6

**NOTE:** For STAT requests, please call Serology at 945-7582 or 945-7634.

During outbreaks or increased disease activity, testing may occur more frequently. The published testing frequencies represent baseline service.

### **3.2 SAMPLE REQUIREMENTS**

- Collect 5-10 ml of blood as early as possible after onset of illness. Tests requiring serum require 1 red top serum tube. For plasma collect 1 EDTA tube (purple topped tube). Refer to Serology test list on previous page. Acute and convalescent or paired sera are:
  - **Acute: 1-3 days after onset of illness.**
  - **Convalescent: 21 days after onset of illness.**

#### **Instructions for Viral Load (PCR) Testing and Genotyping**

Specimen type and patient history are required to determine the most suitable test to perform. Following the instructions will reduce specimen rejection and unnecessary phone calls.

##### Requisition Requirements

- Use the general requisition (MG696) for Hepatitis B/C Viral Load (quantitative)
- Use 'Retrovirus Nucleic Acid' testing requisition form (MG5126) for HIV Viral Load
- Proper requisition and specimen information are required. Indicate specimen collection date, time, and patient's code or name on the requisition.
- Patient history, i.e. initial assessment, follow up, patient on treatment etc. is required to perform the proper test.
- Please ensure physician number, facility number, and the name and address for report distribution is filled in.

##### Specimen Type

- 10 cc whole blood in EDTA tube (purple top tube) **or**
- EDTA plasma: label the tube "EDTA plasma", and include patient's name or code and requisition number.

##### Transport Requirements

- Deliver whole EDTA blood to CPL within 4 hours and no later than 4:00 p.m. during working days.
- Keep plasma refrigerated and deliver it within 24 hours, from collection time, on cold pack. For longer transportation time, store plasma frozen at -20°C and deliver frozen to CPL.

#### **Guidelines for Submitting Specimens for HIV1/2 or HTLV I/II Provirus Testing** Specimen Collection - Whole Blood Only (Anticoagulant Required)

- Collect blood in EDTA tubes (lavender). Heparin (green) tubes will not be accepted. It is not necessary to include a red top (no anticoagulant) tube for serological testing.



- A minimum volume of 2 mL is required; 5 mL is preferred.
- Blood should be kept at room temperature at all times.

#### Shipping

- Record date blood was drawn. Specimen should be received at CPL Monday-Wednesday AM only.
- Ship by courier directly to CPL to ensure receipt within 4 hours of collection.
- Ship at ambient temperature; do not freeze or cool.

For further information, please contact 945-3183 or 945-7545.

#### **Collection for syphilis DFA**

##### Lesions

- Remove any scab or crust covering the suspected chancre.
- Secondary infection exudate, if any, should be removed with a gauze sponge.
- If necessary, compress the base of the lesion or apply a suction cup to the lesion to promote the accumulation of tissue fluid on the ulcer surface.
- Apply a glass slide to the oozing lesion or use a sterile bacteriological loop to transfer the fluid from the lesion to the glass slide.
- Allow to air dry, do not apply any fixative. Label the slide.
- Place in slide holder and transport to CPL with a CPL requisition.

#### **Serologic Antibody Tests Requiring Acute And Convalescent Bloods or Paired Sera**

- Legionella

#### **Serologic Antibody Tests Performed On A Single Sample; A Second Sample May Be Requested**

- |                         |              |
|-------------------------|--------------|
| • Arbovirus             | • Mycoplasma |
| • <i>Chlamydia spp.</i> | • Parvovirus |
| • Lyme disease          |              |

### **Serologic Tests Generally Requiring Only A Single Sample**

- Cytomegalovirus (CMV) IgM or IgG
- Epstein Barr virus (EBV) Antibodies
- *H. pylori* Antibody
- Hepatitis A, B, C, or other hepatropic viruses
- Herpes simplex IgM or IgG
- HTLV - 1/2 Antibody
- Human immunodeficiency virus (HIV) Antibody
- Measles Antibody
- Mumps Antibody
- Parasitic serology
- Parvovirus IgM or IgG
- Rubella Antibody
- Syphilis Serology
- Toxoplasma Antibody
- Varicella zoster virus IgM or IgG

### **3.3 REQUISITIONS**

Use the general requisition (MG696) for all testing except non-nominal HIV. Fill out requisition completely, see Section 1.1.1 for requirements.

**HIV Ab requisition (MG13396)** must be completed for HIV antibody testing. Include complete code, epidemiologic data, physician billing number, facility number, phone number, and check off informed consent box. Results will be delayed if code and epidemiological data are not complete. Do not submit consent form to CPL.

**HIV viral load** - use the Retrovirus Nucleic Acid Testing Requisition (MG5126): see previous section *Instructions for Viral Load (PCR) Testing and Genotyping*.

### **3.4 TRANSPORT**

Ship serum or plasma in tightly capped polypropylene tubes. Place specimens individually in leak proof bags with requisition on the outside of the bag. **DO NOT USE STAPLES.**

### 3.5 REFERRED OUT SEROLOGY TESTS

Amoeba Ab	HIV-genotyping
Anti-delta Ab	HIV1/2 provirus
Anti-hepatitis E virus	HTLV I/II provirus
Arbovirus Ab	Human granulocytic ehrlichiosis
Aspergillus Ab	Hydatid/Echinococcus Ab
Babesia Ab	Leishmania Ab & PCR
Bartonella henselae Ab	Leptospira species Ab
Blastomycoses Ab	Lyme PCR
Botulinum Ab titre	Lymphocytic choriomeningitis Ab
Botulinum toxin	Lymphogranuloma venereum Ab
Brucella Ab	Malaria Ab
Chlamydia Ab	Meningococcus Ab
Coccidioides Ab	Pneumococcus Ab
Cryptococcus Ag & Ab	Rabies antibody (RFFIT)
Cysticercosis Ab	Rickettsia Ab
Diphtheria Ab	Schistosoma Ab
EBV early Ag	Strongyloides Ab
EBV nuclear Ag	Tetanus Ab
Endemic typhus Ab	Toxocara Ab
Escherichia	Trichinella Ab
Fascioliasis Ab	Trypanosoma Ab
Filaria Ab	Tularemia Ab
Haemophilus influenzae B Ab	Yellow fever Ab
Hantavirus Ab	Yersinia species Ab
Histoplasma Ab	

**Note:** This is not an exclusive list.

### 3.6 PARASITOLOGY TESTING

#### **Specimen Collection:**

- Reliable screens for enteric ova and parasites require 2-3 stool samples collected on different dates.
- Sample must be thoroughly emulsified at the time of collection. One part of stool into three parts SAF fixative.
- Samples not in fixative, or not in proper clinical containers, will not be processed (exception: *Cryptosporidium* speciation).

- Ensure patients are free of barium, cathartics or antibiotics as these substances may interfere with the examination.
- For information regarding other parasitology services offered at CPL, please call the Parasitology lab directly at 945-7825.
- For scabies skin scrapings: Place a single drop of mineral oil over unexcoriated burrow. Scrape lesion 6-7 times with a 15 scalpel blade until tiny specks of blood appear. The mineral oil will emulsify the scrapings. Transfer the emulsified scrapings with the blade to a clean glass slide and cover with a cover slip. Repeat several times. Package securely and forward to CPL expediently. CPL will do a microscopic exam to look for any stage or sex of mite, feces, eggs and egg casings.

**NOTE:** For *Cryptosporidium* speciation in outbreak cases, submitted stool must be preservative-free.

#### **Requisition:**

- Use the CPL General Requisition (MG696) for ordering parasite investigations.
- Ensure all relevant clinical information is given; i.e. symptoms, history of travel, etc. This will affect test selection at CPL. The more information the better.

<b>CAUSAL AGENT</b>	<b>SPECIMEN REQUIRED</b>	<b>TEST PERFORMED</b>
<i>Note: Processed daily.</i>		
<i>Acanthamoeba</i> species	Contact Parasitology laboratory	
<i>Ancylostoma duodenale</i> <i>Necator americanus</i>	Feces in SAF	Microscopy
<i>Angiostrongylus cantonensis</i>	Contact Parasitology laboratory	
Arthropods (mites, ticks, fleas, lice, fly maggots, etc.)	Dead: submit dry or in 70% alcohol Alive: submit with slightly moistened cotton	Microscopy, Gross ID
<i>Ascaris lumbricoides</i>	Feces in SAF	Microscopy
	Worm passed in feces Submit unpreserved in 0.85% NaCl, or if there is a delay in transit of three or more days, submit in 5% formalin or 70% alcohol.	Gross ID Serology

**CAUSAL AGENT****SPECIMEN REQUIRED****TEST PERFORMED***Note: Processed daily.*

<i>Babesia</i> species	Thick and thin blood films Blood with anticoagulant (EDTA)	Microscopy
<i>Balantidium coli</i>	Feces in SAF	Microscopy
<i>Blastocystis hominis</i>	Feces in SAF	Microscopy
<i>Clonorchis sinensis</i> (Chinese liver fluke) <i>Opisthorchis felineus</i> <i>Opisthorchis viverrini</i> <i>Metorchis conjunctus</i>	Feces in SAF	Microscopy
<i>Cryptosporidium</i> species	Feces in SAF	Microscopy
<i>Cyclospora cayetanensis</i>	Feces in SAF	Microscopy
Cysticercosis (pork tapeworm, larval stage)	Serum	Serology
<i>Demodex folliculorum</i> <i>Demodex brevis</i>	Skin scrapings including hair follicles and sebaceous glands Submit dry or mounted between two slides. Prior consultation is preferable.	Microscopy
<i>Dientamoeba fragilis</i>	Feces in SAF	Microscopy
<i>Diphyllbothrium</i> species (broad fish tapeworm)	Feces in SAF Worm segments Submit unreserved in 0.85% NaCl, or if there is a delay in transit of three or more days, submit in 5% formalin or 70% alcohol.	Microscopy
<i>Echinococcus granulosus</i> (dog tapeworm) <i>Echinococcus multilocularis</i>	Aspirated fluid from cyst Contact Parasitology laboratory Cyst, excised Serum	Ab detection
<i>Entamoeba histolytica</i>	Feces in SAF Serum	Microscopy Ab detection



CAUSAL AGENT	SPECIMEN REQUIRED	TEST PERFORMED
<i>Note: Processed daily.</i>		
<i>Enterobius vermicularis</i> (pinworm)	Pinworm paddle applied to perianal region  Vaseline paraffin anal swab or cellulose (transparent NOT translucent or opaque) tape preparations	Microscopy
<i>Fasciola gigantica</i> <i>Fasciola hepatica</i>	Feces in SAF	Microscopy
<i>Fasciolopsis buski</i>	Feces in SAF	Microscopy
<i>Giardia lamblia</i> (duodenales)	Feces in SAF Duodenal drainage	Microscopy
<i>Heterophyes heterophyes</i> <i>Metagonimus yokogawai</i>	Feces in SAF	Microscopy
<i>Leishmania tropica</i>	Smear from edge or base of lesions Contact Parasitology laboratory Serum	Micro (indicate)
<i>Leishmania brasiliensis</i> <i>Leishmania mexicana</i>	Smear from edge or base of lesions Contact Parasitology laboratory Serum	Micro (indicate)
<i>Leishmania donovani</i> (infantum)	Thick and thin blood films Contact Parasitology laboratory Biopsy material (spleen, liver, lymph nodes) Blood with anticoagulant (EDTA) Serum	Microscopy (indicate)  Ab detection
<i>Loa loa</i> (African eye worm)	Thick and thin blood films Blood with anticoagulant (EDTA)	Microscopy (indicate) Ab detection
Microsporidia	Feces in SAF	Microscopy (indicate)
Maggots	Maggots (clinical specimens only) Dead: submit dry or in 70% alcohol Alive: submit with slightly moistened cotton	Microscopy Gross ID
<i>Naegleria</i> species <i>Hartmannella</i> species <i>Acanthamoeba</i> species others	Contact Parasitology laboratory	
<i>Onchocerca volvulus</i> <i>Mansonella streptocerca</i>	Skin biopsy Contact Parasitology laboratory Aspirated material from skin nodules Excision of nodule	Microscopy (indicate)

CAUSAL AGENT	SPECIMEN REQUIRED	TEST PERFORMED
<i>Note: Processed daily.</i>		
<i>Paragonimus species</i> (lung fluke)	Feces in SAF Sputum	Microscopy Microscopy
<i>Pediculus humanus capitis</i> (head louse) <i>Pediculus humanus corporis</i> (body louse) <i>Phthirus pubis</i> (crab louse)	Adults, nymphs, or eggs ("nits") Submit dry or in 70% alcohol.  Infested hairs	Microscopy
<i>Plasmodium vivax</i> <i>Plasmodium malariae</i> <i>Plasmodium ovale</i> <i>Plasmodium falciparum</i> <i>Plasmodium knowlesi</i>	Thick and thin blood films from finger blood (at height of paroxysm and 8-16 hours later) Blood with anticoagulant (EDTA)	Microscopy (indicate)  Ab detection
<i>Sarcoptes scabiei</i>	Skin scrapings at end of tracks, fresh Submit dry or mounted mineral oil scrapings between two slides. Prior consultation is preferable.	Microscopy
<i>Schistosoma haematobium</i> (bladder blood fluke)	Urine	Submit mid-stream to terminal urine Microscopy Serology
<i>Schistosoma japonicum</i> (oriental blood fluke)	Feces in SAF	Microscopy Serology
<i>Schistosoma mansoni</i> <i>Schistosoma intercalatum</i>	Feces in SAF	Microscopy Serology
<i>Strongyloides stercoralis</i>	Feces in SAF Duodenal contents by intubation serum Sputum	Microscopy Microscopy Ab detection
<i>Taenia saginata</i> (beef tapeworm) <i>Taenia solium</i> (pork tapeworm)	Feces in SAF  Worm segments Submit unpreserved in 0.85% NaCl, or if there is a delay in transit of three or more days, submit in 5% formalin or 70% alcohol. Gross ID	Microscopy
<i>Toxoplasma gondii</i>	CSF Contact Parasitology laboratory Biopsy material Whole blood in anticoagulant	Microscopy  Serology Ab detection
<i>Trichinella spiralis</i>	Serum	Ab detection

CAUSAL AGENT	SPECIMEN REQUIRED	TEST PERFORMED
<i>Note: Processed daily.</i>		
<i>Trichostrongylus</i> species	Feces in SAF	Microscopy
<i>Trichuris trichiura</i> (human whipworm)	Feces in SAF	Microscopy
<i>Trypanosoma rhodesiense</i> <i>Trypanosoma gambiense</i>	Blood films, thick and thin Lymph aspirated from nodes CSF	Microscopy
<i>Trypanosoma cruzi</i>	Blood films, thick and thin Lymph aspirated from nodes CSF Serum	Microscopy Microscopy  Ab detection
<i>Wuchereria bancrofti</i> <i>Brugia malayi</i> <i>Mansonella perstans</i> <i>Mansonella ozzardi</i> <i>Loa loa</i>	Blood smear, thick and thin Blood with anticoagulant Aspiration from lymph vessels and nodes	Microscopy (indicated) Ab detection) Microscopy Serology

**Note: Not an exclusive test list.**

## 4.0 VIRUS DETECTION

Clinical virology services involve the isolation or detection and identification of human viral pathogens from clinical specimens using established procedures such as:

- Cell culture – many viruses are grown and identified in established cell lines.
- Rapid diagnostics – sensitive and specific procedures that provide accurate results within hours to aid in patient management.
  - Electron microscopy (EM)
  - Immunofluorescent antibody techniques
  - Latex agglutination
  - Enzyme linked immunosorbent assays (ELISA)
  - Molecular-based diagnostics
- Viral strain identification – subtyping for epidemiological and public health purposes, i.e., outbreak management, etc.
  - Neutralization
  - Immunofluorescence
  - Hemagglutination inhibition
  - Immunoelectron microscopy
  - Molecular-based typing
  - ELISA
  - Embryonic egg inoculation

### Emergency on-call service (paging system)

Call 945-6655 and the on-site Security Guard will refer the call to the medical staff on call. A technologist is always available for STAT testing requests.

### Other services

- Virology studies and projects generated from outside Manitoba Health.
- Education activities including participation in University of Manitoba post-graduate medical education programs.
- International Reference Centre for EM Virology.

### Transplant program support

Surveillance for, and diagnosis of viral infections common in immunocompromised patient populations and monitoring of response to antiviral therapy.

### Public Health program support

- Participation in various surveillance programs (national and international).
- Viral strain characterization.

- Meningoencephalitis exanthematous, respiratory and enteric outbreak investigations.
- Setting public health policy regarding viral disease.

#### Referral services

- Low volume or esoteric test requests are forwarded to reference laboratories.
- Level 3 and 4 pathogen investigations and prion investigations are forwarded to appropriate reference facilities.

**Note: This is a general description of services and not meant to be exclusive.**

### **4.1 SPECIMEN REQUIREMENTS**

- Most viruses do not survive well at room temperature.
- Feces should be sent fresh, unpreserved and unfrozen.
- EDTA blood must be kept at room temperature and must reach the laboratory within six hours of collection if it is being sent as whole blood. See 4.2.1 for further instructions.
- Other specimens for virus isolation should be refrigerated and transmitted to the laboratory as quickly as possible (e.g. cold pack). Swabs should be sent in virus transport medium (VTM), see Transport Media (Section 1.6). Ensure each specimen is properly numbered.

If in doubt, always consult with the laboratory before sending the specimen. A brief clinical history, DATE OF ONSET of symptoms, the DATE OF COLLECTION and the TYPE of specimen must accompany each specimen. Be sure to include presence of outbreak and a contact phone number.

### **4.2 SPECIMEN COLLECTION**

Specimens should be collected as early in the disease as possible. **Do not send dry swabs or swabs in bacterial transport medium.**


**4.2.1 Blood for quantitative CMV NAAT:** 5 ml purple stoppered EDTA tube (without separator). Must be kept at room temperature and must reach CPL within six hours of collection (during working hours) if it is being sent as whole blood. If the specimen cannot be forwarded to CPL within these parameters, the specimen should be spun at 800xg for 20 minutes and the plasma removed. The plasma should be kept cold until it is sent and transported with a cold pack.



- 4.2.2 Blood for BK virus:** Requires a minimum 5 ml purple stoppered EDTA tube (without separator).
- 4.2.3 Blood for quantitative EBV:** Requires a purple stoppered EDTA tube (without separator) taken and delivered to CPL on Monday or Wednesday before 12:00 p.m., or Tuesday before 4:30 p.m.
- 4.2.4 Bone marrow aspirate:** Collect in purple stoppered EDTA tube (without separator). Immediately invert several times to mix properly. Must be kept at room temperature and deliver to CPL within six hours.
- 4.2.5 Bronchial Alveolar Lavage (BAL):** 10 ml in a sterile container. If < 3 mls, add to a VTM bijou. Keep at 4°C.
- 4.2.6 CSF:** See Section 2.1.5, items 1-4 for collection instructions. A minimum of 0.5 ml is required for each test requested. Keep at 4°C.
- 4.2.7 Genital or mouth swabs for Herpes:** Use dacron or cotton-tipped swabs. Calcium alginate swabs not recommended. Swab the affected area, break off swab into VTM. Consult laboratory for rapid HSV testing. Keep at 4°C.
- 4.2.8 Lesion:** Expose and clean base of lesion with sterile gauze and saline. Scrape epithelial cells from base vigorously with a sterile swab. If dry, moisten swab in sterile saline, swab lesion, break off into VTM. Use a cotton or rayon swab, not calcium alginate. Keep at 4°C.
- 4.2.9 Lesion Smear for Molluscum contagiosum only:** Touch a glass microscope slide directly to an unroofed lesion. Air dry. Do not use cover slip or fixative. Place slide in slide-mailer and secure with an elastic band or clip. Keep at 4°C.
- 4.2.10 Nasal Swab:** Swab anterior nares as far back as possible. Break off into VTM. Keep at 4°C.
- 4.2.11 Nasopharyngeal aspirate:** Place a flexible plastic catheter gently into the posterior nasopharynx. Apply gentle suction with a syringe or wall suction, collect sample into a trap device, flush with 2.0 ml of VTM, then transfer to a sterile bijou bottle; **do not send the trap or tubing.** Keep at 4°C.

#### 4.2.12 Nasopharyngeal swab:

- a. Per nasal method: Remove any mucous from the patient's nose. Tilt the patient's head back slightly (about 70°) to straighten the passage from the front of the nose to the nasopharynx to make insertion of the swab easier. Gently insert the flocked swab into the nasopharynx (half the distance from the corner of the nose to the front of the ear). In adults, this distance is usually about 4 cm, in children this distance is less. Gentle rotation of the swab may be helpful. Rotate the swab several times to dislodge the columnar epithelial cells. Place swab(s) in VTM and break/cut shaft short enough to fit in bottle. Keep at 4°C.

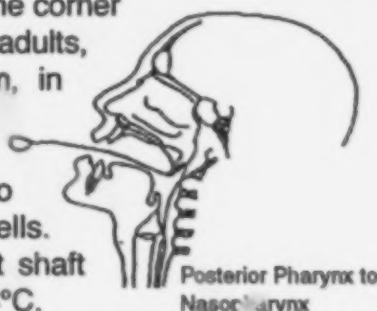

- b. Per ora method: Insert a flocked swab (if small-tipped, wire shaft used then bend shaft to give a slight curve) into the nasopharynx by passing the swab up behind the soft palate (see Figure). Vigorous swabbing will be more likely to collect the needed nasoepithelial cells. Place swab in VTM and break/cut shaft short enough to fit in bottle. Keep at 4°C.

**4.2.13 Rectal Swab:** If stool is unobtainable, a rectal swab may be submitted. Break off into VTM. Keep at 4°C.

**4.2.14 Stool (Raw):** Submit raw material in sterile container (no more than one-half full), without any preservatives or transport media. **DO NOT FREEZE.** Keep at 4°C.

**4.2.15 Throat Swab (Oropharyngeal):** Use dacron or rayon-tipped swab. Swab back of throat vigorously. Break off swab into VTM. Keep at 4°C.

**4.2.16 Tissue Biopsy:** A minimum specimen diameter of 2 mm is required. Tissue should be suspended in VTM for transport. Keep at 4°C.

**4.2.17 Tracheal Secretion:** Add specimen to VTM. Keep at 4°C.

**4.2.18 Urine:** Approximately 15 to 20 ml is required. Place in sterile container. Keep at 4°C.

**4.2.19 Vesicle Fluid:** Disinfect area with alcohol swabs (except if vesicle is located on mucous membrane.) **Remove 1.5 ml** of VTM from Bijou bottle. Fluid is collected by piercing the vesicle with a sterile needle attached to a tuberculin syringe, and aspirating as much material as possible. Rinse needle and syringe in the 0.5 ml VTM remaining in Bijou bottle. Keep at 4°C.

LABORATORY TESTS	SPECIMEN REQUIRED	TEST METHODS	FREQUENCY OF TESTING	TURN-AROUND TIME
Adenovirus	Respiratory, Throat NPA, Eye, Fecal Fecal	Tissue Culture	Daily	3-14 days
		Electron Microscopy	Daily	1-2 days
Coronavirus	Fecal	Electron Microscopy	Daily	1-2 days
Coxsackievirus	Respiratory, Fecal	Tissue Culture	Daily	3-14 days
CytomegaloVirus (CMV)	Urine, amniotic fluid, Respiratory, Biopsy, Bone Marrow Aspirate EDTA Blood CSF (on request)	Tissue Culture	Daily	7-21 days
		NAAT	Mon.-Wed.	1-5 days
Echovirus & Enterovirus	Respiratory, Fecal CSF (meningitis)	Tissue Culture	Daily	3-14 days
		NAAT	Referred out	2-7 days
Epstein Barr (EBV)	CSF EDTA blood Should be collected on Sun. or Mon. and received by noon Mon., or should be collected on Tues. or Wed. and received by Wed. noon	NAAT	Referred out Mon. & Wed. p.m.	5-14 days
Herpes simplex	Lesion swabs, vesicle fluid Lesion swab only CSF	Tissue Culture	Daily	1-8 days
		IFA NAAT	Referred out	STATS-3 hrs. 3-7 days
Herpes	Vesicle fluid	Electron Microscopy	Daily	1-2 days
		Tissue Culture		1-10 days
Human Herpesvirus 6, 7, 8	CSF	NAAT	Referred out	5-10 days
Influenza	Respiratory	Tissue Culture	Daily	3-10 days
		EIA	Seasonal	STAT-2 hrs.
Measles	Resiratory	Tissue Culture	Referred out	1-5 days
Molluscum contagiosum	Pustular Slide Swab in VTM	Electron Microscopy	Daily	1-2 days

LABORATORY TESTS	SPECIMEN REQUIRED	TEST METHODS	FREQUENCY OF TESTING	TURN-AROUND TIME
Mumps Virus	Respiratory	Tissue Culture NAAT	Daily	7-15 days
			Referred out	7-10 days
ORF	Pustular	Electron Microscopy	Daily	1-2 days
Papilloma Virus	Tissue, Biopsy	NAAT	Referred out	5-10 days
Polyoma (BK,JC)	Urine EDTA Blood	Electron Microscopy NAAT	Daily	1-2 days
			Referred out	
Parainfluenza	Respiratory	Tissue Culture	Daily	5-14 days
Poliovirus	Respiratory, Fecal Fecal only	Tissue Culture Electron Microscopy	Daily	3-14 days
				1-2 days
Poxvirus	Vesicle, Pustule (call lab before sending)	Electron Microscopy PCR	Referred out	1-2 days
Reovirus	Fecal	Tissue Culture Electron Microscopy	Daily	5-14 days
				1-2 days
Respiratory syncytial virus	NPA, NPS, Trach secretions	Tissue Culture EIA	Daily	3-14 days
			Seasonal testing	1-4 days
Rhinovirus	Respiratory	Tissue Culture	Daily	7-14 days
Rotavirus	Fecal	Electron Microscopy	Daily	24-48 hrs.
Rubella	Products of conception Urine	Tissue Culture	Daily	7-14 days
SARS	Respiratory, fecal (call MOH before sending)	Culture PCR	Referred out	2-7 days
Small Round Enteric Virus	Stool	Electron Microscopy	Daily	1-2 days
Varicella Zoster	Lesion swab	Culture IFA	Daily	6-10 days
			Thursday	1-5 days
				STAT - 3 hrs.

**NOTE: For STAT requests, please call Virology at 945-6858.  
Not an exclusive test list.**

## 5.0 NEWBORN SCREENING AND PUBLIC HEALTH CHEMISTRY

### Newborn Screening

The Newborn Screening Section of CPL screens all newborn babies in Manitoba for inherited disorders of metabolism and endocrine dysfunction. The screening program is guided by the Manitoba Perinatal Screening Committee. Screening is performed at the biochemical and genetic levels using dried blood spot specimens collected following birth.

The Manitoba Perinatal Screening Committee defines neonates as less than 28 days of age; therefore any unscreened children older than this (including refugee and immigrant children) should have plasma amino acid and TSH studies requested through HSC Clinical Chemistry DSM. Further laboratory studies relevant to ethnic background or country of origin should also be considered (e.g., sickle cell screen or hemoglobinopathy) when conducting screening for children over 28 days of age.

Babies are screened for inborn errors of carbohydrate metabolism (Galactosemia), fatty acid oxidation defects (Carnitine Palmitoyl Transferase-type 1 Deficiency), and amino acid metabolism defects (Glutaric Acidemia-type 1 and Phenylketonuria). Other screening detects disorders of endocrine function (Congenital Hypothyroidism and Congenital Adrenal Hyperplasia) and multiple carboxylase deficiency.

### Maternal Serum Screening

CPL, in collaboration with the Department of Human Genetics of the University of Manitoba, provides Maternal Serum Screening (MSS) to pregnant women in Manitoba as part of their prenatal care. This test provides an estimation of the risk for **fetal open neural tube defects**, Down syndrome, trisomy 18 and Smith-Lemli-Opitz syndrome (SLOS). CPL tests 4 biochemical markers (Quad Test) in the mother's blood that are produced by the fetus and placenta. The biochemical markers are alpha-fetoprotein (AFP), unconjugated estriol (uE3) and human chorionic gonadotropin (hCG), and Inhibin A (DiA).



### **Emergency or On-call**

By special arrangement only. Contact the Chief Technologist to arrange coverage.

### **Uninsured Services**

- Out-of-Province Neonatal Screening and Genetic Family studies.

## **5.1 NEWBORN SCREENING PROGRAM**

Cadham Provincial Laboratory supplies blood collection cards to birthing facilities, information pamphlets about the program to new parents, and program information to health care professionals upon request. Specimen collection instructions are given on the back of the blood collection card and in the newborn screening guidelines for health care providers. Call (204) 945-7458 for supplies.

### **Specimen collection**

**Collection:** Clean skin of heel with alcohol pledget, wipe dry with a sterile gauze pad. Puncture with disposable, 2.0 mm Lancet. If bleeding is slow, it is helpful to hold limb dependant for a short period of time before spotting blood on filter paper. Do not layer blood sample.

**Full Term Baby:** Take sample at time of hospital discharge, regardless of age. Babies sampled at less than 24 hours of age will require a repeat sample.

**Premature or ill Baby:** Take first specimen at five days of age and second specimen at two to three weeks of age or at time of hospital discharge, whichever comes first. Mark second specimen "repeat".

**Home Birth:** Take sample at two to three days of age.

**Blood Card Handling:** Fill **all** circles with blood, apply from **one** side only. Let blood soak through to the periphery of each circle. Allow to dry on a clean, dry surface at room temperature, *before* placing in the plastic cover (minimum 4 hours). Do not handle or contaminate blood spot area. Keep the card out of direct sunlight and away from heat sources while drying.

Deliver cards immediately to Cadham Provincial Laboratory.

All **requested information** must be supplied. Do not use regular postal services to mail the cards as this may cause unnecessary delays. Forward on the day of collection, if possible.

**Samples collected at <24 hours of age**

Infants whose first sample is taken at <24 hours of age, will be referred to the appropriate public health office for the collection of a second sample.

**Twins**

At 2 weeks of age, twins will be referred to the appropriate public health office for the collection of a second sample for a repeat TSH test. This is to rule out congenital hypothyroidism of one infant masked by the healthy twin on the first screen due to *in utero* twin-to-twin transfusion. Repeat samples will be requested from all same-gender twins. In cases of male and female twins (fraternal twins) the public health office will be contacted to verify the gender of the twins, but further follow-up will not be required.


**Results**

A negative report will be sent to the Medical Records Department of the birthing facility or to the practitioner if no facility is indicated. Immediate referral to a pediatric consultant is made in cases of significantly abnormal findings or critically elevated results. In cases of moderately elevated results, a request for a repeat sample is made to the infant's follow-up physician, midwife, Public Health Nurse, or the Nurse-In-Charge for newborns living in remote communities.



## Neonatal Screening: Blood Specimen Collection and Handling Procedure


1. Equipment: sterile lancet with tip less than 2.0 mm, sterile alcohol prep, sterile gauze pads, soft cloth, blood collection card, gloves.

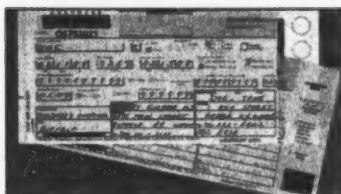

2. Complete ALL information. Do not contaminate filter paper circles by allowing the circles to come in contact with spillage or by touching before or after blood collection.



3. Hatched area (  ) indicates safe areas for puncture site.

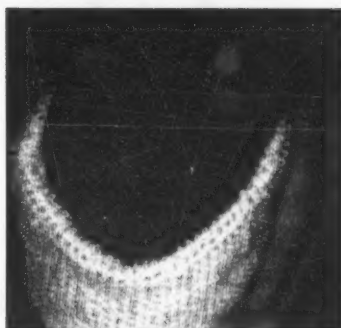

4. Warm site with soft cloth, moistened with warm water up to 41°C, for three to five minutes.




5. Cleanse site with alcohol prep. Wipe DRY with sterile gauze pad.


Reference: Schleicher & Schuell




6. Puncture heel. Wipe away first blood drop with sterile gauze pad. Allow another **LARGE** blood drop to form.

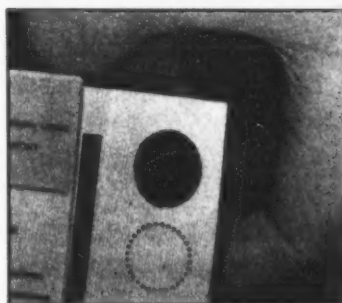

7. Lightly touch filter paper to **LARGE** blood drop. Allow blood to soak through and completely fill circle with **SINGLE** application to **LARGE** blood drop. (To enhance blood flow, **VERY GENTLE** intermittent pressure may be applied to area surrounding puncture site). Apply blood to one side of filter paper only.

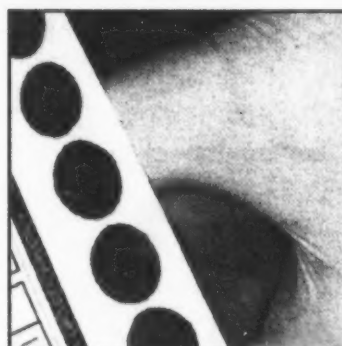


8. Fill remaining circles in the same manner as step 7, with successive blood drops. If blood flow is diminished, repeat steps 5 through 7. Care of skin puncture site should be consistent with your institution's procedures.

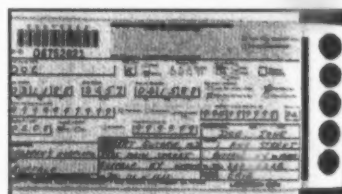

9. Dry blood spots on a dry, clean, flat non-absorbent surface for a minimum of four hours. (**Note:** This is not a current sample of the CPL form).



10. Deliver completed card for testing to CPL within 24 hours of collection if possible.

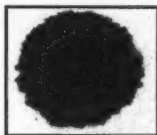


VALID SPECIMEN


## Simple Spot Check

Allow a sufficient quantity of blood to soak through to completely fill the preprinted circle on the filter paper. Fill all required circles with blood. Do not layer successive drops of blood or apply blood more than once in the same collection circle. Avoid touching or smearing spots.


### Poor Specimens - Possible Causes:

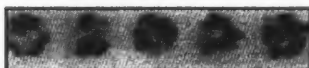

1. Specimen quantity insufficient for testing.

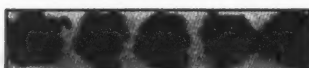


2. Specimen appears scratched or abraded.

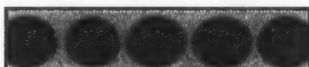

3. Specimen not dry before mailing.



4. Specimen appears supersaturated.

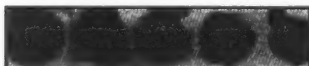

5. Specimen appears diluted, discolored or contaminated.


- Removing filter paper before blood has completely filled circle or before blood has soaked through to second side.
- Applying blood to filter paper with a capillary tube.
- Allowing filter paper to come in contact with gloved or ungloved hands or substances such as hand lotion or powder, either before or after blood specimen collection.
- Applying blood with a capillary tube or other device.
- Mailing specimen before drying for a minimum of two hours.
- Applying excess blood to filter paper, usually with a device.
- Applying blood to both sides of filter paper.
- Squeezing or "milking" area surrounding the puncture site.
- Allowing filter paper to come in contact with gloved or ungloved hands or substances such as alcohol, formula, antiseptic solutions, water, hand lotion or powder, etc., either before or after blood specimen collection.
- Exposing blood spots to direct heat.





6. Specimen exhibits serum rings.

- Not wiping alcohol from puncture site before making skin puncture.
- Allowing filter paper to come in contact with alcohol, hand lotion, etc.
- Squeezing area surrounding puncture site excessively.
- Drying specimen improperly.
- Applying blood to filter paper with a capillary tube.

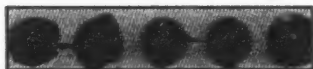

7. Specimen appears clotted or layered.

- Touching the same circle on filter paper to blood drop several times.
- Filling circle on both sides of filter paper.

## MANITOBA NEWBORN SCREENING PROGRAM

Condition	Congenital Hypothyroidism (CH)	Phenylketonuria (PKU)	Galactosemia	Multiple carboxylase deficiency	Congenital adrenal hyperplasia (CAH)	Glutaric Acidemia Type 1 (GA-1)	Carnitine palmitoyltransferase Type 1 (CPT-1)
<b>Test</b>	Thyroid stimulating hormone (TSH)	Phenylalanine	Galactose and galactose-1-phosphate	Biotinidase	17-OH-progesterone (17-OHP)	Allele Specific PCR for the Island Lake mutation	PCR amplification restriction enzyme digestion
<b>Follow-up Test</b>			Beutler spot test uridyl transferase	Semi-quantitative biotinidase		N/A	N/A
<b>Reference Range</b>	<30 mIU/L at <=72 hours of age; <20 at >72 hours of age	<200 $\mu$ mol/L	<400 $\mu$ mol/L	O.D.>0.06	<25 nmol/L	<b>Not</b> homozygous for the Island Lake	<b>Not</b> homozygous for mutation
<b>Action: Moderately elevated results</b>						N/A	N/A
<b>Action: Critically elevated results</b>							
<b>Time tested</b>	Twice weekly: Monday & Thursday	Daily: Monday - Friday	Daily: Monday - Friday	Twice weekly: Tuesday & Friday	Daily: Monday - Friday	Bi-weekly	Bi-weekly

## 5.2 MATERNAL SERUM SCREENING (QUAD TESTING) PROGRAM

There are two components to Maternal/Serum Screening. The first is alpha-fetoprotein (AFP) screening, and the second component includes unconjugated estriol (uE3), dimeric inhibin A (DIA) and human chorionic gonadotropin (hCG) with the AFP. This combination is known as the quad test or four-marker.

### Specimen requirements

A minimum of 0.5 ml of serum is required at 16 - 18 weeks of gestation for testing. Collect a minimum of 5.0 ml of blood in a serum separator tube; a full 10 ml serum separator tube preferred (minimum half-full required). Within two hours of collection, centrifuge to separate the serum from cells (Plasma is not suitable). Send centrifuged primary tube, (no aliquots) properly labeled to CPL as soon as possible, **DO NOT SEND ALIQUOTS IN PLASTIC TUBES**. Delay in sending sample should be avoided. If delay is unavoidable, do not freeze. Store the centrifuged specimen at 4°C until shipment. Samples seven or more days in transit may be compromised for analytes tested.

### Requisition requirements

Please fill in all clinical information on the requisition as completely and accurately as possible. This is essential to ensure correct calculation of the risk. Turn around time is significantly improved when all required information is included on the requisition. The following patient information is **absolutely required** for accurate interpretation:

1. Name of patient
2. PHIN
3. Date of birth
4. Gestational age:
  - Ultrasound date and measurements (BPD, CRL or composite gestational age), (most accurate for interpretation).
  - Last menstrual period (LMP) date.
  - Expected Date of Confinement (EDC).

### Results

Samples obtained between 16-18 weeks of gestation are optimal. However, samples received between 15 weeks and 20 weeks 6 days can be interpreted. A limited MS AFP interpretation can be made up to 23 weeks 6 days. Follow-up of abnormal results is the responsibility of the Human Genetics Department of the University of Manitoba, CSB - FE211, 820 Sherbrook Street, Winnipeg, Manitoba; telephone 787-2097 or 787-2098.

## 6.0 INFORMATION MANAGEMENT

The Information Management section is responsible for the provision of service in five main areas for microbiology, serology and virology:

- Data entry services
- Patient inquiry services
- Results reporting (paper/electronic/fax)
- Data requests
- Request for data retrieval

### **Data Entry Services**

- Data entry of all incoming requisitions (excluding Newborn Screening and Maternal Screening) will match patient demographics with the Client Registry to ensure the integrity of the CPL database in matching test requests and maintaining a comprehensive patient profile.
- All requests for laboratory testing must be accompanied by a completed CPL requisition. The integrity of the data entered into the CPL databases is compromised by incomplete or illegible requisition information, and may result in rejection or delayed reporting. Requisition entry requires extensive searching on the Client Registry to ensure correct patient matching.
- **Equally important is the clarity and completeness of the "return report to" portion of the requisition. This will ensure prompt reporting to the ordering practitioner and if requested on requisition, a copy will be forwarded to another practitioner. If no practitioner or address is provided, the sample will be rejected.**
- Request to amend requisition information - changes to requisition information will require the completion of this form (see section 8.0).
- Data entry and verification of laboratory results.
- Ensuring reportable results are flagged and reported to CDC in Manitoba and Nunavut.
- Creating and maintaining outbreak codes and their respective reporting, as requested.

### **Patient Inquiry Services**

- **Patient Inquiry - 945-6611**

Hours of operation:      Monday - Friday 0800 - 1630 hrs.  
   Saturday 0900 - 1600 hrs.

When making a telephone inquiry for results you will be asked to provide the following information to ensure the authenticity of the requester and also to ensure correct patient and/or results are being given. A log of all patient inquiry calls is maintained at CPL.

- a. Your name, telephone number and institution you are calling from.
- b. The patient PHIN is the most efficient way of searching patient results from both the online system and/or the archived tape files.

When a PHIN is not readily available, you will be asked for the patient name, gender and date of birth. As many people have aliases or change names (married, etc.) this information may be required to verify the correct patient file.

- c. The clerk will then verify with the caller the corresponding demographics and pertinent requisition information before providing results. If the request is for HIV results, or if any type of interpretation or explanation of results is required, the call will be transferred to the appropriate section.
- d. Results (verbal and/or hard copy) may be provided with a valid request. **Reporting to multiple practitioners remains the responsibility of the ordering practitioner as the CPL system is limited to only one "copy report to" address per requisition.**
- e. The requester must be verified before confidential information is provided by phone. Results will be released only to physicians, midwives, public health practitioners and/or their designates.



- f. In some cases it may be deemed necessary to have the requested results telephoned back to the physician and/or facility (certain results, results requiring interpretation or to validate the requester, preliminary results, etc).
- g. Calls made to Patient Inquiry (204) 945-6611 after regular hours will receive a voice message and an alternate number to call for **emergent** requests - (204) 945-6655. CPL Security will record the information required and relay it to the physician on call. The computer system is available to the physician on call every day between 0600 and 2200 hours.

#### **Results Reporting Provincial Data Network (PDN), Paper and Fax**

- Based on the return address information, result reports will be produced either on paper or electronically on the PDN. The electronic report files are available for downloading by 0900 hours, the day following results entry and is accessible 24 hours a day, 7 days a week. Clients have access to the most current 3 weeks of laboratory results and 30 days on the PDN. Private access to the PDN reporting is available to all requesters who have a computer, communication software and a printer. Inquiries for information or access may be made to the Information Co-ordinator at (204) 945-6845. \*The PDN reporting will be discontinued as of May 1, 2010.
- Some paper reports printed at CPL are mailed or couriered out. Laboratory results will be auto-faxed to a confirmed secure fax.

#### **Requests for Data Retrieval**

- Co-ordinate and respond to all incoming requests for the retrieval of mainframe data back to 1982 including statistical requests, patient profiles, reprinting of previous reports, etc.

Requesters may be asked to complete a CPL Data Request Form according to the CPL Data Request Guidelines (see below).

A charge for service will be applicable to agencies outside Manitoba Health. Once the Data Request is received an estimate of cost will be provided.

#### **Data Request Guidelines**

1. Requests for CPL data (internal and external), will be made using the CPL Data Request Form (see Section "Forms and Requisitions" and forwarded to the Information Co-ordinator at 750 William Avenue (fax # 786-4770) or by e-mail – [laura.girden@gov.mb.ca](mailto:laura.girden@gov.mb.ca).

2. Once a request is assessed for priority and an estimated project time (internal) is given, the requestor will be notified with the estimated date of completion. Please provide as much lead time as possible.

**Note:** Research project requests for data may require ethics approval from an appropriate institutional review board (i.e., Faculty Committee on the Use of Human Subjects in Research, U. of M., etc.) prior to approval at CPL.

3. All relevant correspondence must accompany the CPL Data Request form. In addition, all external research projects will require approval from the Health Information Privacy Committee (HIPC) if identifiable personal health information is requested.
4. The data report is sent to the requestor with the specialist's name(s) and contact number. If data includes personal health identifiers, Manitoba Health Information Systems Branch will receive the data and provide it to the requester. Any report generated using CPL data must be copied to the CPL specialist for review prior to publication to ensure appropriate interpretation.
5. The data request procedure will be reviewed at 6 months and then annually to ensure appropriateness.

#### **Requests by Individuals or Their Representatives for Personal Health Information**

- Requests to CPL by individuals or legal representatives for personal health information must be in writing, except in an emergency (involving an immediate threat to the mental/physical health or safety of an individual the information is regarding). Whenever possible, the report will be forwarded to the appropriate attending physician.
- A standard request for information form will be used (see section 8.0). The form will be submitted to the Privacy Officer.
- CPL must respond within thirty (30) days to the requester.
- A fee may be charged for acquisition of personal health information. A charge-back policy consistent with the fee regulated under *The Freedom of Information and Protection of Privacy Act* will apply (see section 8.0).

## 7.0 ALPHABETICAL INDEX OF TESTING

### SECTIONS:

CM = Clinical Microbiology  
NBS/PHC = Newborn Screening &  
Public Health Chemistry  
PA = Parasitology

SE = Serology

VD = Virus Detection

Disease or Syndrome (Causal Agent(s))	Specimen Requirements	Tests / Examinations	Section
Actinomycosis ( <i>Actinomyces israelii</i> )	Pus (preferably with granules) Bronchial washing Intrauterine device	Microscopy and special culture	CM

NOTES: Organism takes 3 or more days to grow in cultures.  
If suspected, *Actinomyces* culture must be specifically requested.

Adenovirus infections (Adenoviruses types 1- 41) upper respiratory tract; pneumonia; acute respiratory disease syndrome (ARD); Infections of conjunctiva	Throat swab in VTM Conjunctival swab Lung aspirate or biopsies in VTM.	Viral culture	VD
Gastroenteritis	Feces	EIA (Types 40/41), EM, Culture	VD
<i>Aeromonas hydrophila</i> group	Stool or rectal swab	Culture	CM

NOTES: Causal agent of diarrhea.

AIDS  (HIV Virus)	Clotted blood or serum  EDTA blood	Serology - EIA WB confirmatory Viral load Genotyping	SE
-------------------------	--	--	----

NOTES: HIV viral load increases up to 1000 copies/ml after several below detection results have been reported by several laboratories. These increases are not typically indicative of the development of drug resistance. Follow-up specimen in 4 weeks may help in resolving this issue.  
Viral load and genotyping only done on confirmed positives. Use special HIV antibody and viral load requisitions only. See section 3.2.

Amebiasis (see  
Dysentery, amebic,  
amebic encephalitis  
amebic hepatitis)

Disease or Syndrome (Causal Agent(s))	Specimen Requirements	Tests / Examinations	Section
Amebic encephalitis ( <i>Hartmannella</i> , <i>Naegleria</i> , <i>Acanthamoeba</i> species)	CSF	Microscopy for trophozoites	PA
NOTES: DO NOT REFRIGERATE CSF specimens. Infection acquired from water by swimming or bathing, even in chlorinated pools. Suspicion of amebic disease must be indicated on the requisition.			
Amebic hepatitis ( <i>Entamoeba histolytica</i> )	Clotted blood or Serum	Serology	PA
Anaerobic infections (see also Gas Gangrene)	Pus from deep abscesses, brain, lung or pelvic region or body cavities in TM. See also pus. Lung aspirates or biopsies.	Microscopy and culture	CM
NOTES: Many anaerobic strains are slow to grow and reports cannot be expected for 4 days or more. Many body sites have anaerobic normal flora. Culture of such sites is unprofitable (eg. skin, mouth, throat, sputum, vagina and bowels).			
Ancylostomiasis (see Hookworm Disease)			
Anthrax ( <i>Bacillus anthracis</i> )	Isolate  Swab or pustular fluid from skin lesion in TM Sputum in rare instances of pulmonary infection Blood for culture	Identification, toxin testing Microscopy and special culture  Culture	CM
NOTES: Use gloves and mask for collection. If lesion is dry, moisten swab in sterile water, saline or broth and rotate beneath the edge of the eschar. Notify the laboratory when anthrax is suspected and mark the requisition clearly "suspect anthrax." Consider notification of regional MOH.			
Arbovirus infections (Arboviruses)	Clotted blood or serum	Serology (Referred out)	SE
NOTES: Consult with laboratory, Virology Section, for submission of brain biopsies. Requests for other than West Nile Virus are referred out to the National Microbiology Laboratory.			

Disease or Syndrome (Causal Agent(s))	Specimen Requirements	Tests / Examinations	Section
Ascariasis ( <i>Ascaris lumbricoides</i> )	Feces in SAF Worm (passed in feces or vomit)	Microscopy for ova Identification of worm	PA
NOTES: Worms may be submitted in sterile urine container with or without formalin or alcohol preservative. Infected persons must be treated.			
Aseptic meningitis (see Meningitis, viral)			
Aspergillosis ( <i>Aspergillus fumigatus</i> <i>Aspergillus species</i> )	First morning respiratory secretion Biopsy material Pus	Microscopy and fungal culture	CM
	Clotted blood or serum	Serology (Referred out)	SE
NOTES: May be a contaminant in sputum and isolation does not necessarily indicate infection. Please indicate on requisition form if aspergillosis is suspected.			
Atypical mycobacteria (see Mycobacteria)			
Atypical pneumonia (see <i>Mycoplasma</i> infections; also Pneumonia, viral and other non-bacterial)			
Balanitis (Various bacteria and yeast)	Swab in TM	Microscopy and culture	CM
Balantidiasis ( <i>Balantidium coli</i> )	Feces in SAF Scrapings of ulcerated bowel (sigmoidoscopic)	Microscopy for trophozoites and cysts	PA
Bilharziasis (see Schistosomiasis)			
Biotinidase deficiency	Newborn screening card	Qualitative spot test and semi-quantitative spectrophotometry	NBS/ PHC



Disease or Syndrome (Causal Agent(s))	Specimen Requirements	Tests / Examinations	Section
Blastomycosis (North American) ( <i>Blastomyces dermatitidis</i> )	Scrapings from skin lesions Purulent exudate from base of skin lesion Aspirated material from abscesses Sputum Cerebrospinal fluid	Microscopy and fungal culture	CM
	Serum	Serology (Referred out)	SE

NOTES: A dimorphic fungus found in the Kenora area, and in Southeast Manitoba. Consult with the laboratory and mark requisition clearly "blastomycosis suspected."

---

Bornholm's Disease  
(see Coxsackievirus infections)

---

Botulism ( <i>Clostridium botulinum</i> types A, B, and E)	Isolate Feces, tissue exudate in TM, gastric contents unfixed	Confirmation Culture for <i>C. botulinum</i> neurotoxin detection in blood, food, vomit or gastric contents (Referred out)	CM
	Clotted blood or Serum (3 X 10 ml)		SE

NOTES: Consult laboratory ASAP! Collect in sterile screw-capped jars. Transport immediately. Take blood early in illness (30 ml if possible) BEFORE giving antitoxin. Very rarely *C. botulinum* causes wound infections.

---

Bronchiolitis and viral respiratory disease (Respiratory syncytial virus; adenovirus, parainfluenza, influenza)	NPA, NPS in VTM	Viral culture EIA	VD
--	-----------------	----------------------	----

---

NOTES: Throat swabs will not be rejected, but sensitivity reduced.  
Pack specimens for virus isolation in cold packs and send by fastest possible means.  
DO NOT FREEZE (or RSV will not survive).

---

Disease or Syndrome (Causal Agent(s))	Specimen Requirements	Tests / Examinations	Section
Brucellosis ( <i>Brucella abortus</i> , <i>B. melitensis</i> , <i>B. suis</i> )	Isolate Blood for culture	Identification and Typing Culture	CM
	Clotted blood or Serum	Serology (Referred out)	SE
NOTES: Multiple cultures are recommended; they are given prolonged incubation before being reported as negative. Interpretation of serologic findings is difficult in some chronic infections, consult the laboratory.			
Campylobacter infections ( <i>Campylobacter jejuni</i> , <i>C. coli</i> )	Feces Rectal swab in TM	Culture	CM
NOTES: Major cause of sporadic gastroenteritis, occurring in summer and early fall. Highest incidence in infants and young children, associated with ingestion of contaminated milk and water or improperly handled or cooked food - primarily poultry products.			
Candidiasis ( <i>Candida albicans</i> <i>Candida spp.</i> )	Mouth, throat, cervical vaginal swabs, urethral in TM Skin and nail scrapings Urine	Microscopy and culture	CM
NOTES: If pulmonary involvement suspected, submit transtracheal or lung aspirates. Refrigerate samples if delay in transport.			
Carditis (see Coxsackie virus infections)			
Catscratch fever ( <i>Bartonella henselae</i> )	Lymph node biopsy	Special pathology (Referred out) Special culture PCR	CM
	Tissue - unfixed		
	Clotted blood or Serum	Serology (Referred out)	SE
NOTES: Tissue to be delivered on ice ASAP. Please phone lab.			

Disease or Syndrome (Causal Agent(s))	Specimen Requirements	Tests / Examinations	Section
Cercopithecine herpes virus (simian herpes virus, Herpes B virus)	Wound swab in VTM Clotted blood or serum, 3-4 ml CSF, biopsy or necropsy tissue in VTM	Viral culture (Referred out) PCR (Referred out)	VD

NOTE: In all cases, contact the lab before collecting or submitting specimens.

Cervicitis (see  
Gonorrhoeae, Chlamydial  
Infections)

Chagas disease ( <i>Trypanosoma cruzi</i> )	Blood films (thick and thin; unstained) Lymph aspirated from nodes or chagoma	Microscopy PCR (Referred out)	PA
	Clotted blood or Serum	Serology (Referred out)	SE

NOTES: Endemic occurs in Central and South America.

Chancroid (Soft chancere) ( <i>Haemophilus ducreyi</i> )	Swab of pus or scrapings from lesions Dacron swab in 2SP CTM	Special culture  PCR (Referred out to National Microbiology Lab)	CM
--	---	--	----

NOTES: Culture requires special medium and is available only by special request. Please phone lab: 945-7204. There is no optimal transport medium available, but Amies charcoal TM may be used. Transport to CPL immediately. Requests for molecular testing must include clinical background leading to suspicion of *H. ducreyi* infection in a patient, including any additional factors that may increase the probability of *H. ducreyi* infection. Swab samples for molecular testing require dacron swab and 2SP CTM. Phone lab at 945-7204 to obtain 2SP CTM.

Chickenpox (Varicella) Shingles (Zoster) (Varicella-Zoster virus)	Vesicle fluid Base of lesion swabbed vigorously In VTM	Electron microscopy and Viral culture Rapid test - DFA	VD
	Clotted blood or serum or plasma	Serology	SE

NOTES: Immune status: detection of IgG.  
Diagnosis: detection of specific IgM, presence indicates recent infection.

Disease or Syndrome (Causal Agent(s))	Specimen Requirements	Tests / Examinations	Section
Chlamydia infections: respiratory <i>Chlamydia psittaci</i> , <i>C. pneumoniae</i> Ornithosis Psittacosis Pneumonitis TWAR	Sputum or nasopharyngeal aspirate preferred Nasopharyngeal swab in 2SP Chlamydia TM is acceptable	PCR (Referred out)	CM

NOTES: *C. pneumoniae* is tested for IgM on only acute sample.

Chlamydia infections: STI and others <i>Chlamydia trachomatis</i> Trachoma Inclusion conjunctivitis N.G.U. Infantile pneumonitis P.I.D. Prepubertal vaginitis	Cervical swab Urethral swab Urine-first void (20-30 ml) (see 2.2) Tracheal secretions Nasopharyngeal secretions Rectal swab (anal columns) Throat swab Conjunctival swab	NAAT (GenProbe Aptima)  DFA (Microtrak)  NAAT (genProbe Aptima) DFA (Microtrak)	CM
	Serum	IgM in neonatal infections Serology (Referred out)	SE
Lymphogranuloma venereum (LGV) (see Lymphogranuloma Venereum (LGV))			

NOTES: For NAAT testing by GenProbe Aptima, only cervix swabs, urethral swabs and urine are acceptable. If eye swabs are submitted for NAAT, results will be reported as "for investigational purposes only." Vaginal swabs are not appropriate. Please use the swabs provided in the kit, and place only the blue swab in the tube. Do not discard the liquid preservative in the tube. Positive GenProbe samples are retained for 3 weeks in event further testing is required.  
The presence of IgM-antibody is diagnostic.

Disease or Syndrome (Causal Agent(s))	Specimen Requirements	Tests / Examinations	Section
Cholera ( <i>Vibrio cholerae</i> including the El Tor biotype)	Isolate Feces	Typing (Referred out) Microscopy and culture	CM
NOTES: Specify on requisition if cholera is suspected.			
Cholinesterase deficiency (acquired)	Clotted blood or serum	Qualitative	NBS/ PHC
Chromoblastomycosis ( <i>Phialophora species</i> , <i>Cladosporium carrionii</i> )	Scales from skin lesions Pus Sputum Biopsy	Microscopy and fungal culture	CM
<i>Clostridium difficile</i> (antibiotic-associated diarrhea, pseudomembranous colitis) ( <i>Clostridium difficile</i> toxin)	Feces (10 ml)	Cytotoxin testing Special culture  Rapid toxin test	CM
NOTES: Request <i>C. difficile</i> testing. Inappropriate specimens include swabs, stool in transport medium or fixative, and formed stools. <i>C. difficile</i> culture from stool is by special request only.			
Colitis (See <i>Clostridium difficile</i> )			
CJD (Creutzfeldt-Jakob disease)	CSF (1 ml)	Immune blot (Referred out)	VD
NOTES: Please contact Lab before sending.			
Clonorchiasis ( <i>Opisthorchis</i> ) ( <i>Clonorchis sinensis</i> ), the Chinese liver fluke <i>Metorchis conjunctus</i> (Canadian liver fluke)	Feces in SAF	Microscopy for ova	PA
NOTES: Occurs in the Far East. May persist twenty years or more. In North America, mainly in Aborigines. It is difficult to differentiate the ova of these and related species and when in doubt they are reported as 'opisthorchid' ova.			



Disease or Syndrome (Causal Agent(s))	Specimen Requirements	Tests / Examinations	Section
Clostridial infections (see Anaerobic Infections and Gas Gangrene)			
CMV (see Cytomegalovirus)			
Coccidiomycosis ( <i>Coccidioides immitis</i> )	Sputum Aspirated material Pleural fluid Biopsy material Cerebrospinal fluid	Microscopy and fungal culture	CM
	Clotted blood or Serum	Serology (Referred out)	SE
NOTES: Infection is generally contracted in arid areas of North America, especially California, therefore a travel history is helpful. Mark requisition clearly "coccidioides suspected."			
Common cold or minor respiratory illness (rhinovirus)	NPA, NPS in VTM Throat swab in VTM Nasopharyngeal aspirate in VTM	Viral culture	VD
NOTES: Laboratory confirmation rarely necessary. Testing only recommended as part of an outbreak investigation.			
Congenital adrenal hyperplasia (CAH)	Newborn screening card	Fluoroimmunoassay followed by extraction of free 17-OHP	NBS/ PHC
Congenital primary hypothyroidism	Newborn screening card	Fluoroimmunoassay for TSH	NBS/ PHC
Congenital infections (Various including rubella, herpes, cytomegalovirus, toxoplasma, Listeria, Chlamydia, Group B <i>Streptococcus</i> )	EDTA plasma Body fluids (urine, CSF, amniotic fluid, etc.) Tissues (biopsies)	NAAT Viral culture Culture for bacteria	VD CM
	Clotted blood or serum	Serology	SE
NOTES: See under the individual causative agents.			

Disease or Syndrome (Causal Agent(s))	Specimen Requirements	Tests / Examinations	Section
Conjunctivitis (Many bacteria including Staphylococcus, Haemophilus spp., Streptococcus pneumoniae, Moraxella spp, Neisseria gonorrhoeae, Chlamydia trachomatis, and several viruses, including Herpes and Adenovirus)	Conjunctival swab in Amies charcoal TM Swab for chlamydia  Swab in VTM	Microscopy and culture  NAAT (Gen Probe Aptima) DFA (Microtrak)  Viral culture	CM    VD

NOTES: Use Amies Charcoal transport medium for bacterial causal agents.

Use NAAT or DFA collection kits for chlamydia.

If specimen submitted for NAAT, results will be reported as "for  
investigational purposes only."

Swabs to be cultured for viruses must be sent in VTM. Take specimens  
before using topical anesthetics which may be anti-microbial.

Serology for adenoviruses, Herpes.

---

Contagious pustular  
dermatitis (see Poxvirus Infections)

---

Coryza (see Common colds)

---

Cowpox (see Poxvirus Infections)

---

Coxsackievirus infections (including aseptic meningitis, pleurodynia (Bornholm disease), febrile illness often with rash; Hand, foot and mouth disease; myocarditis, pericarditis, etc.) Herpangina (Coxsackievirus A, types 1-24 Coxsackievirus B, types 1-6)	Feces Cerebrospinal fluid, vesicle fluid, throat swab in VTM	Viral culture NAAT	VD
---	---	-----------------------	----

NOTES: Specimens may be mailed unfrozen as these viruses are relatively resistant.  
Each virus type requires a separate test.

---

Disease or Syndrome (Causal Agent(s))	Specimen Requirements	Tests / Examinations	Section
Creutzfeldt-Jakob disease (See CJD)			
Croup (Parainfluenza virus Influenza virus, Respiratory syncytial virus (RSV) and <i>Bordetella pertussis</i> )	NPS, NPA or Throat swabs in VTM	Viral culture EIA (where applicable)	VD
	NPS, NPA or throat swab in Amies charcoal medium	Bacterial culture	CM
Cryptococcosis ( <i>Cryptococcus neoformans</i> )	CSF Sputum Aspirated material from abscesses Biopsy material (ulcers, lymph nodes)	Microscopy and culture  Antigen detection (Referred out)	CM
	Serum	Serology (Referred out)	SE
Cryptosporidium Speciation	Feces in SAF Feces without Preservative	Microscopy (Referred out)	PA
NOTES: Cryptosporidium must be specifically requested on the requisition. Upon approval, speciation will be referred out for public health purposes.			
<i>Cyclospora cayentanensis</i>	Feces in SAF	Microscopy	PA
Cysticercosis ( <i>Taenia solium</i> )	Clotted blood or serum	Serology (Referred out)	SE
Cystitis (see Urinary Tract Infections)			

Disease or Syndrome (Causal Agent(s))	Specimen Requirements	Tests / Examinations	Section
Cytomegalovirus infections (Cytomegalovirus)	Urine Throat wash, biopsy Autopsy material in VTM Bone marrow aspirate EDTA blood* *specimen must reach laboratory same day or plasma removed and sent to lab.	Viral culture PCR	VD
	Clotted blood, serum or plasma	Serology	SE

NOTES: Send urine in sterile bottle, URGENTLY. (Do not freeze, 4°C is optimal). In urine, CMV is very labile at 25°C.  
It is also slow growing, culture may take 20 days or more.  
Prolonged excretion of virus in urine and saliva may occur in both congenital and acquired infection.  
Blood may be collected at any time during the illness.

---

Dengue (see  
Arbovirus infections)

---

Dermatophytosis ( <i>Epidermophyton floccosum</i> <i>Microsporum spp.</i> <i>Trichophyton equinum</i> , <i>T. mentagrophytes</i> , <i>T. rubrum</i> , <i>T. tonsurans</i> , <i>T. verrucosum</i> , <i>T. violaceum</i> )	Skin scrapings Nail clippings Hair	Microscopy and fungal culture	CM
---	--	----------------------------------	----

NOTES: Send complete affected hair.  
Wrap specimen in black paper.  
DO NOT USE slides, cotton, wool or tubes.  
Results of fungal culture may not be available until 2 to 3 weeks after receipt of specimen.  
Never use transport medium for skin scrapings.

---

Disease or Syndrome (Causal Agent(s))	Specimen Requirements	Tests / Examinations	Section
Diarrhea, bacterial (see also Cholera; *Food poisoning; Paratyphoid fevers; Typhoid fever) (Numerous bacterial species including <i>Salmonella</i> , <i>Shigella</i> <i>Escherichia coli</i> (enteropathogenic and verotoxin producing strains), <i>Yersinia</i> <i>enterocolitica</i> , <i>Campylobacter</i> , <i>Plesiomonas</i> <i>shigelloides</i> , <i>Aeromonas spp.</i> and <i>Vibrio spp.</i>	Feces - 30 ml/10g. Fill container minimum 1/3 full and no more than 1/2 full. No fixative.	Culture <i>E. Coli</i> Verotoxin, <i>C. Difficile</i> Toxin  FBI Investigation, Outbreak Investigation	CM
	Clotted blood or serum	Serology (Referred out for <i>Yersinia</i> )	SE

NOTES: Send material containing blood or mucus if these are present.  
Feces is always preferable to a rectal swab but if fecal specimen unobtainable, send rectal swab in TM.  
Enteropathogenic *E. coli* are reported only in feces of children less than 3 years of age except on special request.  
Acute diarrhea requires one stool. No transport media or preservative.  
Chronic diarrhea: submit stool on three separate days.  
Outbreaks should be noted on the requisition. Contact the regional MOH to establish outbreak or FBI status and obtain outbreak code for clinical specimens.

Diarrhea, Infantile (see also Diarrhea, viral) (Specific serotypes of <i>Escherichia coli</i> , <i>Campylobacter spp.</i> <i>Aeromonas spp.</i> )	Feces - fill container minimum 1/3 full and no more than 1/2 full. No fixative.	Culture Verotoxin testing	CM
---	--	------------------------------	----

NOTES: Commonly due to certain specific serotypes of *Escherichia coli* termed enteropathogenic serotypes, eg. 0111, 055, 0119, etc.  
Screening for these organisms is done on all fecal specimens from children under 3 years of age.



Disease or Syndrome (Causal Agent(s))	Specimen Requirements	Tests / Examinations	Section
Diarrhea, parasitic (See also Giardiasis, <i>Cryptosporidium</i> , <i>Cyclospora</i> , Balantidiasis, Ascariasis, Hookworm disease, Dysentery, amebic Diphyllobothriasis Microsporidiosis Teniasis, Trichostrongyloidiasis, Trichinosis, Trichuriasis, Worm infections)	Feces in SAF Segments of worm in feces	Microscopy for ova and parasites Antigen detection test	PA
	Clotted blood or serum	Serology (Referred out)	SE

**NOTES:** Multiple stool specimens are usually required to adequately rule out parasitic diarrhea. Three specimens collected on separate days is recommended.

Diarrhea, viral (Rota, small round enteric virus Enteroviruses (echo- and coxsackieviruses)), Enteric Adeno (Type 40-41), reoviruses	Feces (Fill container no more than 1/2 full). No fixative.	Viral culture EM EIA	VD
--	--	----------------------------	----

**NOTES:** When Rotavirus is suspected in infants, electron microscopy should be requested. This virus is probably the most common causal agent of winter vomiting and diarrhea in infants under 3 years of age. Fecal swab is a suboptimal specimen. Outbreaks should be noted on requisition. Contact the regional MOH to establish outbreak or FBI status and obtain code for clinical specimens.

Diphyllobothrium ( <i>Diphyllobothrium</i> species (fish tape worm))	Feces, in SAF Segments of worm in feces	Microscopy for ova and segments	PA
--	---	------------------------------------	----

**NOTES:** Treatment of infected individuals is recommended.

Disease or Syndrome (Causal Agent(s))	Specimen Requirements	Tests / Examinations	Section
Diphtheria ( <i>Corynebacterium diphtheriae</i> )	Isolate  Throat swab NPS, Ear swab Swab from skin lesion in TM	Confirmation and Toxigenicity Culture	CM
	Clotted blood or serum	Serology (Referred out)	SE

NOTES: If symptoms are suggestive of *C. diphtheriae* infection, please indicate on requisition "suspect *C. diphtheriae*". Toxigenicity tests are performed on all *C. diphtheriae* isolated.

Duchenne Muscular Dystrophy	Newborn screening card		NBS/ PHC
Duodenal and gastric ulceration	Gastric Biopsy Clotted blood or serum	Culture Serology	CM SE

NOTES: Submit biopsy material in a wide-mouthed, screw-capped, sterile container with sufficient sterile saline to keep moist. Transport ASAP. If a delay is anticipated, transport with ice packs.  
For serology, please note on requisition if patient is on treatment. Pre- and 6-12 months post-bloods will be tested together.

Dysentery, amebic ( <i>Entamoeba histolytica</i> )	Feces in SAF	Microscopy for trophozoites and cysts	PA
	Clotted blood or serum	Serology (Referred out)	SE

NOTES: IHA serology relevant only for extraintestinal amebiasis.

Dysentery, bacillary (see also Shigellosis) ( <i>Shigella</i> spp.)	Isolate Feces	Speciation/typing Culture	CM
---	------------------	------------------------------	----

NOTES: Submit feces or rectal swab. Send material containing blood or mucus if these are present. Outbreak should be noted on the requisition. Contact the Regional MOH to establish outbreak or FBI status and obtain outbreak code for clinical specimens. Refrigerate if delay in transport.

Disease or Syndrome (Causal Agent(s))	Specimen Requirements	Tests / Examinations	Section
Ear infections Otitis externa Acute otitis media, otomycosis (Several species of bacteria and fungi)	Ear swab in TM Pus aspirated through intact eardrum	Microscopy and culture	CM
Respiratory viruses	Ear swab in VTM	Culture	VD
NOTES: Place swab in transport medium. Mixed cultures generally occur; only pathogenic bacteria or fungi are reported unless a special request is made.			
Eastern equine encephalitis (see Arbovirus infections)			
EBV (Epstein Barr Virus) See Infectious Mononucleosis			
Echinococcosis ( <i>Echinococcus</i> <i>granulosus</i> , <i>E.</i> <i>multilocularis</i> )	Serum	Serology (Referred out)	PA
NOTES: Sensitivity of test will vary depending on size, integrity and location of cyst.			
Echovirus infections Aseptic meningitis, Rash, Diarrhea, Upper respiratory infection (Echoviruses, types 1-34)	Throat swab in VTM Feces CSF	Viral culture NAAT	VD
NOTES: See Coxsackie infections.			
Ectoparasites ( <i>Arthropods</i> ) See also Scabies	Parasite Hair with nits	Identification Microscopy and/or visual	PA
NOTES: Send in alcohol, if possible.			

Disease or Syndrome (Causal Agent(s))	Specimen Requirements	Tests / Examinations	Section
Encephalitis, viral, including epidemic, sporadic, and post- infectious types (Many viruses including arbo-, herpes- myxo-, paramyxo-, entero-, and poxviruses)	Throat swab in VTM Biopsy material (brain) Autopsy material in VTM Feces CSF	Viral culture   NAAT	VD
NOTES: Biopsy material of brain may be sent in suspected herpes encephalitis. Consult with the laboratory before submitting brain biopsies. Never place the brain biopsy in formalin.			
Enteric fever (see Typhoid fever)			
Entamoeba histolytica (see Dysentery, amebic, Amebic encephalitis, Amebic hepatitis)			
Enterobiasis ( <i>Enterobius vermicularis</i> )	Clear sticky tape applied to peri-anal region	Microscopy for ova	PA
NOTES: Pinworm ova are collected by pressing sticky side of clear cellophane (Scotch) tape against the peri-anal skin first thing after waking. The tape is then pressed onto a glass slide and sent for microscopic examination. Repeat exams may be necessary.			
Enterococcus colonization Enterococcus (VRE) (See below)	Swabs from rectum or ostomy, wounds, open skin lesions and/or line or device sites in TM Isolate	VRE Screen  Typing	CM
NOTES: VRE or suspected VRE isolates may also be forwarded to CPL for genetic typing (Van A, B, C), and/or pulsed-field gel electrophoresis (PFGE).			

Disease or Syndrome (Causal Agent(s))	Specimen Requirements	Tests / Examinations	Section
Enterococcus infection ( <i>Enterococcus faecalis</i> , <i>E. faecium</i> , and other species) (See above)	Swab from wound in TM Infected body fluid Urine	Culture	CM

Enterocolitis (*Yersinia*)  
(see Diarrhea,  
bacterial and  
*Yersinia* Infections)

Enteropathogenic <i>E. coli</i> See also diarrhea, bacterial; diarrhea, infantile	Feces Isolate	Culture Typing and toxin	CM
---	------------------	-----------------------------	----

NOTES: These organisms are sought in all feces specimens from children less than 3 years of age, and cases of traveller's diarrhea.  
Specimens from other cases are examined for these sero-types on special request only.

Enterovirus infections (see also Coxsackievirus infections; Echovirus infections; Poliomyelitis)	Throat swab in VTM Feces CSF	Culture NAAT	VD
--	------------------------------------	-----------------	----

Epidemic  
keratoconjunctivitis  
(see Adenovirus  
infections and  
Conjunctivitis)

Epidemic myalgia or  
pleurodynia (see  
Coxsackie virus  
infections)

Epiglottitis, acute ( <i>Haemophilus</i> <i>influenzae</i> )	Nasopharyngeal swab in TM Throat swab in TM	Culture	CM
--	---	---------	----

NOTES: Rapidly progressive, often fatal disease.  
Take respiratory tract specimens only after intubation.



Disease or Syndrome (Causal Agent(s))	Specimen Requirements	Tests / Examinations	Section
Epstein Barr Virus (EBV) (See Infectious Mononucleosis)			
Equine encephalitis, Eastern, Western or Venezuelan forms (see Arbovirus infections)			
Erysipelas (see Streptococcal infections)			
Erysipeloid ( <i>Erysipelothrix</i> <i>rhusiopathiae</i> )	Inject saline into lesion and re-aspirate Biopsy	Culture	CM
NOTES: History of animal or fish contact is usual in this cutaneous disease, usually of the hand and fingers. Clean and disinfect skin before sampling.			
Erythema infectiosum (5th disease) (Parvovirus B <sub>19</sub> )	Clotted blood, serum or plasma	Serology IgM & IgG	SE
ESBL (Extended- spectrum Beta lactamase ( <i>E.coli</i> , <i>Klebsiella</i> spp.)	Isolate Screening Rectal swab	ESBL confirmation Culture	CM
NOTES: Screen swabs are to be requested only on consultation with the CPL Director.			
Farmer's Lung ( <i>thermophilic fungi</i> )	Clotted blood or serum	Serology (Referred out)	SE
NOTES: The fungi, from mouldy hay, etc., act as allergens rather than infective agents.			
Fasciola ( <i>Fasciola gigantica</i> <i>F. hepatica</i> )	Feces in SAF	Microscopy	PA
Favus (see Dermatophytosis)			

Disease or Syndrome (Causal Agent(s))	Specimen Requirements	Tests / Examinations	Section
Fetal neural defects, trisomy 21 or 18 and SLOS (Maternal serum screening, AFP) (see also Triple Test)	Clotted blood or serum between 15-18 weeks gestation Amniotic fluid (AFP only)	Chemiluminescent immunoassay	NBS/ PHC
NOTES: Submit with the fully completed maternal serum screening requisition.			
Filariasis (see also Loiasis; Onchocerciasis) ( <i>Wuchereria bancrofti</i> )	Blood smear, thick or thin or Blood with anti- coagulant Clotted blood or serum	Microscopy for microfilariae  Serology (Referred out)	PA  SE
NOTES: Some specimens are periodic so optimal collection times may vary. Consult with Parasitology at CPL at 945-7825.			
Food poisoning, acute bacterial, viral and toxic forms (see also Botulism) ( <i>Staphylococcus aureus</i> , <i>Salmonella spp.</i> , <i>Clostridium perfringens</i> , <i>Listeria spp.</i> , <i>Shigella spp.</i> , <i>Vibrio spp.</i> <i>Bacillus</i> <i>cereus</i> <i>Yersinia spp.</i> , Verotoxin Small round and other enteric viruses)	Feces - fill container Minimum 1/3 full and no more than 1/2 full. No fixative Isolate  Clotted blood or serum	Culture, tests for toxin where indicated  Typing as indicated  Serology (Referred out)  Culture EM	CM    SE  VD
NOTES: Outbreaks should be noted on the requisition. Contact Regional MOH to establish outbreak or FBI status and obtain outbreak code for clinical specimens.			
Fungal infections (see Dermatophytosis, and individual fungal infections)			
Galactosemia and galactosemia variants	Newborn screening card	Spectrophotometry	NBS/ PHC

Disease or Syndrome (Causal Agent(s))	Specimen Requirements	Tests / Examinations	Section
Gas Gangrene ( <i>Clostridium</i> <i>perfringens</i> , <i>C.</i> <i>septicum</i> , <i>C. novyi</i> and other species)	Swabs from lesions in TM Necrotic tissue	Microscopy and culture	CM
NOTES: Anaerobic streptococci may cause similar lesions.			
Gastritis, peptic ulcers (see Duodenal and Gastric Ulceration)			
Gastroenteritis (see Diarrhea - bacterial, infantile, parasitic and viral. Food poisoning)			
German measles (see Rubella)			
Giardia ( <i>Giardia lamblia</i> )	Feces in SAF	Microscopy for trophozoites and cysts Antigen detection test	PA
NOTES: Examination of duodenal aspirates may be helpful.			
Gonorrhea ( <i>Neisseria</i> <i>gonorrhoeae</i> )	GenProbe Aptima swab of urethra, cervix, prepubertal vagina, conjunctiva Urine (first void 20-30 ml) (see section 2.2) Swab of throat, rectum, ovaries and fallopian tubes, vagina in charcoal TM Joint aspirate Isolate	NAAT  NAAT C&S  Typing and sensitivity	CM
NOTES: For GenProbe Aptima unisex swab collection kits, immerse blue swab in the provided preservative. Do not discard liquid or swab. Positive GenProbe samples are retained for 3 weeks in the event further testing is required. Conjunctival swabs in charcoal transport media for culture are preferred, if transport time is less than 48 hours. An air dried smear for chlamydia should accompany any NAAT specimens from pre-pubertal children (one streak approx. 1.5 inches or 3 cm. long on a clean glass slide). If eye swab submitted for NAAT, results will be reported as "for investigational purposes only."			

Disease or Syndrome (Causal Agent(s))	Specimen Requirements	Tests / Examinations	Section
Group B Streptococcus agalactiae (prenatal screen) (see Streptococcal infections)			
Guillain-Barre syndrome (Echo, Coxsackie, EBV, WNV, <i>Campylobacter</i> )	Throat swab in VTM Feces CSF	Viral culture  NAAT	VD
	Feces	Culture	CM
Hand-foot-and- mouth disease (Coxsackie group A viruses, especially types 16 and 9)	Vesicle fluid or swabs Feces Throat swabs in VTM	Viral culture	VD
Hantavirus (Muert canyon virus) (Sin Nombre virus)	Clotted blood or serum Tissue	Serology Referred out to Federal lab	SE VD
<i>Helicobacter pylori</i> (see Duodenal and gastric ulcerations)			
Hemolytic uremic syndrome (HUS) (Verotoxin producing <i>Escherichia coli</i> , <i>Shigella</i> )	Feces - fill container minimum 1/3 full no more than 1/2 full.	Direct fecal Verotoxin test (FVT) VT from Colony Sweeps (VT/PECS)	CM
NOTES: Send refrigerated stool as soon as possible <u>without</u> transport medium or fixative.			
Hepatitis A (Infectious hepatitis) (Hepatitis A virus)	Clotted blood Serum or (EDTA) plasma	Serology	SE
NOTES: Anti-HAV IgM is present in patients with acute hepatitis A infections. Presence of Anti-HAV IgG indicates immunity.			

Disease or Syndrome (Causal Agent(s))	Specimen Requirements	Tests / Examinations	Section
Hepatitis B (Serum hepatitis) (Hepatitis B Virus)	Clotted blood Serum or (EDTA) plasma	Serology  Viral load	SE
NOTES: HBsAg is present in patients with acute or chronic hepatitis. For interpretation of other tests, contact the laboratory.			
Hepatitis C (Hepatitis C Virus)	Clotted blood, serum or EDTA plasma	HCV antibody, HCV, RNA, RIBA, Viral load Genotyping	SE
NOTES: Hepatitis C accounts for a large proportion of cases of what was previously known as Non A-Non B Hepatitis.			
Hepatitis D, E	Clotted blood, serum or EDTA plasma	Serology (Referred out)	SE
Herpangina (Coxsackie A viruses)	Swab from lesions Throat swabs in VTM Feces	Viral culture	VD
	Paired clotted blood	Serology	SE
NOTES: Infrequently isolated in tissue culture.			
Herpes B Virus ( <i>simian herpes</i> ) (see Cercopithecine herpes virus)			
Herpes simplex virus infections (including herpes encephalitis, neonatal herpes, eczema herpeticum, genital herpes) (Herpes virus type I and II)	Vesicle fluid Base of lesion swabbed vigorously in VTM Throat swab, Biopsy material, Autopsy material or Urethral swab in VTM Cerebrospinal fluid	Electron microscopy examination of vesicular fluid Viral culture DFA  NAAT	VD
	Clotted blood and serum or plasma	Serology	SE
NOTES: For herpes encephalitis, brain biopsy may be sent after consultation with CPL Virology Section. A positive serum HSV IgM test is indicative of a recent infection.			

Disease or Syndrome (Causal Agent(s))	Specimen Requirements	Tests / Examinations	Section
Herpes zoster (see Chickenpox)			
Heterophyiasis ( <i>Heterophyes</i> <i>heterophyes</i> <i>Metagonimus</i> <i>yokogawai</i> <i>Opisthorchis</i> )	Feces in SAF	Microscopic examination for ova	PA
NOTES: Occurs in Middle and Far East and Southern Europe.			
Histoplasmosis ( <i>Histoplasma</i> <i>capsulatum</i> )	Bronchial washings Sputum Swab or scrapings from ulcer Biopsy material (lymph nodes, marrow) Cerebrospinal fluid	Microscopy and fungal culture	CM
	Clotted blood serum	Serology (Referred out)	SE
NOTES: <i>Histoplasma capsulatum</i> may be present in soil contaminated by birds and bats. Consult with the Lab & clearly mark requisition "Histoplasma suspected".			
Hookworm disease, ancylostomiasis (see also Trichostrongyliasis) ( <i>Ancylostoma</i> <i>duodenale</i> , <i>Necator americanus</i> , <i>Trichostrongylus</i> <i>species</i> )	Feces in SAF	Microscopy examination for ova and larvae	PA
Human herpes type 6 (HHV6) Roseola	Clotted blood or serum	Serology (Referred out)	SE
Hydatid Disease (see Echinococcosis)			
Hymenolepiasis ( <i>Hymenolepis nana</i> <i>H. diminuta</i> )	Feces in SAF	Microscopy examination for ova	PA



Disease or Syndrome (Causal Agent(s))	Specimen Requirements	Tests / Examinations	Section
Impetigo ( <i>Streptococcus pyogenes</i> , <i>Staphylococcus aureus</i> )	Swabs from lesions in charcoal TM	Culture	CM
	Clotted blood or serum	Serology (streptococcal)	SE
NOTES: Nephritis is occasionally associated with impetigo. Culture is most indicated during outbreaks.			
Infantile diarrhea (see Diarrhea, Infantile)			
Infectious hepatitis (see Hepatitis, A, B, C, D, E)			
Infectious mononucleosis (Epstein Barr virus, EBV)	Clotted blood or serum	Serology	SE
	Should be collected on Sun. or Mon. and be received by noon Mon. or should be collected Tues. or Wed. and received by Wed. noon	NAAT (Referred out)	VD
Influenza (Influenza viruses types A and B)	Nasopharyngeal swab or aspirate in VTM Throat swab or washing in VTM Autopsy material (lung) in VTM	Viral culture EIA NAAT	VD
NOTES: Outbreaks should be noted on the requisition. Contact the regional MOH to establish outbreak and obtain code for clinical specimens.			
Kala Azar (see Leishmaniasis, visceral form)			
Keratoconjunctivitis, viral (see Adenovirus infections; Herpes simplex infection)			
Keratomycosis (Many fungi)	Corneal scrapings	Fungal culture	CM
NOTES: Consult the Mycology Section before sending specimens.			

Disease or Syndrome (Causal Agent(s))	Specimen Requirements	Tests / Examinations	Section
Laryngitis, bacterial and acute laryngotracheo-bronchitis (croup) ( <i>Corynebacterium diphtheriae</i> , <i>Haemophilus influenzae</i> , <i>Streptococcus pyogenes</i> and other organisms)	Throat swab in TM Aspirated respiratory secretion	Culture	CM

NOTES: The majority of cases are caused by viruses.

Laryngitis, viral and acute laryngotracheo-bronchitis (croup) (Several viruses, including adeno-, parainfluenza, measles, respiratory syncytial, influenza, rhino-, and echoviruses)	NPA, NPS Throat swab in VTM	Viral culture EIA	VD
--	--------------------------------	----------------------	----

NOTES: Send specimens for viral isolation on a cold pack. DO NOT FREEZE.

Legionnaires' disease ( <i>Legionella pneumophila</i> )	Sputum Lung biopsy Bronchoscopy Specimens Tracheal secretions or aspirates Isolate Urine	Culture DFA    Typing Antigen detection	CM
	Paired clotted blood or serum 21 days apart	Serology	SE

NOTES: May mimic viral pneumonia. DFA not done on sputum. Transport specimens in sterile dry containers. Add a small amount of sterile non-bacteriostatic, distilled water to prevent desiccation if necessary. Do not add saline due to its inhibitory effect. Refrigerate in transit.

Disease or Syndrome (Causal Agent(s))	Specimen Requirements	Tests / Examinations	Section
Leishmaniasis, cutaneous form ( <i>Leishmania tropica</i> )	Biopsy from edge or base of lesions in skin Smear from edge of base of lesion	Microscopy Culture/PCR (Referred out)	PA
	Clotted blood or serum	Serology (Referred out)	SE
Leishmaniasis, mucocutaneous form (Espundia) ( <i>L. brasiliensis</i> , <i>L. mexicana</i> )	Skin biopsy	Culture/PCR (Referred out)	PA
	Smear from edge or base of lesion	Microscopy	
	Clotted blood or serum	Serology (Referred out)	SE
NOTES: Occurs in Mexico, Central and South America.			
Leishmaniasis, visceral form (Kala Azar) ( <i>L. donovani</i> )	Bone marrow films	Microscopy	PA
	Biopsy material (spleen, liver, lymph nodes)	PCR (Referred out) Culture (Referred out)	SE
	Clotted blood or serum	Serology (Referred out)	SE
Leprosy ( <i>Mycobacterium leprae</i> )	Biopsy of tissue affected, usually skin nodes Nasal scrapings	Microscopy	CM
NOTES: <i>M. leprae</i> cannot usually be cultured in vitro. Diagnosis is essentially a clinical one, supported by demonstration of acid-fast bacilli in the specimen. Referred to HSC TB lab. For further information, call HSC TB lab at 787-7652.			
Leptospirosis ( <i>Leptospira ictero- haemorrhagiae</i> , <i>L. canicola</i> , <i>L. pomona</i> and others)	Blood, Urine Autopsy material (liver, kidneys), CSF	PCR (Referred out)	CM
	Clotted blood or serum	Serology (Referred out)	SE
NOTES: Consult the laboratory before sending for PCR. Most infections are diagnosed serologically.			

Disease or Syndrome (Causal Agent(s))	Specimen Requirements	Tests / Examinations	Section
Listeriosis ( <i>Listeria monocytogenes</i> )	Blood, CSF, Vaginal swab, Amniotic fluid, placenta Isolate	Microscopy and culture  Typing	CM
NOTES: May cause meningitis, or granulomatous disease in the newborn or fetal death. Please indicate "possible Listeriosis" on requisition.			
Loiasis ( <i>Loa loa</i> )	Thick and thin blood films, Blood with anti-coagulant	Microscopy examination for microfilariae	PA
Lung fluke disease (see Paragonimiasis)			
Lyme Disease ( <i>Borrelia burgdorferi</i> )	Biopsy of tissue, CSF  Clotted blood or serum	Molecular testing (Referred out)	CM  SE
NOTES: Molecular testing by prior arrangement only. Serology may be negative during first stage erythema chronicum migrans.			
Lymphocytic choriomeningitis (LCM) (Lymphocytic choriomeningitis virus)	Blood (early), Urine Cerebrospinal fluid (late)  Serum	Referred out   Referred out	VD   SE
NOTES: Provide patient history and onset of illness			
Lymphogranuloma venereum (LGV) (Chlamydia trachomatis Serovars L1, L2, L3) (see also Chlamydia Infections)	Dacron swab of: Bubo Anogenital ulcers (rectal, vaginal, urethral) If no ulcers – cervical, urethral, rectal swabs Fluid aspirate  Serum	PCR (Referred out)     Referred out	CM     SE
NOTES: Place swabs in 2SP Chlamydia TM and transport as soon as possible. Also submit GenProbe Aptima urine, urethral or cervical swab for routine Chlamydia testing (see: Chlamydia Infections). Routine Chlamydial testing will not specifically confirm LGV as is not a substitute for the required specimens outlined above. Relevant clinical information is required before testing will be performed. Consult CPL at 945-7184.			
Maduromycosis ( <i>Madura foot</i> ) (see Mycetoma)			

Disease or Syndrome (Causal Agent(s))	Specimen Requirements	Tests / Examinations	Section
Malaria ( <i>Plasmodium vivax</i> <i>P. malariae</i> , <i>P. ovale</i> , <i>P. falciparum</i> ) <i>P. Knowlesi</i>	Thick and thin blood films on separate slides Clotted blood or serum	Microscopy examination for parasites Referred out (only if microscopy is negative)	PA  SE
NOTES: Take films at 6-18 hour intervals for 3 days. Air-dry blood films without heat.			
Measles, including diseases associated with the measles virus; giant cell pneumonia; encephalitis; subacute sclerosing panen- cephalitis (SSPE) (see also Panencephalitis) (Measles virus)	Throat swab in VTM Autopsy material (lung, brain) in VTM Cerebrospinal fluid Urine  Clotted blood or serum	Culture (Referred out)     Serology	VD     SE
NOTES: Presence of IgM is diagnostically significant. Consult the CPL Virology section prior to sending culture specimens.			
Melioidosis ( <i>Burkholderia</i> <i>pseudomallei</i> )	Sputum Swab of abscesses in TM	Culture	CM
NOTES: Endemic in South-East Asia and Northern Australia. Travel history is helpful as it is found in tropical and sub-tropical areas worldwide.			
Meningitis, bacterial (see also Meningococcal infections) ( <i>Neisseria</i> <i>meningitidis</i> , <i>Haemophilus</i> <i>influenzae</i> , <i>Streptococcus</i> <i>pneumoniae</i> , <i>Listeria</i> <i>monocytogenes</i> , and in the newborn, coliform organisms and Group B streptococci)	CSF, Skin lesions Swab in TM Isolate	Microscopy and culture   Typing	CM
NOTES: Send specimens by most rapid means of transport. If delay in transport is anticipated, send a smear and up to 3 ml of CSF in pediatric blood culture bottle. Transport ASAP.			

Disease or Syndrome (Causal Agent(s))	Specimen Requirements	Tests / Examinations	Section
Meningitis, viral ( <i>Aseptic meningitis</i> ) coxsackie A and B, echo, mumps, measles, Less common causes are: lymphocytic choriomeningitis, poliomyelitis, herpes simplex, EBV, varicella-zoster, rubella, arboviruses, influenza, adenoviruses	Throat swab in VTM Autopsy material in VTM (brain, spinal cord, intestinal contents) Feces Cerebrospinal fluid	Viral culture NAAT	VD
NOTES: Consult with the CPL Virology Section prior to collection of brain biopsy. Clinical and epidemiological history must be provided.			
Meningococcal infections, including meningitis and meningococcemia ( <i>Neisseria meningitidis</i> )	CSF Blood culture Swabs from petechial lesions in TM Isolate	Microscopy and culture PCR (consult with laboratory) (Referred out) Typing	CM
NOTES: Send specimens by most rapid means of transport. If delay in transport is anticipated, send a smear and up to 3 ml of CSF in pediatric blood culture bottle. Transport ASAP. PCR should not replace culture.			
Methicillin-resistant <i>Staphylococcus aureus</i> (MRSA) (See Staphylococcus colonization and Staphylococcus infections)			
Microsporidiosis	Feces in SAF	Microscopy for spores	PA
NOTE: Requires special request on requisition for special staining.			
Molluscum contagiosum (see Poxviruses)			
Moniliasis (see candidiasis)			



Disease or Syndrome (Causal Agent(s))	Specimen Requirements	Tests / Examinations	Section
Mononucleosis, infectious (see infectious mononucleosis)			
MRSA (See Staphylococcus colonization and Staphylococcus infections)			
Mucormycosis (see Phycomycosis)			
Mumps, including complicating meningoencephalitis, pancreatitis or orchitis (Mumps virus)	Throat swab in VTM CSF, Saliva Urine if orchitis present  Clotted blood or serum	Viral culture NAAT (Referred out)  Serology	VD  SE
NOTES: Presence of IgM is diagnostically significant.			
Myalgia, epidemic (see Coxsackievirus infections)			
Mycetoma (see also Maduromycosis) ( <i>Actinomyces israelii</i> , <i>Exophiala jeanselmei</i> , <i>Nocardia species</i> , <i>Streptomyces species</i> , <i>Pseudoallescheria</i> <i>boydii</i> , <i>Madurella spp.</i> , and other filamentous fungi)	Pus with or without granules Aspirated material from fluctuant areas Biopsy material using a wide-mouthed, screw- capped sterile container	Microscopy and fungal culture	CM
NOTE: Transport ASAP and refrigerate if a delay in transport is anticipated.			

Disease or Syndrome (Causal Agent(s))	Specimen Requirements	Tests / Examinations	Section
Mycobacteria, atypical (see also Tuberculosis) ( <i>Mycobacterium</i> <i>kansasii</i> , <i>M.</i> <i>scrofulaceum</i> , <i>M.</i> <i>avium-intracellulare</i> <i>M. marinum</i> ( <i>balnei</i> ), <i>M. fortuitum</i> , and others)	Sputum Swabs of skin lesions or pus in TM Stool, Blood, Body fluids, Tissue, Bone marrow	Referred out to HSC	CM
NOTES: Please specify that examinations for these organisms are required. Referred to the HSC TB lab. For further information call the HSC TB lab at 787-7652.			
Mycoplasma infections ( <i>Mycoplasma</i> <i>pneumoniae</i> )	Sterile fluids, tissue Respiratory secretions	PCR (Referred out)	CM
Mycoplasma pneumonia IgM	Clotted blood or serum	Serology – CF	SE
Notes: <i>M. pneumoniae</i> causes primary atypical walking pneumonia. <i>M.</i> <i>pneumoniae</i> is tested for IgM on only acute samples.			
Mycoplasma Infections – Genital Mycoplasmas ( <i>Mycoplasmas hominis</i> , <i>Ureaplasma urealyticum</i> )	Sterile fluids and tissue – neonates and children Placental swab, amniotic fluid Respiratory secretions (neonates and children) – NOT SPUTUM Urethral/cervical swab* Urines/semen**	Culture/PCR (Referred out)	CM
* Relevant history required for culture – physician must consult St. Boniface Microbiology lab at 237-2484 to arrange for testing.			
** Relevant history required for culture, e.g., prostatitis, unresolving UTI following treatment – consult St. Boniface Microbiology lab at 237-2484.			
Notes: For culture, specimens must arrive at CPL same day as collected, early enough to allow for shipment and arrival at reference lab on same day. Consult CPL at 945-7184.			
Myocarditis (Coxsackie B and other enteroviruses)	Throat swab in VTM Feces Pericardial fluid	Viral culture	VD

Disease or Syndrome (Causal Agent(s))	Specimen Requirements	Tests / Examinations	Section
Necrotizing fasciitis (Group A <i>Streptococcus</i> )	Swab in TM Tissue Isolate	Culture  Typing (Referred out)	CM
Nephritis, acute glomerulo- (see also Streptococcal infections) (Sequelae of <i>Streptococcus</i> <i>pyogenes</i> infections)	Nose, throat or skin swabs in TM  Clotted blood or serum	Microscopy and culture  ASOT Anti-DNase B	CM  SE
NOTES: Certain M serotypes (e.g. type 12) of <i>S. pyogenes</i> are associated with nephritis; strains isolated from cases will generally be typed but results may not be available for weeks.			
Nocardiosis ( <i>Nocardia asteroides</i> <i>Nocardia spp.</i> )	Sputum Pleural fluid Pus in TM	Microscopy and culture	CM
NOTES: Grows slowly in culture. Specify on the requisition if no cardiosis is suspected.			
Non-specific urethritis (see Trichomoniasis, <i>Herpes simplex</i> , <i>Candidiasis</i> , <i>Mycoplasma</i> , <i>Chlamydia</i> <i>trachomatis</i> )			
Onchocerciasis (see also Filariasis; ( <i>Onchocerca volvulus</i> , <i>Mansonella</i> <i>streptocerca</i> )	Biopsy of skin Aspirated material from skin nodules Excision of nodule	Microscopy examination for microfilariae and search for adult worm	PA
Onychomycosis (see Dermatophytosis)			
Ophthalmia neonatorum (see Conjunctivitis)			
Orchitis, viral (see Mumps)			

Disease or Syndrome (Causal Agent(s))	Specimen Requirements	Tests / Examinations	Section
Orf (see Poxvirus infections)			
NOTES: Usually transmitted from sheep to man.			
Ornithosis (see Chlamydia infections: respiratory)			
Osteomyelitis, acute (see also Staphylococcal infections) ( <i>Staphylococcus aureus</i> and other bacterial species)	Blood for culture Purulent discharge from skin or other lesions in TM Aspirated pus in TM Swab from primary lesion in TM	Microscopy and culture	CM
Otitis media (see Ear infections)			
Otomycosis (see Ear infections)			
Pancreatitis, viral (Coxsackie B virus Mumps virus)	Stool for Coxsackie Throat swab in VTM and urine for mumps	Viral culture NAAT (Referred out)	VD
Panencephalitis, subacute sclerosing (see also Measles) (Probably associated with measles virus infections)	Clotted blood or serum CSF, Brain biopsy culture Postmortem specimen	Serology Viral culture (Referred out)	SE VD
NOTE: Consult the laboratory for brain biopsy culture.			
Papilloma virus (Warts-Epidermal (genital) Uterine Cervical Dysplasia, Carcinoma)	Tissue biopsy Exfoliated cervical cells from transformation zone	NAAT (Referred out)	VD
Paracoccidiomycosis ( <i>Paracoccidioides brasiliensis</i> )	Mouth (or lip) swab and scrapings, Pus Sputum Biopsy material	Microscopy and fungal culture	CM
	Clotted blood or serum	Serology (Referred out)	SE
NOTE: Please indicate presumptive paracoccidioidomycosis on requisition.			

Disease or Syndrome (Causal Agent(s))	Specimen Requirements	Tests / Examinations	Section
Paragonimiasis ( <i>Paragonimus westermani</i> )	Feces in SAF	Microscopy examination for ova	PA
Parainfluenza virus infections including colds, pharyngitis, laryngitis, bronchiolitis, pneumonia (Parainfluenza viruses)	Nasopharyngeal aspirate Throat swab in VTM	Viral culture	VD
NOTE: Infection may give an anamnestic response to other paramyxoviruses, e.g. mumps.			
Paralytic illnesses caused by viruses (see also Encephalitis, viral, and individual viruses) (Several viruses, especially polio- coxsackie-, echo-, and herpesviruses, and as part of encephalomyelitis, or ascending myelitis syndromes) (Western Equine Encephalitis) WEE, (St. Louie Encephalitis) SLE (West Nile Virus) WNV	Throat swab in VTM Biopsy material (brain, spinal cord) Autopsy material in VTM (brain or spinal cord) Cerebrospinal fluid Feces	Viral culture NAAT	VD
NOTES: Consult the CPL Virology section for culture of biopsy material. No serology test available for polio.			
Paratyphoid fever ( <i>Salmonella</i> <i>paratyphi</i> A, B, or C)	Feces Urine Blood Isolate	Culture  Typing	CM
Paronychia, mycotic ( <i>Candida albicans</i> )	Nail scrapings (base of nail)	Microscopy and fungal culture	CM
Parvovirus B <sub>19</sub> (see Erythema infectiosum)			

Disease or Syndrome (Causal Agent(s))	Specimen Requirements	Tests / Examinations	Section
<i>Pasteurella</i> infections (see also <i>Yersinia</i> infections; Plague) ( <i>Pasteurella multocida</i> )	Swab from wound in TM Sputum Blood for culture	Culture	CM
NOTES: Often transmitted by bites of animals; also causes chronic respiratory infection especially in persons with prolonged contact with animals. Relevant clinical information very helpful.			
Pediculosis ( <i>Pediculus humanus capitus</i> , <i>Pediculus humanus corporis</i> , <i>Phthirus pubis</i> )	Parasites or ova in hair or under-clothing	Microscopy identification	PA
NOTE: May be submitted in alcohol.			
Pemphigus ( <i>Staphylococcus aureus</i> )	Swab of vesicle fluid in TM	Microscopy and culture	CM
Pericarditis, viral (Coxsackie B viruses)	Feces Pericardial fluid	Viral culture	VD
Pertussis ( <i>Bordetella pertussis</i> , <i>B. parapertussis</i> )	Nasopharyngeal swab in TM or Auger suction	Culture PCR	CM
	Clotted blood or serum	Serology	SE
NOTES: Organisms are very fastidious and difficult to isolate. Transport ASAP. Serology for outbreaks or immune status only.			
Pharyngitis, bacterial (see also streptococcal infections; Diphtheria) (Several species of bacteria, especially <i>Streptococcus pyogenes</i> , <i>Corynebacterium diphtheriae</i> , <i>C. ulcerans</i> )	Throat swab in TM	Culture	CM
NOTE: Indicate any suspicion of diphtheria on requisition.			



Disease or Syndrome (Causal Agent(s))	Specimen Requirements	Tests / Examinations	Section
Pharyngitis, viral (Many viruses)	Throat swab in VTM Feces	Viral culture	VD
Phenylketonuria (PKU)	Newborn screening card	Fluorometry	NBS/ PHC
Phycomycosis (Fungi in the class of Phycomycetes)	Scrapings from lesion in black paper or dry sterile container.	Microscopy and fungal culture	CM
NOTES: Most frequently associated with diabetes mellitus and other severe debilitating diseases, or during treatment with corticosteroids and antimetabolites.			
Pinta ( <i>Treponema carateum</i> )	Clotted blood or serum	Serology tests for syphilis	SE
NOTES: Antibodies to the treponemes of pinta are indistinguishable from those to the treponemes of syphilis by all the diagnostic tests in current use.			
Pinworm infection (see Enterobiasis)			
Pityriasis versicolor ( <i>Malassezia furfur</i> )	Skin scrapings or scales	Microscopy	CM
NOTE: Culture confirmation is not practical.			
Plague ( <i>Yersinia pestis</i> )	Pus from buboes in TM	Culture	CM
	Throat swab in TM, Sputum Isolate Blood culture	Confirmation	
	Clotted blood or serum	Serology (Referred out)	SE
NOTES: Consult the laboratory if plague, which is endemic in the southwestern USA, is suspected. Label all specimens 'SUSPECT PLAGUE'.			
Pleurodynia, epidemic (see Coxsackievirus infections)			
Pneumonia, primary atypical (see Mycoplasma infections)			

Disease or Syndrome (Causal Agent(s))	Specimen Requirements	Tests / Examinations	Section
Pneumonia, bacterial (Several bacterial species, especially <i>Streptococcus pneumoniae</i> , <i>Streptococcus pyogenes</i> , <i>Staphylococcus aureus</i> , <i>Klebsiella pneumoniae</i> , and <i>Haemophilus influenzae</i> )	Sputum Auger suction or transtracheal aspirate Throat swab in TM Lung aspirate Biopsy	Microscopy and culture	CM
Pneumonia, viral (Numerous viruses)	Nasopharyngeal swab or aspirate in VTM Throat swab in VTM Autopsy material (lung) in VTM	Viral culture RSV EIA NAAT	VD
Poliomyelitis (Polioviruses, types 1, 2, and 3)	Feces Cerebrospinal fluid Autopsy material (brain, spinal cord, intestinal contents) in VTM	Viral culture NAAT	VD
Poxvirus infections (Poxviruses: cowpox, vaccina, contagious pustular dermatitis (Orf). Milker's nodes (paravaccinia), and molluscum contagiosum))	Vesicle fluid, Exudate from skin lesions, Skin crusts, Scrapings from skin lesion in VTM Lesion smear	Electron microscopy Viral culture  Electron Microscopy only	VD
NOTES: The differentiation between poxviruses and varicella (herpes virus) can be made within hours by electron microscopy.			
Psittacosis (see Chlamydia infections: respiratory)			

Disease or Syndrome (Causal Agent(s))	Specimen Requirements	Tests / Examinations	Section
Puerperal fever (Usually haemolytic streptococci)	High vaginal swab in TM Nose and throat swab in TM from mother Umbilical swab in TM from baby	Culture	CM

NOTES: If haemolytic streptococcal infection is suspected, take nose and throat swabs from the mother, and umbilical swab from baby.

---

Pyelitis/Pyelonephritis  
(see Urinary tract  
infections)

---

Pyoderma ( <i>Streptococcus</i> <i>pyogenes</i> <i>Staphylococcus</i> <i>aureus</i> )	Swab in TM	Culture	CM
---	------------	---------	----

---

Pyrexia of unknown origin (PUO) (Infection due to various bacteria and other agents)	Feces, Urine	Culture	CM
--	--------------	---------	----

Clotted blood or serum	Serologic tests for enteric fever, brucellosis, tularemia and other infections	SE
------------------------	---	----

---

Q fever ( <i>Coxiella burnetii</i> )	Clotted blood or serum	Serology (Referred out)	SE
---	---------------------------	----------------------------	----

---

Rabies (Rabiesvirus)	Clotted blood or serum	Serology (Referred out)	SE
	Biopsy material	Referred out	

---

NOTE: Rabies test for immune status testing. Consult with CPL to arrange biopsy collection and transport.

---

Disease or Syndrome (Causal Agent(s))	Specimen Requirements	Tests / Examinations	Section
Reiter's syndrome (see Non-specific urethritis)			
Relapsing fever, louse or tick-borne ( <i>Borrelia spp.</i> )	Blood films Blood in citrate  Serum	Microscopy EIA  Serology	PA  SE
Reovirus infections including upper respiratory infections and diarrhea (Reoviruses, types 1, 2, and 3)	Feces Throat swab in TM	Viral culture EM	VD
Respiratory infections, acute bacterial (see also under individual syndromes and diseases) (Several bacterial species; especially <i>Streptococcus pyogenes</i> , <i>Corynebacterium diphtheriae</i> , <i>Streptococcus pneumoniae</i> , and <i>Haemophilus influenzae</i> )	Sputum Auger suction Throat swab in TM  Clotted blood or serum	Microscopy and culture  Serology (Referred out)	CM  SE
NOTES: Please refer to note for "Pneumonia, bacterial".			
Respiratory infections, acute viral (see also under individual syndromes and diseases) (Numerous viruses, especially adeno-, rhino-, coxsackie-, influenza, and respiratory syncytial virus)	NPA, NPS Throat swab, Autopsy material (lung) in VTM	Viral culture EIA NAAT	VD

Disease or Syndrome (Causal Agent(s))	Specimen Requirements	Tests / Examinations	Section
Respiratory syncytial virus infections (Respiratory syncytial virus)	Nasopharyngeal aspirate in VTM NPS, ETT Autopsy material (lung) in VTM	Viral culture EIA	VD
NOTES: Virus unstable and specimens should NOT BE FROZEN, but should be kept cool (e.g. on cold pack). Transport to CPL as soon as possible after collection.			
Rheumatic fever (see also Streptococcal infections) (A sequelae to infection with <i>Streptococcus</i> <i>pyogenes</i> )	Throat swab in TM Clotted blood or serum	Culture Serology	CM SE
Rhinovirus infections (see common cold)			
Rickettsial infections, Louse-borne typhus, Rocky Mountain Spotted Fever, Rickettsialpox Scrub typhus (see also Q fever) ( <i>Rickettsia prowazekii</i> <i>R. typhi</i> , <i>R. rickettsii</i> <i>R. akari</i> , <i>Orientia</i> <i>tsutsugamushi</i> Various other species)	Clotted blood or serum	Serology (Referred out)	SE
Ringworm (see Dermatophytosis)			
Rocky Mountain Spotted Fever (See Rickettsial infections)			
Rotavirus	Feces	Electron microscopy	VD

Disease or Syndrome (Causal Agent(s))	Specimen Requirements	Tests / Examinations	Section
Rubella, German measles (Rubella virus)	Clotted blood, serum or plasma	Serology	SE
	Aborted material Placenta Throat swab in VTM Urine	Viral culture	VD
Rubella, congenital rubella syndrome (Rubella virus)	Urine Nasal swab, Throat swab, Autopsy material (all organs) in VTM Cerebrospinal fluid	Viral culture	VD
	Clotted blood or serum	Serology	SE
NOTES: Infants with congenital rubella infection may excrete virus in the urine for many months after birth.			
Rubeola (see Measles)			
Salmonella (see also Typhoid fever, paratyphoid fever, and food poisoning) ( <i>Salmonella</i> spp., over 2000 named serotypes)	Feces, Blood for culture	Culture	CM
	Urine Isolate	Typing	
NOTES: In suspected typhoid/paratyphoid fever, cultures of blood and urine are indicated. Outbreaks should be noted on requisition. Contact the regional MOH to establish outbreak or FBI status and obtain outbreak code for clinical specimens.			
Salpingitis (see Gonorrhea)			
SARS (Severe acute respiratory syndrome) (SARS coronavirus)	Clotted blood or serum	Serology (Referred out)	SE
	Nasopharyngeal aspirate Stool	Viral Culture RT-PCR Electron microscopy	VD
NOTES: SARS investigation referred to NML. Please contact CPL prior to submission of specimens.			
Scabies ( <i>Sarcoptes scabiei</i> )	Scrapings of skin at edge of tracks	Microscopy examination for mites	PA



Disease or Syndrome (Causal Agent(s))	Specimen Requirements	Tests / Examinations	Section
Scarlet fever, scarlatina (see Streptococcal infections)			
Schistosomiasis ( <i>Schistosoma</i> <i>haematobium</i> <i>Schistosoma</i> <i>japonicum</i> <i>Schistosoma</i> <i>mansoni</i> )	Urine, for <i>S.haematobium</i> Feces in SAF Clotted blood or serum	Microscopy examination for ova  Serology (Referred out)	PA
Schistosomal dermatitis (Swimmer's itch) ( <i>Trichobilharzia</i> species)		No useful test	
NOTES: Common in lakes in North America. Notify regional MOH if suspected.			
Scrub typhus (see Rickettsial infections)			
Septicemia (Numerous bacteria)	Swabs of septic lesions in TM Urine	Culture	CM
NOTES: Cultures should be taken from any suspected focus of infection.			
Serum hepatitis (see Hepatitis B)			
Severe acute respiratory syndrome (SARS) (see SARS)			
Shigellosis ( <i>Shigella flexneri</i> , <i>Shigella sonnei</i> , <i>S. dysenteriae</i> , <i>S. boydii</i> )	Feces Rectal swab in TM Isolate	Culture  Typing	CM
Shingles, zoster (see Chickenpox)			
Sore throat (see Pharyngitis)			

Disease or Syndrome (Causal Agent(s))	Specimen Requirements	Tests / Examinations	Section
--	--------------------------	----------------------	---------

South American  
Blastomycosis (see  
Paracoccidioidomycosis)

Sporotrichosis ( <i>Sporothrix</i> , <i>Sporotrichium</i> <i>schenckii</i> )	Pus from ulcerated lesions or aspirated from subcutaneous abscesses in TM	Microscopy and fungal culture	CM
---	--	----------------------------------	----

Spotted fever (see  
Rickettsial infections)

St. Louis encephalitis  
(see Arbovirus  
infections)

Staphylococcus colonisation (MRSA, Methicillin resistant <i>Staphylococcus</i> <i>aureus</i> ) (See also Staphylococcus infections)	Swabs from nares, throat, rectum or ostomy, wounds, line or device sites in TM Isolate	MRSA Screen  Confirmation	CM
---	--	---------------------------------	----

NOTES: MRSA isolates may also be forwarded to CPL for pulsed-field gel  
electrophoresis (PFGE).

Staphylococcus infections ( <i>Staphylococcus</i> <i>aureus</i> and some members of the coagulase-negative Staphylococci) (See also Food poisoning) (See also Staphylococcus colonisation)	Wound or nasal swab in TM Infected body fluids Urine CSF	Microscopy and culture	CM
--	--	------------------------	----

NOTES: *Staphylococcus aureus* is commonly present in the nose and is reported in  
nasal swabs only if the requisition indicates that there is a lesion of the  
nasal passages or in the search for a carrier in relation to food poisoning  
or hospital outbreaks.

Disease or Syndrome (Causal Agent(s))	Specimen Requirements	Tests / Examinations	Section
Streptococcal infections including scarlet fever (scarlatina), erysipelas, epidemic streptococcal sore throat and streptococcal-related illnesses such as acute glomerulonephritis, rheumatic fever. Necrotising fasciitis, toxic shock syndrome ( <i>Streptococcus pyogenes</i> and other groups of streptococci)	Throat swab in TM Exudate from infected area  Clotted blood or serum	Culture  Serology	CM  SE
NOTES: The current serologic tests are the ASOT and the ADB titres. A result of 200 or greater is diagnostically significant.			
Group B <i>Streptococcus agalactiae</i> (prenatal screen)	Combination vaginal/rectal swab	Culture	CM
NOTES: Prenatal screen recommended at 35-37 weeks gestation. Culture performed only with clinical information indicating pregnancy. As this is a transient coloniser, a negative culture does not guarantee absence of GBS at time of labour.			
Strongyloidiasis ( <i>Strongyloides stercoralis</i> )	Feces in SAF  Clotted blood or serum	Microscopy examination for larvae  Serology	PA  SE
NOTES: Single blood is adequate.			
Subacute sclerosing panencephalitis (See Panencephalitis)			
Swimmer's itch (see Schistosomal dermatitis)			

Disease or Syndrome (Causal Agent(s))	Specimen Requirements	Tests / Examinations	Section
Syphilis ( <i>Treponema pallidum</i> )	Serous exudate from primary lesion (chancre) or secondary lesion of skin or mucous membrane as a lesion slide Cerebrospinal fluid Clotted blood or serum	Immuno-fluorescent examination for <i>T.</i> <i>pallidum</i>  Serology	SE

NOTES: Note that certain other treponemal infections, notably yaws, cause positive reactions in ALL the serologic tests used in the diagnosis of syphilis.

---

Tapeworms (see  
Taeniasis and  
Diphyllobothriasis,  
Hymenolepiasis,  
Echinococcosis)

Taeniasis ( <i>Taenia saginata</i> <i>Taenia solium</i> )	Feces in SAF Worm, including segments	Microscopy examination for ova, and identification of segments	PA
---	---	---	----

NOTES: *T. saginata* - beef tapeworm.  
*T. solium* - pork tapeworm.

Tetanus ( <i>Clostridium tetani</i> )	Swabs from wounds and other lesions in TM	Microscopy and culture	CM
	Clotted blood or serum	Serology (Referred out)	SE

NOTES: See anaerobic culture for specimen submission details. Serology for immune status testing only.

---

Throat infections  
(see Pharyngitis)

---

Thrush (see  
Candidiasis)

---

Toxocariasis (see  
Visceral larva migrans)

Toxoplasmosis	Clotted blood or serum Biopsy Contact Parasitology Lab	Serology  Microscopy	SE  PA
---------------	--	----------------------------	--------------

Disease or Syndrome (Causal Agent(s))	Specimen Requirements	Tests / Examinations	Section
Trachoma (see Chlamydial Infections <i>Chlamydia trachomatis</i> )			
Trichinosis ( <i>Trichinella spiralis</i> )	Biopsy of muscle	Examination for larvae	PA
	Serum (Preferred)	Serology (Referred out)	SE
NOTES: A single serum sample is adequate.			
Trichomoniasis ( <i>Trichomonas vaginalis</i> )	Air dried slide from Vagina or Urethra Prostatic secretions after massage	Microscopy	CM
NOTES: Secretions dried on glass slides can be examined by microscopy. Do not apply fixative.			
Trichostrongyliasis (see also Hookworm disease) ( <i>Trichostrongylus species</i> )	Feces in SAF	Microscopy examination for ova	PA
NOTES: Occurs in Russia and the Orient.			
Trichuriasis ( <i>Trichuris trichiura</i> , whipworm)	Feces in SAF	Microscopy examination	PA
Triple Test (AFP/uE3/hCG)	Clotted blood or serum	Chemiluminescent immunoassay	NBS/ PHC
Trypanosomiasis, African ( <i>Trypanosoma rhodesiense</i> , <i>T. gambiense</i> )	Blood films, thick and thin Lymph node aspirated	Microscopy examination	PA
	Clotted blood or serum	Serology (Referred out)	SE
NOTES: Occurs in Tropical Africa.			
Trypanosomiasis, American (see Chagas' disease)			

Disease or Syndrome (Causal Agent(s))	Specimen Requirements	Tests / Examinations	Section
Tuberculosis ( <i>Mycobacterium tuberculosis</i> , <i>M. Bovis</i> )	Sputum Fluids from body cavities, joints, etc. Biopsy material (tissue, lymph glands, uterine curettings) Gastric washings Cerebrospinal fluid Purulent exudate in TM Urine, three early morning specimens	Referred out to HSC	CM

NOTES: Send early morning sputa on three consecutive days. All specimens are referred to HSC TB lab for culture. Positive results are sent as soon as available, but negatives are held for 8 weeks prior to being reported. For patients without spontaneous sputum, induction of cough and sputum by inhalation of a warm sterile aerosol of saline is preferred to gastric aspiration. Submission of both induced sputum and gastric contents gives the best results. For further information call the HSC TB lab at 787-7652.

Tularemia ( <i>Francisella tularensis</i> )	Clotted blood or serum	Serology (Referred out)	SE
	Swab from ulcer in TM Aspirated material from lymph nodes Blood for culture Suspected isolate	Microscopy and culture	CM

NOTE: Culture examination only by special arrangement with the laboratory. Requisition should indicate "Suspect tularemia".

Typhoid fever (Enteric fever) ( <i>Salmonella typhi</i> )	Feces Urine Blood	Culture	CM
Typhus fever (see Rickettsial infections)			
Undulant fever (see Brucellosis)			
Urethritis (see Gonorrhea, Candidiasis, Chlamydia, Trichomoniasis, Herpes simplex, Mycoplasma)			



Disease or Syndrome (Causal Agent(s))	Specimen Requirements	Tests / Examinations	Section
Urinary tract infections (Numerous bacterial species, including <i>E. coli</i> , <i>Klebsiella spp.</i> , <i>Proteus species</i> , <i>Pseudomonas aeruginosa</i> , <i>Enterococcus spp.</i> , <i>Staphylococcus aureus</i> , <i>S. saprophyticus</i> )	Midstream urine on Dip slides	Culture and sensitivity (semi-quantitative count)	CM
NOTES: Use dip slides. It is essential to indicate whether the patient is or recently has been on antibiotic therapy in order to assist interpretation of the findings. In general, counts of 10 <sup>8</sup> CFU/L or more are significant.			
Vaccinia (see Poxvirus infections)			
Vaginitis (see Gonorrhea, Trichomoniasis Streptococcal infections, Candidiasis, Chlamydia infections)			
Varicella (see Chickenpox)			
Vancomycin-resistant Enterococcus (VRE) (See Enterococcus colonization and Enterococcus infections)			
Venezuelan equine encephalitis (see Arbovirus infections)			
Vibrio infections (Cholera is listed separately) (see also Food poisoning) ( <i>Vibrio fetus</i> and related Vibrios)	Cerebrospinal fluid	Microscopy and culture	CM

Disease or Syndrome (Causal Agent(s))	Specimen Requirements	Tests / Examinations	Section
Vincent's angina ( <i>Fusobacterium fusiforme</i> concomitant with <i>Borrelia vincenti</i> )	Air dried smear from gums	Microscopy	CM
Viral hepatitis (see Hepatitis A, B, C, D, E)			
Visceral larva migrans ( <i>Toxocara spp.</i> )	Clotted blood or serum	Serology (Referred out)	SE
NOTE: There are considerable cross-reactions in serologic tests with <i>Toxocara</i> and <i>Ascaris</i> antigens.			
Weil's disease (see Leptospirosis)			
Western equine encephalitis (see Arbovirus infections)			
West Nile Virus (WNV) infections (See also Arbovirus infections)	Clotted blood or serum	Serology	SE
NOTE: Serum in all cases is the specimen of choice.			
Whooping cough (see Pertussis)			
Worm infections (see also under individual parasites)	If possible send whole worm in saline Fix in 10% formalin	Identification	PA
Wound infections (Different species of aerobic and anaerobic bacteria)	Wound swab in TM	Microscopy and culture	CM
NOTE: See wound specimen collection.			
Yaws ( <i>Treponema pertenuae</i> )	Clotted blood or serum	Serological tests for syphilis	SE
NOTE: Antibodies to the treponeme of yaws are indistinguishable from those to the treponeme of syphilis by all the diagnostic tests in current use.			
Yellow fever (see also Arbovirus infections)	EDTA blood	PCR (Referred out)	VD
NOTE: For post-vaccine illness only.			

Disease or Syndrome (Causal Agent(s))	Specimen Requirements	Tests / Examinations	Section
Yersinia infections (see also Diarrhea) ( <i>Yersinia Enterocolitica</i> , <i>Y. pseudotuberculosis</i> )	Blood for culture Feces, Swab from abscess in TM Excised mesenteric lymph nodes	Culture	CM
	Clotted blood or serum	Serology (Referred out)	SE

NOTE: May cause enteritis, terminal ileitis and mesenteric lymphadenitis or chronic enteritis.

---


Zoster (see  
Chickenpox)

---

Zygomycosis (see  
Phycomycosis)

---

## 8.0 FORMS AND REQUISITIONS

**Manitoba** 

Health

Cadham Provincial Laboratory  
Provincial Programs & ServicesP.O. Box 8450  
750 William Avenue  
Winnipeg MB R3C 3Y1  
PH: (204) 945-6806  
Fax: (204) 786-4770**SUPPLIES REQUEST FORM (PLEASE FAX ONLY)**

ORDERED BY: \_\_\_\_\_

FACILITY NAME: \_\_\_\_\_ DATE: \_\_\_\_\_

ADDRESS: \_\_\_\_\_ PHONE: \_\_\_\_\_

CITY: \_\_\_\_\_ PROVINCE: ☐ MB ☐ SK ☐ ON POSTAL CODE: \_\_\_\_\_SEND BY: PICK-UP ☐ BUS ☐ AIR ☐ BLUE BOX ☐ MAIL ☐PUROLATOR ☐ MEDICAL CARRIERS ☐ OTHER \_\_\_\_\_

<u>FORMS</u>	<u>QUANTITY</u>	<u>SWABS</u>	<u>QUANTITY</u>	<u>OTHER</u>	<u>QUANTITY</u>
General Requisition	_____	A.C.M.	_____	Stool Container	_____
Address Labels	_____	Rayswab (100/Box)	_____	Urine Container	_____
Newborn Screening Cards	_____	GenProbe (yellow) urine collection kit	_____	VTM (Viral Transport Medium)	_____
HIV Requisition	_____	GenProbe (purple) unisex swab kit	_____	LGV Supplies (Dacron swab & 2SP CTM)	_____
Retrovirus Requisition	_____	Micro Trak	_____		
Guide to Services	_____	Nasopharyngeal	_____		
Protecting Your Baby pamphlets	_____	Flocked swab (adult)	_____		
		Flocked swab (pediatric)	_____		

**MEDIA**

BAsp	_____	O/F Glucose	_____	SIM	_____
BAw	_____	Urea Slants	_____	Orn. Decarb.	_____
BBE	_____	BA Slants	_____	Mineral Oil	_____
BCA	_____	CHOC Slants	_____	Lysine Decarb	_____
BMH	_____	GC Glucose	_____	Nitrate Slants	_____
BHI	_____	GC Maltose	_____	Skim Milk	_____
CHOC	_____	BE	_____	Calf Serum	_____
MH	_____	NaCl	_____	Cooked Meat	_____
NM	_____	Todd Hewitt Scr.	_____	Maintenance Media	_____
OM	_____	BHI for Anaerobes	_____	McFarland Standards (0.5)	_____
PEA	_____	SIMCIT	_____		
TM	_____				

**REAGENTS**

Gram's Crystal Violet	_____	10% KOH	_____	Kovac's	_____
Gram's Safranin	_____	40% KOH	_____	Ehrlich's	_____
Gram's Iodine	_____			SAF Fixative	_____
Decolorizer 50/50	_____				

**PLEASE FAX THE FORM TO : 786-4770. Allow 48 hrs to complete the order.**

**Cadham Provincial Laboratory**  
**Diagnostic and Infectious Specimen Transport Guidelines**

<b><u>Transportation Mode</u></b>	<b><u>Diagnostic Specimens</u></b>	<b><u>Infectious Specimens</u></b>
<b>Surface Transport</b>	Packaging : TC-125-1B  Packing Instruction: IATA 650  Identify as Diagnostic Specimens  Documentation : Waybill	Packaging : TC-125- 1A  Packing Instruction: IATA 602  Class 6.2 Dangerous Goods Risk Groups 2, 3, 4 UN 2814 or UN 2900  Documentation : Shipper's Declaration / Waybill  <ul style="list-style-type: none"> <li>- Specimens known or suspected of containing viable micro-organisms</li> <li>- Viral loads</li> <li>-Cultures or isolates</li> <li>- Specimens sent for the diagnosis of an RG 4 organism</li> </ul>
<b>Air Transport</b>	<b><u>Category B Specimens</u></b>  Packaging : TC-125-1B or equivalent  Packing Instruction: IATA 650  UN 3373 Diagnostic Specimens or Clinical Specimens  Documentation : Waybill  -Diagnostic or infectious specimens that do not meet the criteria for inclusion in Category A	<b><u>Category A Specimens</u></b>  Packaging : TC-125-1A  Packing Instruction: IATA 602  Class 6.2 Dangerous Goods UN 2814 or UN 2900  Documentation : Shipper's Declaration / Waybill  <ul style="list-style-type: none"> <li>-Micro-organisms listed in Table 3.6.D, p. 96-97, 46<sup>th</sup> Ed. IATA DGR.</li> <li>- Pathogens meeting the same criteria</li> <li>- Pathogens amplified by culture</li> </ul>

## **BLUE BOX PACKAGING DIRECTIONS**

---

When packaging **DIAGNOSTIC SPECIMENS FOR TRANSPORT:**

### **MICROBIOLOGY and VIROLOGY specimens:**

- Put diagnostic specimen into a sealable ziplock-specimen bag with requisition pouch.
- Put tissue specimens for viral studies in viral transport media vials and into a sealable ziplock-specimen bag and seal.
- Put tissue specimens for micro studies in sterile screw cap containers with a small amount of sterile saline into a sealable ziplock-specimen bag and seal.
- Place requisition in the outside pouch of the sealable ziplock-specimen bag.
- Put bag, with specimen upright, into foam inserts.
- Place foam inserts into the cooler, ensuring that white absorbent pads line the bottom of the cooler. Securely close blue packaging box/cooler when finished.

### **URINE and STOOL specimens:**

- Put tightly closed specimen containers into a sealable ziplock-specimen bag and seal.
- Place requisition in the outside pouch of the sealable ziplock-specimen bag.
- Put bag, with specimen upright, into foam inserts.
- Place foam inserts into the cooler, ensuring that white absorbent pads line the bottom of the cooler. Securely close blue packaging box/cooler when finished.

### **DIAGNOSTIC BLOOD specimens in venipuncture tubes:**

- Put collection tube into foam insert.
- Put requisitions in corresponding order as per specimen location in foam insert.
- Place all requisitions in a large sealable ziplock-specimen bag and place in cooler.
- Place foam inserts into the cooler, ensuring that white absorbent pads line the bottom of the cooler.
- Securely close blue packaging box/cooler when finished.

### **NEWBORN SCREENING COLLECTION CARDS:**

- Ensure each card is placed within its protective plastic sleeve.
- Place cards into a sealable ziplock-specimen bag and seal.
- Place bag in the main body of the blue box and NOT BETWEEN the box and styrofoam liner.
- Securely close blue packaging box/cooler when finished.



## **FROZEN specimens:**

- > Put diagnostic specimen into foam insert.
- > Put requisitions in corresponding order as per specimen location in foam insert and seal in a separate sealable ziplock-specimen bag.
- > Place this bag over coolpack/icepack and put both into a sealable ziplock-specimen bag.
- > Make sure that the ice pack is not in direct contact with the sample or requisition so it does not get wet if it thaws.
- > If crushed ice is used, have the crushed ice inside the bag with the specimen in the outer pocket wrapped around by the ice. Then place this in a second specimen bag with the requisition in the outer pocket.
- > Put into blue packaging box/cooler, ensuring that white absorbent pads line the bottom of the cooler/blue packaging box.
- > Securely close blue packaging box/cooler when finished.

## **CYTOLOGY specimens:**

- > Put diagnostic specimen into 90 ml Starplex container. Depending on specimen type, there may be about 40 ml of 50% ethanol or commercial Cytospin Collection Fluid in the container.
- > Securely close the container and place in a sealable ziplock-specimen bag and seal bag.
- > Place requisition in the outside pouch of the sealable ziplock-specimen bag.
- > Put into blue packaging box, ensuring that white absorbent pads line the bottom of the cooler.
- > Securely close blue packaging box/cooler when finished.

## **HISTOLOGY specimens:**

### **1. Tissue cassettes**

- > Put tissue cassettes in cassette holder.
- > Put cassette holder into a primary water-tight plastic container.
- > Put about 300 ml 10% formalin into primary container
  - JUST enough to keep cassettes covered.
- > Put primary water-tight plastic container into a secondary water-tight plastic container, ensuring that the secondary container is lined with a BLUE mini-fan pad.
- > Put into blue packaging box/cooler, securely close blue packaging box/cooler when finished.

### **2. Slides**

- > Place slides in large cardboard slide carrier. Close cover.
- > Put into blue packaging box, securely close blue packaging box/cooler when finished.

### **3. Pathology specimens**

- > Put diagnostic specimen into 90 ml Starplex container containing about 30 ml of 10% formalin.
- > Securely close the container and place in a sealable ziplock-specimen bag and seal bag.
- > Place requisition in the outside pouch of the sealable ziplock-specimen bag.
- > Put into blue packaging box, ensuring that BLUE fan-pads line the bottom of the cooler.
- > Securely close blue packaging box/cooler when finished.

**REQUEST TO AMEND REQUISITION INFORMATION SUBMITTED TO  
CADHAM PROVINCIAL LABORATORY**

**Please amend the above requisition to reflect the corrected data as provided below:**  
Please print clearly.

**PATIENT:** Surname \_\_\_\_\_ Given Name \_\_\_\_\_ Gender \_\_\_\_\_

DOB (yyyy/mm/dd) \_\_\_\_\_ MB Health No. \_\_\_\_\_ \*\*HOSPITAL/CLINIC # \_\_\_\_\_

Address \_\_\_\_\_

Physician \_\_\_\_\_ Facility \_\_\_\_\_

**REQUEST TO AMEND  
REQUISITION DATA  
AUTHORIZED BY:** \_\_\_\_\_

Print Name \_\_\_\_\_ \*\* Signature - Required \*\*

Please provide contact phone number for confirmation \_\_\_\_\_

**RETURN TO CADHAM PROVINCIAL LABORATORY, INFORMATION MANAGEMENT**

**FAX 948-2124      Attention:** \_\_\_\_\_

Ubjacks\forms\reuest to amend req 948-2124.doc

**Manitoba**



Health  
Healthy Living

Cadham Provincial Laboratory  
Public Health Branch

P.O. Box 8450  
750 William Avenue  
Winnipeg MB R3C 3Y1  
PH: (204) 945-6123  
FAX: (204) 786-4770

**REQUEST FOR ACCESS TO INDIVIDUAL'S  
PERSONAL HEALTH INFORMATION**

I, \_\_\_\_\_  
Name of Applicant Date of Birth

\_\_\_\_\_  
Address

\_\_\_\_\_  
Telephone # PHIN #

hereby authorize \_\_\_\_\_  
Branch / Program Area

of Manitoba Health to disclose \_\_\_\_\_  
Extent or Nature of Information to be Released

to \_\_\_\_\_  
Name of Individual to Receive Information

\_\_\_\_\_  
Date of Consent

\_\_\_\_\_  
Signature of Individual (or Legal Representative)  
whose Information is to be Released

\_\_\_\_\_  
Expiry Date of Consent

\_\_\_\_\_  
Signature of Witness

**Manitoba**



Health  
Healthy Living

Cadham Provincial Laboratory  
Public Health Branch

P.O. Box 8450  
750 William Avenue  
Winnipeg MB R3C 3Y1  
PH: (204) 945-6123  
Fax: (204) 786-4770

[www.gov.mb.ca/health/publichealth/cpl](http://www.gov.mb.ca/health/publichealth/cpl)

**CONSENT FOR DISCLOSURE OF  
PERSONAL HEALTH INFORMATION**

I, \_\_\_\_\_  
Name of Applicant Date of Birth

\_\_\_\_\_  
Address

\_\_\_\_\_  
Telephone #

\_\_\_\_\_  
PHIN #

Request access to \_\_\_\_\_  
Extent or Nature of Information to be Released

\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_

From \_\_\_\_\_ of Manitoba.  
Branch/Program Area

Note: If the information requested involves laboratory test results, the information collected will be provided to your physician in order to provide a correct clinical interpretation to you. Please provide your physician's name and mailing address below:


\_\_\_\_\_  
\_\_\_\_\_

\_\_\_\_\_  
Date of Request

\_\_\_\_\_  
Signature of Individual (or Legal Representative)

If legal representative, what is relationship to individual: \_\_\_\_\_

W:\ELTCPLAdministration\FORMS\PHIA Consent for disclosure .doc

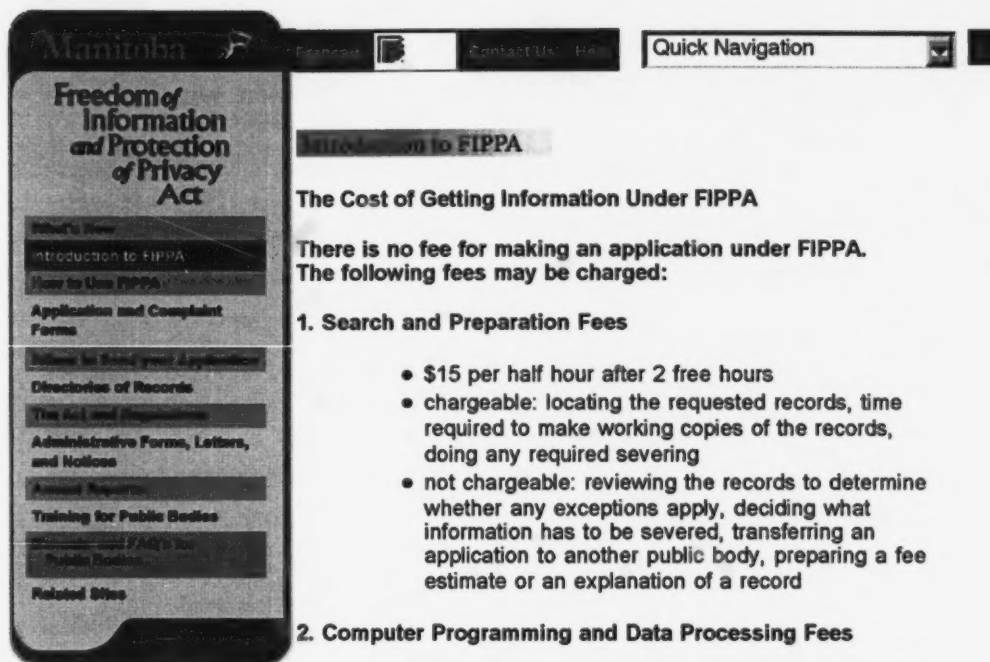

## Introduction to FIPPA

### The Cost of Getting Information Under FIPPA

There is no fee for making an application under FIPPA. The following fees may be charged:

#### 1. Search and Preparation Fees

- \$15 per half hour after 2 free hours
- chargeable: locating the requested records, time required to make working copies of the records, doing any required severing
- not chargeable: reviewing the records to determine whether any exceptions apply, deciding what information has to be severed, transferring an application to another public body, preparing a fee estimate or an explanation of a record

#### 2. Computer Programming and Data Processing Fees

- \$10 for each 15 minutes of internal programming or data processing, or the actual cost incurred when the work is done by an external agency

#### 3. Copying Records (If Applicant Requests a Copy)

- photocopies and computer printouts: 20 cents per page
- prints from microfilm: 50 cents per page
- any other copying method: actual cost

Note: applicants requesting copies of their own personal information are not required to pay for the copies if the total copying charge is less than \$10.

#### 4. Delivery Fees

- regular mail: no charge
- courier delivery: actual cost

### Fee Estimates

If you will be required to pay fees, the head of the public body must send you a fee estimate before providing the service. You have up to 30 days from the date the estimate is given to advise the public body if the estimate is accepted or to modify the request in order to reduce the fee. If you want the public body to proceed with the original request, payment for the total amount should be submitted, along with a signed copy of the Estimate of Costs form.

**CADHAM PROVINCIAL PUBLIC HEALTH LABORATORY DATA REQUEST FORM**

Submit Completed form to: Information Coordinator, Cadham Provincial Public Health Laboratory.

Ph: 545-3334 Fax: 545-3334

**Requested by:****Name:****Telephone:****Organization:****Fax:****Address:****Email:****Data Format Required : Tabulated Data**☐**ASCII File**☐**Other Software Format:** (Please specify:)**Description:** Please be as specific as possible. (Specify the data to be displayed, time period, gender, age, how the data is to be grouped/ordered, what summaries are to be displayed, e.g., counts, totals, subtotals, etc.)**Date of Request:** \_\_\_\_\_ **Date Required:** \_\_\_\_\_**Justification:** (Brief description of proposed use of data) \_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_**Ethics Committee Approval Obtained:** Yes☐

No

☐

(If Yes, please specify name of committee &amp; date of approval)

**CADHAM LABORATORY USE ONLY**

Health Information Privacy Committee Approval (if required) Yes \_\_\_\_\_ No \_\_\_\_\_

Rec'd. By: \_\_\_\_\_ Date: \_\_\_\_\_ Approved by: \_\_\_\_\_ Date: \_\_\_\_\_

Laboratory Manager

Specialist: \_\_\_\_\_ Date of Notification: \_\_\_\_\_


Project Manager (in-house/info systems) \_\_\_\_\_

Estimated completion time (hrs.) \_\_\_\_\_ Actual Completion time (hrs.) \_\_\_\_\_

Date Completed \_\_\_\_\_ Date sent report to client: \_\_\_\_\_

Report received from client Yes \_\_\_\_\_ No \_\_\_\_\_



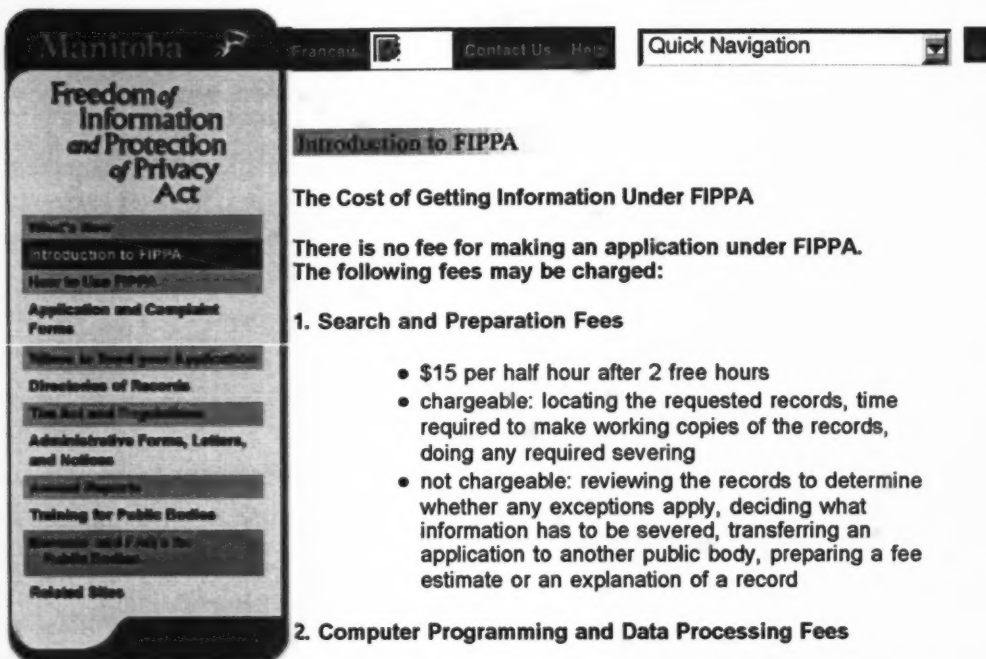

## Introduction to FIPPA

### The Cost of Getting Information Under FIPPA

There is no fee for making an application under FIPPA. The following fees may be charged:

#### 1. Search and Preparation Fees

- \$15 per half hour after 2 free hours
- chargeable: locating the requested records, time required to make working copies of the records, doing any required severing
- not chargeable: reviewing the records to determine whether any exceptions apply, deciding what information has to be severed, transferring an application to another public body, preparing a fee estimate or an explanation of a record

#### 2. Computer Programming and Data Processing Fees

- \$10 for each 15 minutes of internal programming or data processing, or the actual cost incurred when the work is done by an external agency

#### 3. Copying Records (If Applicant Requests a Copy)

- photocopies and computer printouts: 20 cents per page
- prints from microfilm: 50 cents per page
- any other copying method: actual cost

Note: applicants requesting copies of their own personal information are not required to pay for the copies if the total copying charge is less than \$10.

#### 4. Delivery Fees

- regular mail: no charge
- courier delivery: actual cost

### Fee Estimates

If you will be required to pay fees, the head of the public body must send you a fee estimate before providing the service. You have up to 30 days from the date the estimate is given to advise the public body if the estimate is accepted or to modify the request in order to reduce the fee. If you want the public body to proceed with the original request, payment for the total amount should be submitted, along with a signed copy of the Estimate of Costs form.

**CADHAM PROVINCIAL PUBLIC HEALTH LABORATORY DATA REQUEST FORM**

Submit Completed form to: Information Coordinator, Cadham Provincial Public Health Laboratory,  
Ph: 945-6624, Fax: 945-2124

**Requested by:****Name:****Telephone:****Organization:****Fax:****Address:****Email:****Data Format Required : Tabulated Data**☐**ASCII File**☐**Other Software Format: (Please specify:)**

**Description:** Please be as specific as possible. (Specify the data to be displayed, time period, gender, age, how the data is to be grouped/ordered, what summaries are to be displayed, e.g., counts, totals, subtotals, etc.)

**Date of Request:**  **Date Required:** **Justification: (Brief description of proposed use of data)** **Ethics Committee Approval Obtained: Yes**☐**No**☐

(If Yes, please specify name of committee & date of approval)

**CADHAM LABORATORY USE ONLY**Health Information Privacy Committee Approval (if required) Yes ☐ No ☐Rec'd. By: Date: Approved by: Date: 

Laboratory Manager

Specialist: Date of Notification: Project Manager (in-house/info systems) Estimated completion time (hrs.) Actual Completion time (hrs.) Date Completed Date sent report to client 

Report received from client

Yes ☐No ☐

**Cadham Provincial Laboratory**  
**General Requisition**



ONLY ONE SPECIMEN TYPE PER REQUISITION  
All areas of the requisition must be completed (please print clearly)  
See back for requisition/specimen instructions

Cadham Provincial Laboratory  
P.O. Box 9450  
750 Wilton Avenue  
Winnipeg, MB, R3C 3T1

Tel: (204) 945-6123  
Fax: (204) 736-4770  
E-mail: cadham@gov.mb.ca  
Website: www.gov.mb.ca/health/cadham/

ADDRESSOGRAPH

<b>RELEVANT CLINICAL INFORMATION</b>		<b>PATIENT INFORMATION</b>	
Outbreak Code: <input type="checkbox"/> In-Patient <input type="checkbox"/> Out-Patient <input type="checkbox"/> Uninsured Travel/Treatment History: <input type="checkbox"/> Autopsy <input type="checkbox"/> Diabetes <input type="checkbox"/> Food Borne Illness <input type="checkbox"/> Cancer/Chemotherapy <input type="checkbox"/> Dialysis <input type="checkbox"/> Pregnant Signs and Symptoms: <input type="checkbox"/> Bronchitis <input type="checkbox"/> Fever <input type="checkbox"/> Lymphadenopathy <input type="checkbox"/> Conjunctivitis <input type="checkbox"/> Gastrointestinal <input type="checkbox"/> Pneumonia <input type="checkbox"/> Chest Pain <input type="checkbox"/> Headache <input type="checkbox"/> Rash <input type="checkbox"/> Diarrhea <input type="checkbox"/> Influenza-Like Illness <input type="checkbox"/> Sore Throat <input type="checkbox"/> Encephalitis <input type="checkbox"/> Jaundice Other: _____ Reason for Test: <input type="checkbox"/> Immigration <input type="checkbox"/> Occupational <input type="checkbox"/> Other: _____ <input type="checkbox"/> Needlestick <input type="checkbox"/> Sexual Assault <input type="checkbox"/> Mucosal		PHN: _____ MB Health Reg. # _____ Alternate ID: <input type="checkbox"/> RCMP # <input type="checkbox"/> Other Province/Territories <input type="checkbox"/> Military # <input type="checkbox"/> Other _____ Date of Birth: _____ Sex: <input type="checkbox"/> M <input type="checkbox"/> F <input type="checkbox"/> U <input type="checkbox"/> A Chart/Clinic/Lab # _____ Patient Legal Name: Last _____ First Name _____ Street or Other (e.g., General Delivery) _____ Phone # _____ City/Municipality/First Nations Reserve _____ Postal Code _____	
<b>SPECIMEN INFORMATION</b> Specimen Type: _____ Specimen Source: _____ Collected At: _____ Date/Time: _____ YYYY/MM/DD HH:MM		<b>RETURN INFORMATION</b> Ordering Practitioner Last Name _____ First _____ Initial(s) _____ Facility _____ Facility Address _____ City/Town _____ Postal Code _____ Phone # _____ Secure Fax # _____	
<b>COPIES REPORT TO:</b> Other Practitioner Last Name _____ First Name _____ Facility _____ Secure Fax # _____			

<b>SEROLOGY</b>		<b>IMMUNOLOGY</b>	
Serology Test Panels <sup>(1)</sup> <input type="checkbox"/> Blood Borne Pathogen <input type="checkbox"/> Prenatal Panel <input type="checkbox"/> STI Panel <input type="checkbox"/> Prenatal HIV OPT OUT <sup>(2)</sup> <input type="checkbox"/> Post Exposure: Exposed Panel <sup>(1)</sup> <input type="checkbox"/> Post Exposure: Source Panel <sup>(1,3)</sup> Retrovirus <sup>(4)</sup> Syphilis <input type="checkbox"/> Screen <input type="checkbox"/> DFA <input type="checkbox"/> HIV1/2Ab <input type="checkbox"/> HTLV1/2Ab Hepatitis <input type="checkbox"/> HAV IgG (Immunity) <input type="checkbox"/> HBV PCR/Quant <input type="checkbox"/> WNV <input type="checkbox"/> HAV IgM <input type="checkbox"/> HBsAb (Immunity) <input type="checkbox"/> HCV PCR/Quant <input type="checkbox"/> HBcAb (IgM) <input type="checkbox"/> HBsAg <input type="checkbox"/> HCV Genotyping <input type="checkbox"/> HBcAb (Total) <input type="checkbox"/> HCV Ab Miscellaneous Serology CMV <input type="checkbox"/> IgM <input type="checkbox"/> IgG Parvo B19 <input type="checkbox"/> IgM <input type="checkbox"/> IgG EBV <input type="checkbox"/> IgM <input type="checkbox"/> IgG Rubella <input type="checkbox"/> IgM <input type="checkbox"/> IgG HSV <input type="checkbox"/> IgM <input type="checkbox"/> IgG Toxoplasma <input type="checkbox"/> IgM <input type="checkbox"/> IgG Measles <input type="checkbox"/> IgM <input type="checkbox"/> IgG Varicella <input type="checkbox"/> IgM <input type="checkbox"/> IgG Mumps <input type="checkbox"/> IgM <input type="checkbox"/> IgG WNV <input type="checkbox"/> IgM <input type="checkbox"/> Lyme Ab <input type="checkbox"/> Mycoplasma pneumoniae IgM		<input type="checkbox"/> Ova & Parasites <input type="checkbox"/> Pinworm Examination <input type="checkbox"/> Blood Smears <input type="checkbox"/> Skin Scrapings <input type="checkbox"/> Identification <b>MICROBIOLOGY/PATHOLOGY</b> <input type="checkbox"/> C&S <input type="checkbox"/> Chlamydia and Gonorrhea (NAAT) <input type="checkbox"/> MRSA Screen <input type="checkbox"/> Chlamydia DFA (Microtrak) <input type="checkbox"/> VRE Screen <input type="checkbox"/> Fungus Culture <input type="checkbox"/> Clostridium difficile Toxin <input type="checkbox"/> GBS Prenatal Screen <input type="checkbox"/> Verotoxin Testing <input type="checkbox"/> Spore/Sterilizer Testing Referral isolate: <input type="checkbox"/> Identification <input type="checkbox"/> Susceptibility Testing <input type="checkbox"/> Subtyping Isolate Information: _____ <b>PHYS DETECTION</b> <input type="checkbox"/> Viral Studies <input type="checkbox"/> CMV PCR (NAAT) <input type="checkbox"/> HSV PCR (NAAT) <b>OTHER TESTS OR REQUESTS</b>	

Name \_\_\_\_\_  
PHN \_\_\_\_\_

Name \_\_\_\_\_  
PHN \_\_\_\_\_

Specimen label stickers. Where necessary, please fill one in and affix to the accompanying specimen container.

MB-616 (Rev. 05/01)

#### REQUISITION/REMOVAL INFORMATION

**Mandatory Fields:** The specimen will not be tested until all mandatory fields (PHIN, Patient Legal Name, Date of Birth, Sex, Practitioner Name and Address, Source/Type) are provided.

**Alternate ID:** A unique health ID issued by other authorities such as: RCMP, Military, FNIH, Other Canadian Provinces, Great-West Life, etc. If no PHIN this is a mandatory field

**Sex:** M = Male; F = Female; U = Unknown, A = Ambiguous (Transgender)

#### REPORTS

**Secure Fax Number:** The fax machine must be in a secure location accessible ONLY to persons requiring reports.

**Report Address:** The address where the report(s) will be sent. Complete information including facility name is required to ensure delivery. All reports will be sent by fax unless otherwise indicated.

**Copy Report To:** This area can only be filled out or authorized by the ordering practitioner and is intended for another practitioner providing care.

#### REQUISITION TEST ORDERING INFORMATION

**Outbreak Code:** For Infection Control and Public Health Purposes call Outbreak Coordinator (Microbiology Scientist) for code at (204) 945-7473.

**Specimen Type:** The nature of the specimen (e.g., aspirate, blood, tissue/biopsy, stool, swab, urine, sputum, serum, plasma, CSF, etc.)

**Specimen Source:** The anatomical location or site (e.g., throat, right leg wound, etc.) from where the specimen was taken.

- 1) **Test Panels (T):** Prenatal (HBsAg, Rubella IgG, Syphilis, HIV1/2 Ab); Bloodborne (HBsAg, HCV Ab, HIV1/2 Ab); Serology STI (HBsAg, HBsAb, Syphilis, HIV1/2 Ab); Post Exposure – Exposed (HBsAg, HCV Ab, HIV1/2 Ab, HBsAb); Post Exposure – Source (HBsAg, HCV Ab, HIV1/2 Ab).
- 2) **HIV Opt Out Box (P):** When this box is checked off, HIV antibody testing will not be conducted as part of the panel.
- 3) **Post Exposure Panels (P):** If T55 protocol is required, list T55 in the "Other" space under Reason for Test on the front of this form. The "Other" space can also be used if this testing is required due to a bite.
- 4) **HIV (Retrovirus) (R):** For Non-Nominal HIV testing please use requisition MG #13306. For HIV Viral Load and Genotyping use requisition MG #5126 (Retrovirus Nucleic Acid Testing)
- 5) **Nucleic Acid (R):** (Viral load) Send 10 cc EDTA whole blood (must be received within 6 hours at CPL) or EDTA plasma (stored at 2-8°C and received within 3 days at CPL). Please record on the front of this requisition the date and time of collection.

#### SPECIMEN COLLECTION INFORMATION

**Specimen Labelling:** Label specimen (blank stickers found on the front of this requisition may be used for this purpose) with patient's full name and PHIN or alternate ID.

**Serology Specimen Volume Requirement:** 10 mL serum separator tube (full draw).

**Chlamydia and Gonorrhea (NAAT Testing):**

**Endocervical Swab Specimens:** The cervical swab remains the specimen of choice for women. Use the Gen-Probe Aptima Uniswix Swab Collection Kit.

**Urine:** Urine is the specimen of choice for males. It is the only recommended genital specimen for women without a cervix (hysterectomy) or those refusing a complete genital examination. The patient should not have urinated for at least one hour prior to sample collection. Use the Gen-Probe Aptima Urine Specimen Collection Unit.


**Male Urethral Swab Specimens:** Use the Gen-Probe Aptima Uniswix Swab Collection Kit.

The following specimens are unsuitable for processing: 1) Urine specimens received with liquid levels not between the two black lines; 2) Swab specimen transport tubes containing no swab, the cleaning swab (white shaft), two swabs, or a swab not supplied by Gen-Probe; 3) Urine or Uniswix Swab tubes with the foil cap missing or pierced; 4) Urine or swab specimens in tubes other than the Gen-Probe Aptima Collection tubes.

**Chlamydia DFA:** Use for throat, rectal, eye, nasopharyngeal specimens.

For detailed specimen submission requirements and rejection policy, please consult our *Guide to Services* available online at: <http://www.gov.mb.ca/health/publichealth/cpl/documents.html> or call 204-945-6806 to order a copy of the Guide.

Information collected on this form is for public health surveillance purposes.  
Send properly packaged specimen and requisition to:

**Cadham**  
**Provincial Laboratory** **Manitoba**   
750 William Avenue  
Winnipeg, MB R3E 3J7  
Tel: (204) 945-6123  
Fax: (204) 786-4770  
<http://www.gov.mb.ca/health/publichealth/cpl/index.html>

### Non-Nominal HIV Antibody Laboratory Requisition

All areas of this requisition must be filled in for complete lab service.

PATIENT INFORMATION									
PATIENT CODE									
LAST 2 LETTERS OF MOTHER'S MARRIAGE NAME	YEAR OF PATIENT'S BIRTH	DAY OF PATIENT'S BIRTH	SEX M or F	RHA * CODE OF PATIENT'S RESIDENCE (SEE RIGHT)	3 DIGIT PATIENT POSTAL CODE (FIRST 3)				
CHART/CLINIC #:									
Manitoba Health <input type="checkbox"/> INSURED <input type="checkbox"/> UNINSURED									

ORDERING PRACTITIONER	
SEND LABORATORY REPORT TO:	
PRACTITIONER FULL NAME & INITIALS	BILLING #
[Insert and seal this requisition in an envelope for the ordering lab use only]	
RETURN ADDRESS	FACILITY #
	PHONE #
CITY/TOWN	POSTAL CODE

**91234567**

MH 130 MG-133966 (Rev. 12/06)

Fold right  
edge of  
requisition  
to this line  
and seal to  
maintain  
patient  
confidentiality

#### RHA \* CODES

- 10 Winnipeg
- 15 Brandon
- 20 N. Eastman
- 25 S. Eastman
- 30 Interlake
- 40 Central
- 45 Assiniboine
- 60 Parkland
- 70 Norman
- 80 Burnwood
- 90 Churchill
- 95 Unknown
- 99 Non-resident

For help with selecting  
correct RHA, see map  
on back of this  
requisition.

**91234567**

SHADED AREA FOR LABORATORY USE ONLY	
ADMISSION #	DATE RECEIVED
SPECIMEN SOURCE	
DATE COLLECTED	
ANALYSIS REQUESTED	
<input type="checkbox"/> HIV ANTIBODY	
PREVIOUS HIV Ab TEST	
DATE YYYY MM DD RESULT (POS, NEG, INDET)	
REASONS FOR TESTING (Check all appropriate boxes)	
PRE <input type="checkbox"/> PREGNANT SWM <input type="checkbox"/> SEX WITH MEN SWW <input type="checkbox"/> SEX WITH WOMEN STW <input type="checkbox"/> SEX TRADE WORKER STI <input type="checkbox"/> PATIENT HAS A S.T.I. OR IS A CONTACT OF S.T.I. CASE ID <input type="checkbox"/> INJECTION DRUG USE OF <input type="checkbox"/> BORN TO HIV-INFECTED MOTHER ACP <input type="checkbox"/> OCCUPATIONAL EXPOSURE BLT <input type="checkbox"/> RECIPIENT OF BLOOD/BLOOD PRODUCTS PRIOR TO 1985 SAS <input type="checkbox"/> VICTIM OF SEXUAL ASSAULT IMM <input type="checkbox"/> IMMIGRATION REQUIREMENT GLO <input type="checkbox"/> EMIGRATED FROM HIV-ENDemic COUNTRY OTH <input type="checkbox"/> OTHER _____ REF <input type="checkbox"/> REFUSED BY PATIENT	
SHADED AREA FOR LABORATORY USE ONLY	

### Requisition Number

(Please affix to  
specimen container)

**91234567**

**91234567**

**91234567**

**91234567**

**91234567**

## MANITOBA MATERNAL SERUM SCREEN

G \_\_\_\_\_ P \_\_\_\_\_ W \_\_\_\_\_

HOSPITAL OF DELIVERY \_\_\_\_\_

## ETHNIC ORIGIN:

CAUCASIAN ☐ BLACK ☐  
 ABORIGINAL ☐ ASIAN ☐  
 ORIENTAL ☐ OTHER ☐

## FIRST DEGREE RELATIVE WITH:

CHROMOSOMAL ABNORMALITY YES ☐ NO ☐  
 SPINA BIFIDA/ANENCEPHALY YES ☐ NO ☐

EXPLAIN (CIRCLE RELATIONSHIP)

PATIENT / PARTNER / PREVIOUS CHILD / PARENT / SIB

## CURRENT PREGNANCY

CVS OR AMNIOCENTESIS YES ☐ NO ☐

LMP D M Y  
 360902

CYCLE \_\_\_\_\_ DAYS

CERTAIN YES ☐ NO ☐NORMAL YES ☐ NO ☐PILL WITHDRAWAL YES ☐ NO ☐

IF AN ULTRASOUND HAS BEEN PERFORMED:

DATE OF ULTRASOUND: \_\_\_\_\_

IF POSSIBLE ATTACH A COPY OF ULTRASOUND REPORT

SPD \_\_\_\_\_ CRL \_\_\_\_\_ OR

GESTATIONAL AGE BY US \_\_\_\_\_  
 (BPO/CRL MEASUREMENT PREFERRED)

EDC D M Y

EDC CALCULATED BY LMP / US / EXAM  
 (CIRCLE ONE)

CIGARETTE SMOKER YES ☐ NO ☐DIABETIC ON INSULIN OR MEDS YES ☐ NO ☐MULTIPLE PREGNANCY (TWINS/TRIPLETS) YES ☐ NO ☐REPEAT SAMPLE YES ☐ NO ☐

MOST RECENT WEIGHT \_\_\_\_\_ LB. \_\_\_\_\_ KG.

ACCURATE RESULTS CAN ONLY BE REPORTED WHEN ACCURATE INFORMATION IS PROVIDED.

## PATIENT INFORMATION (PLEASE PRINT)

NAME: LAST FIRST

ADDRESS

POSTAL CODE

D.O.B.

ALN.S.C.

P.M.I.N.

CLINIC CHART NO.

## PHYSICIAN:

ADDRESS:

PHONE: FAX:

## SPECIMEN PREPARATION

1. COLLECT 10 CC IN RED SERUM SEPARATOR TUBE
2. CENTRIFUGE AND FORWARD PRIMARY TUBE
3. ATTACH ONE REQUISITION LABEL TO THE TUBE
4. PRINT FIRST AND LAST NAME ON THE TUBE
5. STORE AT 4°C UNTIL SHIPPING

## FORWARD SPECIMEN TO:

CADHAM PROVINCIAL LAB  
 PUBLIC HEALTH CHEMISTRY  
 750 WILLIAM AVENUE  
 WINNIPEG, MANITOBA  
 R3C 3Y1

IF YOU HAVE QUESTIONS RE SAMPLE COLLECTION CALL (204) 945-8002

THIS IS A VOLUNTARY TEST. IF YOU REQUIRE  
 FURTHER INFORMATION REGARDING THIS TEST  
 PLEASE CALL THE CLINICAL CO-ORDINATOR MATERNAL SERUM  
 SCREENING PROGRAM AT (204) 787-2097.

A PROVINCIAL PROGRAM RUN JOINTLY BY CADHAM PROVINCIAL LAB,  
 SECTION OF GENETICS/METABOLISM, AND THE FETAL ASSESSMENT  
 UNIT.

## (TO BE COMPLETED AT TIME OF SAMPLE COLLECTION)

COLLECTED AT \_\_\_\_\_ FACILITY

PHLEBOTOMIST INITIALS

COLLECTION DATE D M Y

GESTATIONAL AGE AT COLLECTION DATE

360902

FORM NO. 854RS 11/88

360902

360902





CPL No. \_\_\_\_\_

### Manitoba Newborn Screening

MHSC \_\_\_\_\_ Mother's PHIN \_\_\_\_\_

Mother's name (first & last) \_\_\_\_\_

Mother's address \_\_\_\_\_

Mother's Permanent Residence (if different) \_\_\_\_\_

Infant's (last) name \_\_\_\_\_

Infant gender Male ☐ Female ☐ First ☐

Supplemental tests CPT1 ☐ GA1 ☐ Repeat ☐  
(Hutterite) (Oji-Cree)

Use 24 hour clock (i.e. 1:00 pm is 13:00 hrs)

Baby's birthdate 

Hr.	Min.	Day	Month	Year					

Date of sample 

Hr.	Min.	Day	Month	Year					

Sample taken by \_\_\_\_\_ (Print name)

Birth weight (grams) \_\_\_\_\_

Gestational age (weeks) \_\_\_\_\_ Premature: Yes ☐

Feeding: Breast milk ☐ Formula ☐ TPN ☐

NAME OF FOLLOW-UP PHYSICIAN/MIDWIFE FOR FAMILY/INFANT

Name \_\_\_\_\_

Address \_\_\_\_\_

Telephone \_\_\_\_\_

Addressograph

Information collected is for CPL Public Health Chemistry  
MG-8017 (Rev. 07/05)



### Instructions

- 1. Full-term baby** – Take sample at time of hospital discharge, regardless of age, but preferably over 24 hours of age. Babies sampled at less than 24 hours of age will require a repeat sample.
- 2. Premature or ill baby** – Take first specimen at five days of age and second specimen at two to three weeks of age or at time of hospital discharge, whichever comes first. Mark second specimen "Repeat".
- 3. Home/midwife assisted birth** – Take sample at two to three days of age.
- 4. Collection** – Sterilize skin of heel and puncture with disposable, no more than 2.0 mm lancet. If bleeding is slow, it is helpful to hold limb dependent for a short period before spotting blood on filter paper.
- 5. Handling** – Fill **all** circles with blood, apply from **one** side only. Let blood soak through and to the periphery of each circle. Allow to dry (4 hours) on a clean, dry surface at room temperature, **before** placing in plastic cover. Do not handle or contaminate blood spots. Ensure all information on newborn screening card is filled in.
- 6. Deliver immediately.**

To: Cadham Provincial Laboratory  
750 William Avenue/P.O. Box 8450  
Winnipeg, Manitoba R3C 3Y1

For more information call: (204) 945-7980

Filter Paper Lot #W041  
MG 8017 (Rev. 07/05)



89084108

Manitoba  
Health

# Retrovirus Nucleic Acid Testing

SEND SPECIMEN TO: CADHAM PROVINCIAL LABORATORY  
750 William Avenue,  
Winnipeg, Manitoba  
R3E 5J7 Tel. (204) 945-6123

PATIENT CODE NO.	LAST 3 LETTERS OF MEDICAL SIGNED NAME	YEAR OF PATIENT'S BIRTH	DAY OF BIRTH	SEX M OR F	REGIONAL HEALTH AUTHORITY (SEE LIST ON BACK) (SEE LIST ON BACK)	3 DIGIT POSTAL CODE
	R					

SPECIMEN DATE  
YR. | MO. | DAY  
2 | 0 |

CLINIC NUMBER

PHYSICIAN NO. | FACILITY NO. | TELEPHONE NO.

PHYSICIAN

FACILITY  
NAME  
AND  
ADDRESS

## REGIONAL HEALTH AUTHORITY (RHA) CODES

10 Winnipeg	40 Central	80 Burntwood
15 Brandon	50 Marquette	90 Churchill
20 North Eastman	55 South Westman	
25 South Eastman	60 Parkland	95 Unknown
30 Interlake	70 Norman	99 Non-resident

**NOTE** - In order to protect the confidentiality of the patient, this form has been designed to seal the Clinical and Epidemiological Data. Once the numbered stubs at the right have been removed, a seal will be revealed on the last copy. Remove the protective strip on this seal, fold and seal along margin at left.

PHIN →

TEAR OFF STUB(S)  
AND STICK  
AROUND TUBE(S)  
CONTAINING  
SPECIMEN(S)

## Type of Test Requested

Plasma Viral Load: ☐ Regular (400 - 750,000 Copies/ml)  
☐ Ultrasensitive (40 - 50,000 Copies/ml)  
Order if previous viral load <400 and patient adhering to therapy  
☐ Genotyping  
☐ Order if primary HIV is suspected and patient in window period

Provirus:

89084108

## Reason for Viral Load Testing:

- ☐ Pregnant  
☐ Initial Assessment  
☐ Three month follow-up  
☐ Monitoring Therapy (Four weeks after start of medication)  
☐ Change of Therapy (4 weeks following change of medication)  
☐ Primary Seroconversion Illness  
☐ Newborn Diagnosis

Date of Last Viral Load Test: \_\_\_\_\_

Most Recent CD4 Count: \_\_\_\_\_

Cells/ $\mu$ l

89084108

Date: \_\_\_\_\_

## Current Antiretrovirals:

### NRTI

- ☐ Zidovudine  
☐ Stavudine  
☐ Lamivudine  
☐ Didanosine  
☐ Zalcitabine  
☐ Abacavir  
☐ Adefovir

### NNRTI

- ☐ Nevirapine  
☐ Delavirdine  
☐ Efavirenz

### PI

- ☐ Saquinavir HGC or SGC  
☐ Ritonavir  
☐ Indinavir  
☐ Nelfinavir  
☐ Amprenavir  
☐ Lopinavir  
☐ Tipranavir

☐ Others: \_\_\_\_\_ (please specify)

89084108

## Specimen Collection:


- Collect 5 CC of Blood in EDTA tube.
- Send the specimen ASAP to the Serology Section of the Cadham Provincial Laboratory
- For viral load testing plasma should be separated and stored at -70C within 4 hr from time of blood collection

## NOTES

---

## NOTES

---

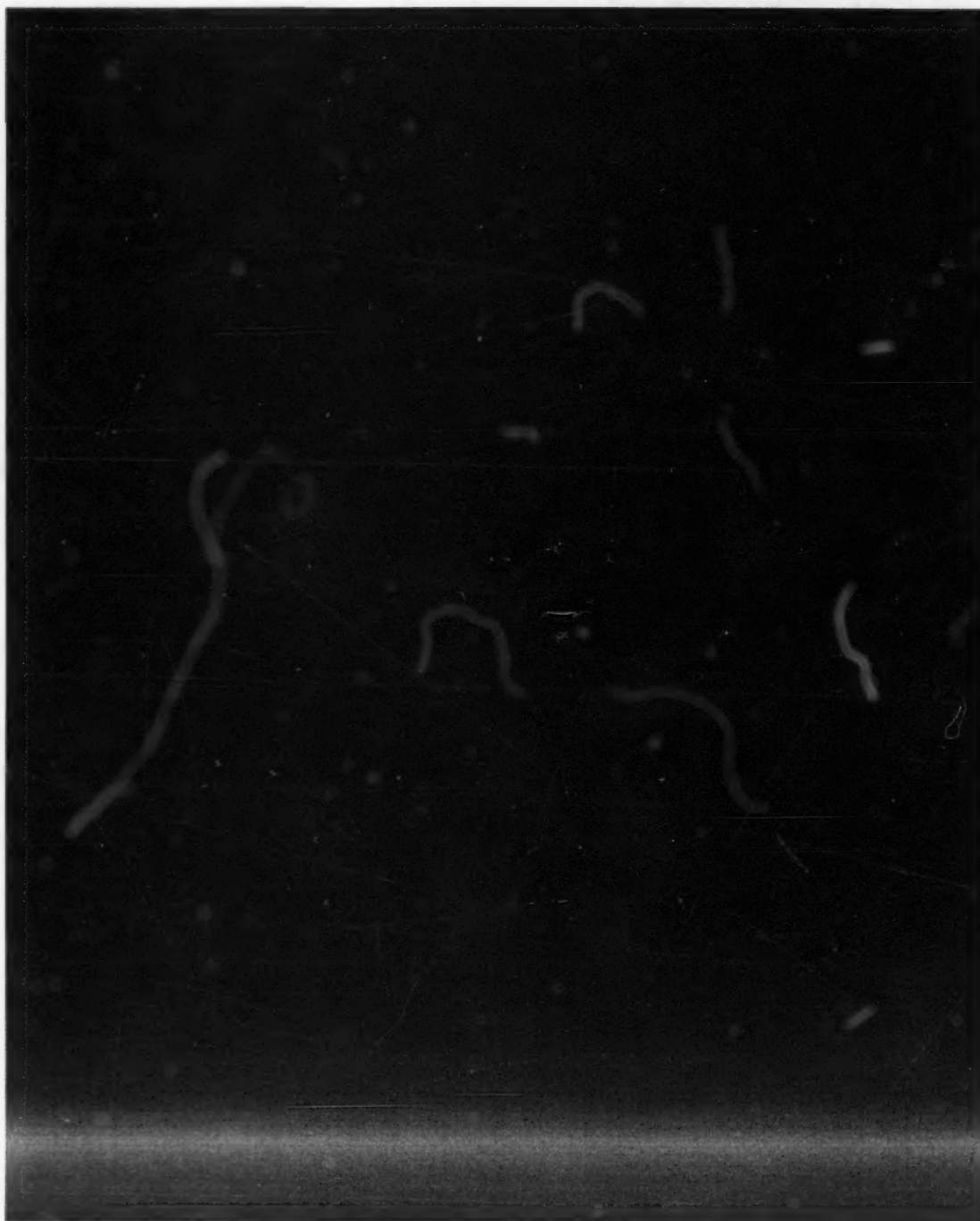

**CADHAM PROVINCIAL LABORATORY**

Manitoba Health